



CERTIFICATION by CUSTODIAN of RECORDS  
Tucson Police Department Crime Laboratory

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002201 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

TS4  
11-4-2008

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 10-23-08 Cal 1 0.099 AC Cal 2 0.097 AC

Between the dates of 9-25-08 and 10-23-08 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.  
The instrument was returned to service on N/A.

Clare Conley

Clare L. Conley

11-4-08

ADPS QAS Permit # 33314

SUBSCRIBED AND SWORN to before me this 4th day of November, 20 08.

[Signature]

Notary Public

\*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY  
STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000  
STANDARD QUALITY ASSURANCE PROCEDURE

1067  
Cu  
TS4  
11.4.2008

QA SPECIALIST Clare Conley AGENCY Tucson Police Department  
DATE 10-23-08 TIME 0943  
INTOXILYZER SERIAL # 80-002201 LOCATION Dui HQ ODW

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 097 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS None.  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE Clare Conley

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-002201  
10/23/2008

09:44:34

2 of 7 cu  
154

Standard Lot#: 730901I  
Last Changed By: J. OLIVARES #36552

QAS: C. CONLEY #34517  
TUCSON PD

Operator: C. CONLEY #34517  
TUCSON PD

Subject: TEST, SQAP  
DOB: 11/11/1111  
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:45:35
Diagnostic Test	Pass	09:46:03
Air Blank	0.000	09:46:30
0.100 Cal Check	0.099	09:46:49
Air Blank	0.000	09:47:18
Subject Test	0.000	09:47:34
Air Blank	ABT*	09:47:41
Air Blank	0.000	09:48:09

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-002201  
10/23/2008

09:48:13

3567  
Cu  
154

---

Standard Lot#: 730901I  
Last Changed By: J. OLIVARES #36552

---

QAS: C. CONLEY #34517  
TUCSON PD

---

Operator: C. CONLEY #34517  
TUCSON PD

---

Subject: TEST, SQAP  
DOB: 11/11/1111  
Sex: M Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	RFI*	09:48:50
Air Blank	IPS**	09:48:58

\*RFI Detect  
\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-002201  
10/23/2008

09:49:04

4567  
Cu  
104

Standard Lot#: 730901I  
Last Changed By: J. OLIVARES #36552

QAS: C. CONLEY #34517  
TUCSON PD

Operator: C. CONLEY #34517  
TUCSON PD

Subject: TEST, SQAP  
DOB: 11/11/1111  
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:49:52
Diagnostic Test	Pass	09:50:19
Air Blank	0.000	09:50:47
0.100 Cal Check	0.098	09:51:05
Air Blank	0.000	09:51:35
Subject Test	0.258*	09:54:43
Air Blank	0.000	09:55:17
Air Blank	0.000	09:55:44
Subject Test	XXX**	09:56:13
Air Blank	0.000	09:56:48

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-002201  
10/23/2008

09:57:51

-----  
Standard Lot#: 730901I  
Last Changed By: J. OLIVARES #36552  
-----

QAS: C. CONLEY #34517  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	09:58:18
Cal Check	0.097	09:58:37
Air Blank	0.000	09:59:06

5567  
cu  
154

*Claire Conley*

\_\_\_\_\_  
Operator's Signature

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6567  
Cu  
TS4 11-4-2008

QA Specialist Clare Conley AGENCY Tucson Police Department  
DATE 10-23-08 TIME 0943  
Intoxilyzer Serial # 80-002201 Location DUI HQ ODW

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 099 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Clare Conley

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-002201  
10/23/2008 09:43:06

-----  
Standard Lot#: 730901I  
Last Changed By: J. OLIVARES #36552  
-----

QAS: C. CONLEY #34517  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	09:43:34
Cal Check	0.099	09:43:52
Air Blank	0.000	09:44:22

7567  
cu  
754

C. Conley  
Operator's Signature



TUCSON POLICE DEPARTMENT  
CRIME LABORATORY

1 of 2  
en

INTOXILYZER SERVICE RECORD  
FOR  
CALIBRATION GAS TANK CHANGE

TSY  
11.4.2008

Date: 10-20-08  
 QAS: OLIVARES  
 Intoxilyzer Serial Number: 80-002201  
 New Lot Number: 730901I  
 New Expiration Date: 11-06-09

- (S) 1. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (S) 2. Remove the old gas tank and install the new one.
- (S) 3. Hit the "ESC" key twice and swipe the QAS card.
- (S) 4. Go to the "M-Maintenance" menu and select "S-Standard Change." Select "D" for dry gas and enter the following information:  
 Enter Target Value: 0.100  
 Enter Control Lot #: From the tank label  
 Enter Exp Date: From the tank label (mm/dd/yyyy)  
 Edit info Y/N? Y or N  
 Update Std Info? Y ENTER  
 "Please wait saving..."
- (S) 5. Hit the "ESC" key twice.
- (S) 6. Go to the "M-Maintenance" menu and select "T-Control Testing." Select "D" for dry control test. After selection is made press ENTER.
- (S) 7. Air blank completed.
- (S) 8. Calibration check completed. Test results 0. 098 AC.
- (S) 9. Air blank completed.
- (S) 10. Remove and sign printed record. Attach the record to the completed checklist.

SIGNATURE [Signature] 34552

2062  
cu

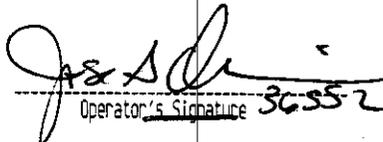
134

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-002201  
10/20/2008 18:48:31

-----  
Standard Lot#: 7309011  
Last Changed By: J. OLIVARES #36552  
-----

QAS: J. OLIVARES #36552  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	18:48:59
Cal Check	0.098	18:49:18
Air Blank	0.000	18:49:47

  
-----  
Operator's Signature 36552