



CERTIFICATION by CUSTODIAN of RECORDS  
 Tucson Police Department Crime Laboratory  
 270 S. Stone Avenue Tucson, AZ 85701-1917

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002201 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

*TS4  
2-12-09*

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 2-12-09 Cal 1 0.098 AC Cal 2 0.097 AC

Between the dates of 1-15-09 and 2-12-09 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.  
 The instrument was returned to service on N/A.

Clare Conley

Clare L. Conley


2-12-09

ADPS QAS Permit # 33314

SUBSCRIBED AND SWORN to before me this 12<sup>th</sup>, day of February, 2009.

A.A. Hernandez

Notary Public

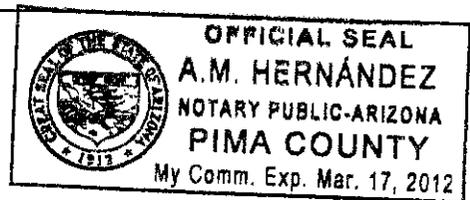

\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

1067  
Cu

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

TSY  
2-12-09

QA SPECIALIST Clare Conley AGENCY Tucson Police Department

DATE 2-12-09 TIME 1036

INTOXILYZER SERIAL # 8D-002201 LOCATION Crime Lab

( ✓ ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- ( ✓ ) 1. Clock time check.
- ( ✓ ) 2. Date check.

OPERATIONAL TESTS

- ( ✓ ) 1. Alcohol-free subject test result 0. 000 AC.
- ( ✓ ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- ( ✓ ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- ( ✓ ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 097 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS None.  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE Clare Conley

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-002201  
02/12/2009

10:37:44

267  
en

154

---

Standard Lot#: 730901I  
Last Changed By: J. OLIVARES #36552

---

QAS: C. CONLEY #34517  
TUCSON PD

---

Operator: C. CONLEY #34517  
TUCSON PD

---

Subject: TEST, SQAP  
DOB: 11/11/1111  
Sex: M Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	10:38:31
Diagnostic Test	Pass	10:38:59
Air Blank	0.000	10:39:26
0.100 Cal Check	0.098	10:39:45
Air Blank	0.000	10:40:14
Subject Test	0.000	10:40:30
Air Blank	ABT*	10:40:37
Air Blank	0.000	10:41:06

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-002201  
02/12/2009

10:41:11

3 of 7

cu

BA

Standard Lot#: 730901I  
Last Changed By: J. OLIVARES #36552

QAS: C. CONLEY #34517  
TUCSON PD

Operator: C. CONLEY #34517  
TUCSON PD

Subject: TEST, SQAP  
DOB: 11/11/1111  
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	10:41:44
Air Blank	IPS**	10:41:55

\*RFI Detect  
\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-002201  
02/12/2009

10:42:00

4 of 7  
Cu  
154

Standard Lot#: 730901I  
Last Changed By: J. OLIVARES #36552

QAS: C. CONLEY #34517  
TUCSON PD

Operator: C. CONLEY #34517  
TUCSON PD

Subject: TEST, SQAP  
DOB: 11/11/1111  
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:42:48
Diagnostic Test	Pass	10:43:15
Air Blank	0.000	10:43:43
0.100 Cal Check	0.097	10:44:02
Air Blank	0.000	10:44:31
Subject Test	0.149*	10:47:40
Air Blank	0.000	10:48:13
Air Blank	0.000	10:48:40
Subject Test	XXX**	10:49:08
Air Blank	0.000	10:49:42

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-002201  
02/12/2009 10:50:11

507  
lu  
TSH

-----  
Standard Lot#: 730901I  
Last Changed By: J. OLIVARES #36552  
-----

QAS: C. CONLEY #34517 -  
TUCSON PD -

Test	g/210L	Time
Air Blank	0.000	10:50:39
Cal Check	0.097	10:50:58
Air Blank	0.000	10:51:28

Clare Conley  
Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

607  
Cu  
134  
2-12-09

QA Specialist Clare AGENCY Tucson Police Department  
DATE 2-12-09 TIME 1036  
Intoxilyzer Serial # 80-002201 Location Crime Lab

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 098 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Clare Conley

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-002201  
02/12/2009 10:36:13

787  
cu

-----  
Standard Lot#: 730901I  
Last Changed By: J. OLIVARES #36552  
-----

TS4

-----  
QAS: C. CONLEY #34517  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	10:36:40
Cal Check	0.098	10:36:59
Air Blank	0.000	10:37:28

*Clare Conley*

Operator's Signature