



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

TS4
6-10-2009

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002201 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 06/04/09 Cal 1 0.099 AC Cal 2 0.099 AC

Between the dates of 05/07/09 and 06/04/09 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on —.
The instrument was returned to service on —.

Cheryl Straub

Cheryl K. Straub


06/10/09

ADPS QAS Permit # 33408

SUBSCRIBED AND SWORN to before me this 10th, day of June, 2009.

A.M. Hernandez

Notary Public

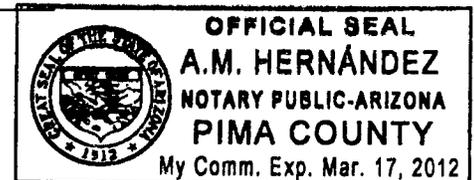

*Includes documents specified in A.A.C. R13-10-104.A.5

178 CG

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

TS4
6-9-2009

QA SPECIALIST Cheryl STRAUB AGENCY Tucson Police Department

DATE 06/04/09 TIME 08:11

INTOXILYZER SERIAL # 80-002201 LOCATION DUF HQ

() 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- () 1. Clock time check.
- () 2. Date check.

OPERATIONAL TESTS

- () 1. Alcohol-free subject test result 0. 000 AC.
- () 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- () 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- () 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 099 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Cheryl Straub

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

278-22

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002201
06/04/2009

08:17:46

BY

Standard Lot#: 804302I
Last Changed By: J. OLIVARES #36552

QAS: C. STRAUB #50772
TUCSON PD

Operator: C. STRAUB #50772
TUCSON PD

Subject: TEST, TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:18:48
Diagnostic Test	Pass	08:19:15
Air Blank	0.000	08:19:43
0.100 Cal Check	0.098	08:20:01
Air Blank	0.000	08:20:31
Subject Test	0.000	08:20:52
Air Blank	0.000	08:21:19
Five Minute Wait	ABT*	
Air Blank	0.000	08:22:11

*Sequence Aborted

Not a Successfully Completed Test Sequence

Pursuant to section 28-1527, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department

388
C8

TS4

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002201
06/04/2009 08:22:48

Standard Lot#: 804302I
Last Changed By: J. OLIVARES #36552

QAS: C. STRAUB #50772
TUCSON PD

Operator: C. STRAUB #50772
TUCSON PD

Subject: TEST,TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	08:23:26
Air Blank	IPS**	08:23:37

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

Pursuant to section 28-1527, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

498
CS

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002201
06/04/2009

08:24:15

134

Standard Lot#: 804302I
Last Changed By: J. OLIVARES #36552

QAS: C. STRAUB #50772
TUCSON PD

Operator: C. STRAUB #50772
TUCSON PD

Subject: TEST,TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:25:07
Diagnostic Test	Pass	08:25:35
Air Blank	0.000	08:26:02
0.100 Cal Check	0.098	08:26:21
Air Blank	0.000	08:26:51
Subject Test	XXX*	08:27:43
Air Blank	0.000	08:28:17

*Invalid Sample
Begin new deprivation period

Not a Successfully Completed Test Sequence

Pursuant to section 28-1527, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

588
CS

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002201
06/04/2009 08:30:59

TS4

Standard Lot#: 804302I
Last Changed By: J. OLIVARES #36552

QAS: C. STRAUB #50772
TUCSON PD

Operator: C. STRAUB #50772
TUCSON PD

Subject: TEST,TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:32:03
Diagnostic Test	Pass	08:32:30
Air Blank	0.000	08:32:58
0.100 Cal Check	0.098	08:33:16
Air Blank	0.000	08:33:46
Subject Test	0.000*	08:36:54
Air Blank	0.000	08:37:22
Air Blank	0.000	08:37:49
Subject Test	DEF**	08:40:58
Air Blank	0.000	08:41:26

*Deficient Sample
**Deficient Sample

Not a Successfully Completed Test Sequence

Pursuant to section 28-1527, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

698
CS

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002201
06/04/2009 08:42:32

TS4

Standard Lot#: 804302I
Last Changed By: J. OLIVARES #36552

QAS: C. STRAUB #50772
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	08:43:00
Cal Check	0.099	08:43:19
Air Blank	0.000	08:43:48

Cheyl Straub

Operator's Signature

788
CS

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

T54
6.9.09

QA Specialist Cheryl Straub AGENCY Tucson Police Department

DATE 06/04/09 TIME 08:11

Intoxilyzer Serial # 80-002201 Location DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 099 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Cheryl Straub

DPS Form Exh G-2 (Rev. 05-01)

848
CS

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002201
06/04/2009 08:11:05

154

Standard Lot#: 804302I
Last Changed By: J. OLIVARES #36552

QAS: C. STRAUB #50772
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	08:11:33
Cal Check	0.099	08:11:51
Air Blank	0.000	08:12:21

Cheryl Straub
Operator's Signature