



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
270 S. Stone Avenue Tucson, AZ 85701-1917

TS4  
6.29.09

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002201 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 6/29/09 Cal 1 0.098 AC Cal 2 0.097 AC

Between the dates of 6/4/09 and 6/29/09 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -  
The instrument was returned to service on -

Sarah Cione

Sarah L. Cione


6/29/09

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 29<sup>th</sup>, day of June, 2009.

A.M. Hernandez

Notary Public

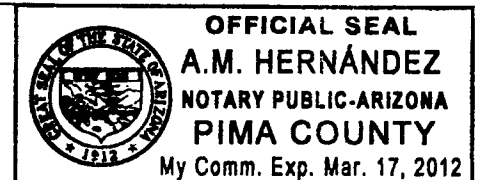

\*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY  
STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

## STANDARD QUALITY ASSURANCE PROCEDURE

1 of 9<sub>sc</sub>  
134  
6-29-09

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department  
DATE 6/29/09 TIME 10:19  
INTOXILYZER SERIAL # 80-002201 LOCATION Crime Lab

( ✓ ) 1. Display reads "PUSH BUTTON TO START"

## DIAGNOSTIC TESTS

( ✓ ) 1. Clock time check.  
( ✓ ) 2. Date check.

## OPERATIONAL TESTS

( ✓ ) 1. Alcohol-free subject test result 0. 000 AC.  
( ✓ ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed  
( ✓ ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.  
( ✓ ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 097 AC

Instrument is operating properly and accurately. YES ☒ NO ☐

## COMMENTS

\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE

Sarah Cione

DPS Form Exh G-4 (Rev. 05-01)

2 of 9<sub>sc</sub>

B4

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-002201  
06/29/2009

10:22:14

---

Standard Lot#: 804302I  
Last Changed By: J. OLIVARES #36552

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP,TEST  
DOB: 11/11/1111  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	10:22:59
Diagnostic Test	Pass	10:23:27
Air Blank	0.000	10:23:54
0.100 Cal Check	0.098	10:24:13
Air Blank	0.000	10:24:42
Subject Test	0.000	10:24:57
Air Blank	0.000	10:25:25
Five Minute Wait	ABT*	
Air Blank	0.000	10:25:59

\*Sequence Aborted

### Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

3 of 9<sub>sc</sub>

BY

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-002201  
06/29/2009 10:26:32

---

Standard Lot#: 804302I  
Last Changed By: J. OLIVARES #36552

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP,TEST  
DOB: 11/11/1111  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	RFI*	10:26:57
Air Blank	IPS**	10:27:05

\*RFI Detect

\*\*Improper Sample

### Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

4 of 9sc  
184

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-002201  
06/29/2009

10:27:11

Standard Lot#: 804302I  
Last Changed By: J. OLIVARES #36552

QAS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SQAP, TEST  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:28:18
Diagnostic Test	Pass	10:28:45
Air Blank	0.000	10:29:13
0.100 Cal Check	0.098	10:29:32
Air Blank	0.000	10:30:01
Subject Test	0.182*	10:33:10
Air Blank	0.000	10:33:42
Air Blank	0.000	10:34:09
Subject Test	XXX**	10:34:40
Air Blank	0.000	10:35:13

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

### Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

5 of 9 sc

INTOXILYZER 8000

Location: ARIZONA

Serial Number: 80-002201

06/29/2009

10:35:30

BH

-----  
Standard Lot#: 804302I

Last Changed By: J. OLIVARES #36552  
-----

QAS: S. CIONE #50616

TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	10:35:58
Cal Check	0.097	10:36:17
Air Blank	0.000	10:36:47

*Sarah Cione*

Operator's Signature

6 of 9 sc

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

TS4  
6-29-09

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department  
DATE 6/29/09 TIME 10:19  
Intoxilyzer Serial # 80-002201 Location Crime Lab

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach  $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 098 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

7 of 9 sc

INTOXILYZER 8000

Location: ARIZONA

Serial Number: 80-002201

06/29/2009

10:19:42

BH

-----  
Standard Lot#: 804302I

Last Changed By: J. OLIVARES #36552  
-----

QAS: S. CIONE #50616

TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	10:20:10
Cal Check	0.098	10:20:28
Air Blank	0.000	10:20:58

*Sarah Cione*

\_\_\_\_\_  
Operator's Signature



TUCSON POLICE DEPARTMENT  
CRIME LABORATORY

INTOXILYZER SERVICE RECORD  
FOR  
CALIBRATION GAS TANK CHANGE

8 of 9 sc

154  
6.29.09

Date: 6/29/09  
QAS: Sarah Cione  
Intoxilyzer Serial Number: 80-002201  
New Lot Number: 826001I  
New Expiration Date: 9/24/2010

- (✓) 1. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 2. Remove the old gas tank and install the new one.
- (✓) 3. Hit the "ESC" key twice and swipe the QAS card.
- (✓) 4. Go to the "M-Maintenance" menu and select "S-Standard Change."  
Select "D" for dry gas and enter the following information:  
Enter Target Value: 0.100  
Enter Control Lot #: From the tank label  
Enter Exp Date: From the tank label (mm/dd/yyyy)  
Edit info Y/N? Y or N  
Update Std Info? Y ENTER  
"Please wait saving..."
- (✓) 5. Hit the "ESC" key twice.
- (✓) 6. Go to the "M-Maintenance" menu and select "Control Testing."  
Select "D" for dry control test. After selection is made press  
ENTER.
- (✓) 7. Air blank completed.
- (✓) 8. Calibration check completed. Test results 0. 097 AC.
- (✓) 9. Air blank completed.
- (✓) 10. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Sarah Cione

9 of 9 sc

INTOXILYZER 8000

Location: ARIZONA

Serial Number: 80-002201

06/29/2009

10:39:56

154

-----  
Standard Lot#: 826001I

Last Changed By: S. CIONE #50616  
-----

QAS: S. CIONE #50616

TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	10:40:24
Cal Check	0.097	10:40:43
Air Blank	0.000	10:41:12

Sarah Cione  
Operator's Signature