



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002201 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 8.21.2009 Cal 1 0.097 AC Cal 2 0.097 AC

Between the dates of 7.23.2009 and 8.21.2009 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on .
The instrument was returned to service on .

Terry S. Gallegos

Terry S. Gallegos

9.2.2009


ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 2nd, day of September, 2009.

Amel Hernandez

Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

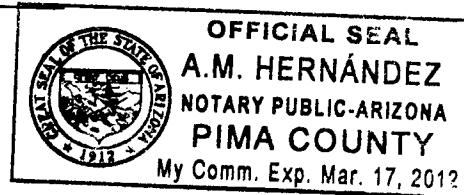

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

CU
9-2-09
10 of 7
TSQA SPECIALIST Terry Gallegos 20847 AGENCY Tucson Police DepartmentDATE 8.21.2009 TIME 1645INTOXILYZER SERIAL # 80-002201 LOCATION Crime Lab

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
(✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
(✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed ✓✓
(✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed ✓
Deficient sample printed. ✓
(✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 097 AC

Instrument is operating properly and accurately. YES ✓ NO

COMMENTS

SIGNATURE

Terry Gallegos

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002201
08/21/2009

16:46:56

Standard Lot#: 826001I
Last Changed By: S. CIONE #50616

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	16:47:44
Diagnostic Test	Pass	16:48:11
Air Blank	0.000	16:48:39
0.100 Cal Check	0.097	16:48:57
Air Blank	0.000	16:49:27
Subject Test	0.000	16:50:03
Air Blank	ABT*	16:50:15
Air Blank	0.000	16:50:43

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002201
08/21/2009

16:50:53

Standard Lot#: 826001I
Last Changed By: S. CIONE #50616

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	16:51:47
Diagnostic Test	Pass	16:52:14
Air Blank	RFI*	16:52:27
Air Blank	IPS**	16:52:36

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002201
08/21/2009

16:52:57

Standard Lot#: 826001I
Last Changed By: S. CIONE #50616

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	16:53:47
Diagnostic Test	Pass	16:54:14
Air Blank	0.000	16:54:41
0.100 Cal Check	0.098	16:55:00
Air Blank	0.000	16:55:30
Subject Test	0.092*	16:58:38
Air Blank	0.000	16:59:11
Air Blank	0.000	16:59:38
Subject Test	XXX**	17:00:07
Air Blank	0.000	17:00:42

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002201
08/21/2009

17:00:59

Standard Lot#: 826001I
Last Changed By: S. CIONE #50616

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Test	g/210L	Time
Air Blank	0.000	17:01:26
Cal Check	0.097	17:01:45
Air Blank	0.000	17:02:15


Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Terry Gallegos 20047 AGENCY Tucson Police Department
DATE 8.21.2009 TIME 1645
Intoxilyzer Serial # 80-002201 Location Crime Lab

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 097 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Terry Gallegos

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000

Location: ARIZONA

Serial Number: 80-002201

08/21/2009

16:45:23

Standard Lot#: 826001I

Last Changed By: S. CIONE #50616

QAS: T. GALLEGOS #28847

TUCSON PD CRIME LAB

Test g/210L Time

Air Blank 0.000 16:45:50
Cal Check 0.097 16:46:09
Air Blank 0.000 16:46:38


Operator's Signature