



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

BS4
5.4.2010

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002201 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 4-30-10 Cal 1 0.101 AC Cal 2 0.099 AC

Between the dates of 4-1-10 and 4-30-10 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.
The instrument was returned to service on N/A.

Clare Conley

Clare L. Conley

4-30-10


ADPS QAS Permit # 33314

SUBSCRIBED AND SWORN to before me this 30th, day of April, 20 10.

A.M. Hernandez

Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

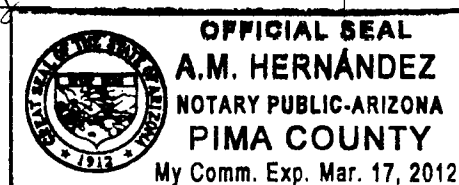

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Clare Conley AGENCY Tucson Police DepartmentDATE 4-30-10 TIME 0735INTOXILYZER SERIAL # 80-502201 LOCATION Crime Lab

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

(✓) 1. Clock time check.

(✓) 2. Date check.

OPERATIONAL TESTS

(✓) 1. Alcohol-free subject test result 0. 000 AC.(✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed(✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.(✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 094 ACInstrument is operating properly and accurately. YES ☒ NO ☐COMMENTS None.SIGNATURE Clare Conley

DPS Form Exh G-4 (Rev. 05-01)

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002201
04/30/2010

07:36:53

Standard Lot#: 826001I
Last Changed By: S. CIONE #50616

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	07:37:40
Diagnostic Test	Pass	07:38:08
Air Blank	0.000	07:38:35
0.100 Cal Check	0.101	07:38:54
Air Blank	0.000	07:39:23
Subject Test	0.000	07:39:41
Air Blank	ABT*	07:39:48
Air Blank	0.000	07:40:16

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002201
04/30/2010

07:40:22

357
CW
TS4

Standard Lot#: 826001I
Last Changed By: S. CIONE #50616

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP

DOB: 11/11/1111

Sex: M

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	07:40:48
Air Blank	IPS**	07:40:59

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002201
04/30/2010

07:41:04

4567
CW
154

Standard Lot#: 826001I
Last Changed By: S. CIONE #50616

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	07:41:51
Diagnostic Test	Pass	07:42:19
Air Blank	0.000	07:42:46
0.100 Cal Check	0.100	07:43:05
Air Blank	0.000	07:43:34
Subject Test	0.467*	07:46:43
Air Blank	0.000	07:47:16
Air Blank	0.000	07:47:43
Subject Test	XXX**	07:48:10
Air Blank	0.000	07:48:43

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000

Location: ARIZONA

Serial Number: 80-002201

04/30/2010

07:49:01


Standard Lot#: 826001I

Last Changed By: S. CIONE #50616

QAS: C. CONLEY #34517

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	07:49:28
Cal Check	0.099	07:49:47
Air Blank	0.000	07:50:17

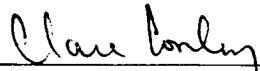

Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6587
CM

TS4
5.4.2010

QA Specialist Clare Conley AGENCY Tucson Police Department
DATE 4-30-10 TIME 0735
Intoxilyzer Serial # 80-002201 Location Crime Lab

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 101 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Clare Conley

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002201
04/30/2010

07:35:32

Standard Lot#: 826001I
Last Changed By: S. CIONE #50616

QAS: C. CONLEY #34517
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	07:35:59
Cal Check	0.101	07:36:18
Air Blank	0.000	07:36:48

Clare Conley

Operator's Signature