



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

SC
6/2/10

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002201 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 5.27.2010 Cal 1 0.100 AC Cal 2 0.099 AC

Between the dates of 4.30.2010 and 5.27.2010 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on _____.
The instrument was returned to service on _____.

Terry S. Gallegos

Terry S. Gallegos

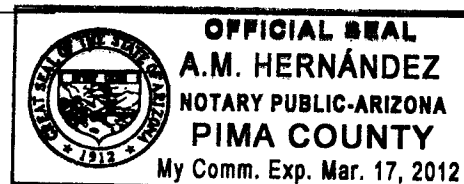
6.1.2010

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 1st, day of June, 2010.

A.M. Hernandez

Notary Public



*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1087
TS4

SC
6/2/10

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department

DATE 5.27.2010 TIME 1540

INTOXILYZER SERIAL # 80-002201 LOCATION Lab

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
(✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
(✓) 2. Error recognition logic system functioning. ✓✓✓
Not a Successfully Completed Test Sequence printed
(✓) 3. Proper sample recognition system. ✓✓
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
(✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 099 AC

Instrument is operating properly and accurately. YES ✓ NO

COMMENTS

SIGNATURE

Terry Gallegos

DPS Form Exh G-4 (Rev. 05-01)

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002201
05/27/2010

15:42:22

Standard Lot#: 826001I
Last Changed By: S. CIONE #50616

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Subject: TEST, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	15:43:14
Diagnostic Test	Pass	15:43:42
Air Blank	0.000	15:44:09
0.100 Cal Check	0.100	15:44:28
Air Blank	0.000	15:44:58
Subject Test	0.000	15:45:17
Air Blank	0.000	15:45:44
Five Minute Wait	ABT*	
Air Blank	0.000	15:46:18

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002201
05/27/2010

15:50:17

30f 7

154

SC

Standard Lot#: 826001I
Last Changed By: S. CIONE #50616

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Subject: TEST, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	15:51:09
Diagnostic Test	Pass	15:51:37
Air Blank	0.000	15:52:04
0.100 Cal Check	0.099	15:52:23
Air Blank	RFI*	15:52:33
Air Blank	IPS**	15:52:41

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002201
05/27/2010

15:52:55

Standard Lot#: 826001I
Last Changed By: S. CIONE #50616

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Subject: TEST, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	15:53:45
Diagnostic Test	Pass	15:54:12
Air Blank	0.000	15:54:40
0.100 Cal Check	0.100	15:54:59
Air Blank	0.000	15:55:28
Subject Test	0.175*	15:58:36
Air Blank	0.000	15:59:08
Air Blank	0.000	15:59:35
Subject Test	XXX**	16:00:00
Air Blank	0.000	16:00:33

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000

Location: ARIZONA

Serial Number: 80-002201

05/27/2010

16:00:56

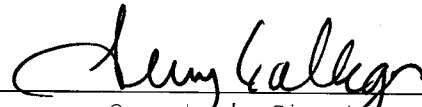
Standard Lot#: 826001I

Last Changed By: S. CIONE #50616

QAS: T. GALLEGOS #28847

TUCSON PD CRIME LAB

Test	g/210L	Time
Air Blank	0.000	16:01:23
Cal Check	0.099	16:01:42
Air Blank	0.000	16:02:12


Operator's Signature

5087
B4 SC

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department
DATE 5.27.2010 TIME 1540
Intoxilyzer Serial # 80-002201 Location Lab

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Terry Gallegos

DPS Form Exh G-2 (Rev. 05-01)

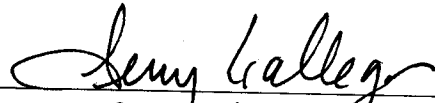
INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002201
05/27/2010

15:40:31

Standard Lot#: 826001I
Last Changed By: S. CIONE #50616

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Test	g/210L	Time
Air Blank	0.000	15:40:59
Cal Check	0.100	15:41:17
Air Blank	0.000	15:41:47


Operator's Signature

7057
TS4
SC