



# CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
270 S. Stone Avenue Tucson, AZ 85701-1917

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002201 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

TSH  
7.6.2010

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 07/01/10 Cal 1 0.098 AC Cal 2 0.098 AC

Between the dates of 05/27/10 and 07/01/10 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.  
The instrument was returned to service on N/A.

Cheryl K. Straub

Cheryl K. Straub


07/01/10

ADPS QAS Permit # 33408

SUBSCRIBED AND SWORN to before me this 1, day of July, 2010.

A.M. Hernandez

Notary Public

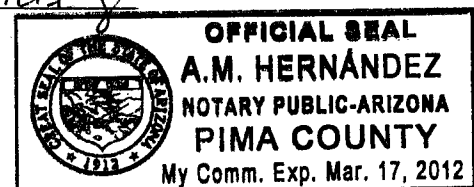

\*Includes documents specified in A.A.C. R13-10-104.A.5

177  
CS

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

TSY  
7.6.2010

QA SPECIALIST Cheryl STRAUB AGENCY Tucson Police Department

DATE 07/01/10 TIME 13:05

INTOXILYZER SERIAL # 80-002201 LOCATION CRIME LAB

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.  
(✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.  
(✓) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed  
(✓) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.  
(✓) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 098 AC

Instrument is operating properly and accurately. YES ☒ NO ☐

COMMENTS

SQAP EXPIRED 06/27/10 WAS NOT BEING  
USED FOR SUBJECT TEST. C8-50772

SIGNATURE

Cheryl Straub

277  
CS

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-002201  
07/01/2010

13:07:10

BY

Standard Lot#: 826001I  
Last Changed By: S. CIONE #50616

QAS: C. STRAUB #50772  
TUCSON PD

Operator: C. STRAUB #50772  
TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	13:08:00
Diagnostic Test	Pass	13:08:27
Air Blank	0.000	13:08:55
0.100 Cal Check	0.099	13:09:13
Air Blank	0.000	13:09:43
Subject Test	0.000	13:10:08
Air Blank	ABT*	13:10:16
Air Blank	0.000	13:10:44

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

387  
CS

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-002201  
07/01/2010

13:11:27

104

Standard Lot#: 826001I  
Last Changed By: S. CIONE #50616

QAS: C. STRAUB #50772  
TUCSON PD

Operator: C. STRAUB #50772  
TPD

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	13:12:27
Air Blank	IPS**	13:12:36

\*RFI Detect

\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

477  
CS

134

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-002201  
07/01/2010

13:13:00

Standard Lot#: 826001I  
Last Changed By: S. CIONE #50616

QAS: C. STRAUB #50772  
TUCSON PD

Operator: C. STRAUB #50772  
TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	13:13:59
Diagnostic Test	Pass	13:14:26
Air Blank	0.000	13:14:54
0.100 Cal Check	0.099	13:15:12
Air Blank	0.000	13:15:42
Subject Test	0.232*	13:18:50
Air Blank	0.000	13:19:22
Air Blank	0.000	13:19:50
Subject Test	XXX**	13:20:16
Air Blank	0.000	13:20:49

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

537  
CS

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-002201  
07/01/2010

13:21:05

-----  
Standard Lot#: 826001I  
Last Changed By: S. CIONE #50616  
-----

134

QAS: C. STRAUB #50772  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	13:21:33
Cal Check	0.098	13:21:52
Air Blank	0.000	13:22:21

  
\_\_\_\_\_  
Operator's Signature

EXHIBIT G-2

687  
CS

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

JSY  
7.6.2010

QA Specialist Cheryl STRAUB AGENCY Tucson Police Department  
DATE 07/01/10 TIME 13:05  
Intoxilyzer Serial # 80-002201 Location Crime Lab

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution N/A AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 098 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Cheryl Straub

DPS Form Exh G-2 (Rev. 05-01)

777  
CS

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-002201  
07/01/2010

13:05:29

SM

-----  
Standard Lot#: 826001I  
Last Changed By: S. CIONE #50616  
-----

QAS: C. STRAUB #50772  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	13:05:57
Cal Check	0.098	13:06:15
Air Blank	0.000	13:06:45

  
Operator's Signature