



CERTIFICATION by CUSTODIAN of RECORDS  
 Tucson Police Department Crime Laboratory  
 270 S. Stone Avenue Tucson, AZ 85701-1917

TSY  
 3-9-2011

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002201 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 03/03/2011 Cal 1 0.104 AC Cal 2 0.101 AC

Between the dates of 02/03/2011 and 03/03/2011 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.  
 The instrument was returned to service on N/A.

Cheryl K. Straub

Cheryl K. Straub


03/08/2011

ADPS QAS Permit # 33408

SUBSCRIBED AND SWORN to before me this 8, day of March, 2011.

A.M. Hernandez

Notary Public

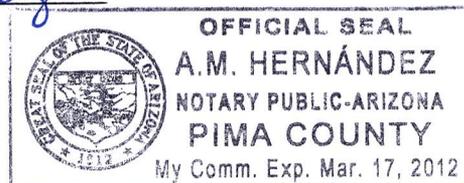

\*Includes documents specified in A.A.C. R13-10-104.A.5

179  
CS

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

TS4  
3-9-2011

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Cheryl STRAUB AGENCY Tucson Police Department

DATE 03/03/2011 TIME 07:54

INTOXILYZER SERIAL # 80-002201 LOCATION TEAM II DUI/HQ

(  ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (  ) 1. Clock time check.
- (  ) 2. Date check.

OPERATIONAL TESTS

- (  ) 1. Alcohol-free subject test result 0. 000 AC.
- (  ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- (  ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- (  ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 101 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS TEST #1; Alcohol-free Subject test, RFI occurred  
during the Air blank After the Calibration Check, Test  
Repeated. TEST #3, Deficient Sample, Ambient fail occurred  
during the Start up of test #3, test Repeated. CS 50772  
03/03/2011

SIGNATURE Cheryl Straub

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

289  
CS

TSM

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-002201  
03/03/2011 07:56:54

-----  
Standard Lot#: 20910100A1  
Last Changed By: C. CONLEY #34517  
-----

OAS: C. STRAUB #50772  
TUCSON PD  
-----

Operator: C. STRAUB #50772  
TUCSON PD  
-----

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111  
-----

15 Minute Deprivation Period? Yes  
-----

Test	g/210L	Time
Air Blank	0.000	07:57:57
Diagnostic Test	Pass	07:58:26
Air Blank	0.000	07:58:54
0.100 Cal Check	0.102	07:59:12
Air Blank	RFI*	07:59:35
Air Blank	0.000	08:00:03

\*RFI Detect

Not a Successfully  
Completed  
Test Sequence

389 CS

754

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-002201  
03/03/2011 08:00:16

-----  
Standard Lot#: 20910100A1  
Last Changed By: C. CONLEY #34517  
-----

QAS: C. STRAUB #50772  
TUCSON PD  
-----

Operator: C. STRAUB #50772  
TUCSON PD  
-----

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111  
-----

15 Minute Deprivation Period? Yes  
-----

Test	g/210L	Time
Air Blank	0.000	08:01:18
Diagnostic Test	Pass	08:01:47
Air Blank	0.000	08:02:14
0.100 Cal Check	0.102	08:02:33
Air Blank	0.000	08:03:02
Subject Test	0.000	08:03:22
Air Blank	ABT*	08:03:29
Air Blank	0.000	08:03:57

\*Sequence Aborted

Not a Successfully  
Completed  
Test Sequence

499  
CS

BY

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-002201  
03/03/2011 08:04:26

-----  
Standard Lot#: 20910100A1  
Last Changed By: C. CONLEY #34517  
-----

QAS: C. STRAUB #50772  
TUCSON PD  
-----

Operator: C. STRAUB #50772  
TUCSON PD  
-----

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111  
-----

15 Minute Deprivation Period? Yes  
-----

Test	g/210L	Time
Air Blank	RFI*	08:05:04
Air Blank	IPS**	08:05:13

\*RFI Detect  
\*\*Improper Sample

Not a Successfully  
Completed  
Test Sequence

589  
CS

TS4

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-002201  
03/03/2011 08:05:19

-----  
Standard Lot#: 20910100A1  
Last Changed By: C. CONLEY #34517  
-----

QAS: C. STRAUB #50772  
TUCSON PD  
-----

Operator: C. STRAUB #50772  
TUCSON PD  
-----

Subject: SOAP, SOAP  
DOB: 11/11/1111  
Sex: F Weight: 111  
-----

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	AMB*	08:06:15
Air Blank	0.000	08:06:43

\*Ambient Fail

Not a Successfully  
Completed  
Test Sequence

689  
CS

TS4

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-002201  
03/03/2011 08:07:18

-----  
Standard Lot#: 20910100A1  
Last Changed By: C. CONLEY #34517  
-----

QAS: C. STRAUB #50772  
TUCSON PD  
-----

Operator: C. STRAUB #50772  
TUCSON PD  
-----

Subject: SOAP, SOAP  
DOB: 11/11/1111  
Sex: F Weight: 111  
-----

15 Minute Deprivation Period? Yes  
-----

Test	g/210L	Time
Air Blank	0.000	08:08:14
Diagnostic Test	Pass	08:08:43
Air Blank	0.000	08:09:10
0.100 Cal Check	0.102	08:09:29
Air Blank	0.000	08:09:58
Subject Test	0.146*	08:13:07
Air Blank	0.000	08:13:39
Air Blank	0.000	08:14:06
Subject Test	XXX**	08:14:43
Air Blank	0.000	08:15:16

\*Deficient Sample  
\*\*Invalid Sample  
Begin new deprivation period

Not a Successfully  
Completed  
Test Sequence

789  
CS

TS4

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-002201  
03/03/2011 08:16:10

-----  
Standard Lot#: 20910100A1  
Last Changed By: C. CONLEY #34517  
-----

QAS: C. STRAUB #50772  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	08:16:38
Cal Check	0.101	08:16:56
Air Blank	0.000	08:17:26



-----  
Operator's Signature  
-----

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY  
STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

TS4  
3.9.2011

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Cheryl STRAUB AGENCY Tucson Police Department  
DATE 03/03/2011 TIME 07:54  
Intoxilyzer Serial # 80-002201 Location TEAM II DUI/HQ

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution N/A AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 104 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Cheryl Straub

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

989  
CS

TS4

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-002201  
03/03/2011 07:54:11

-----  
Standard Lot#: 20910100A1  
Last Changed By: C. CONLEY #34517  
-----

QAS: C. STRAUB #50772  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	07:54:38
Cal Check	0.104	07:54:57
Air Blank	0.000	07:55:27


Operator's Signature