



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002201 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

TS4
6-27-2011

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 6/9/11 Cal 1 0.104 AC Cal 2 0.103 AC

Between the dates of 5/26/11 and 6/9/11 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on 6/9/11
The instrument was returned to service on -

Sarah Cione

Sarah L. Cione


6/27/11

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 27, day of June, 2011.

A.M. Hernandez

Notary Public

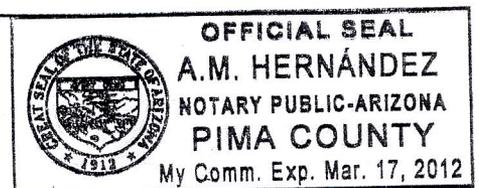

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

1 of 9 sc

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

TS4
6-27-2011

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department
DATE 6/9/11 TIME 1458
INTOXILYZER SERIAL # 80-002201 LOCATION Crime Lab

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 103 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Sarah Cione

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002201
06/09/2011

14:59:26

2 of 9 sc

B4

Standard Lot#: 20910100A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	15:00:12
Diagnostic Test	Pass	15:00:39
Air Blank	0.000	15:01:07
0.100 Cal Check	0.103	15:01:25
Air Blank	0.000	15:01:55
Subject Test	0.000	15:02:12
Air Blank	0.000	15:02:39
Five Minute Wait	ABT*	
Air Blank	0.000	15:03:15

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002201
06/09/2011

15:03:23

3 of 9 sc

BH

Standard Lot#: 20910100A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	15:04:07
Diagnostic Test	Pass	15:04:34
Air Blank	RFI*	15:04:53
Air Blank	IPS**	15:05:04

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002201
06/09/2011 15:05:08

4 of 9_{sc}

BS4

Standard Lot#: 20910100A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	15:05:55
Diagnostic Test	Pass	15:06:22
Air Blank	0.000	15:06:49
0.100 Cal Check	0.104	15:07:08
Air Blank	0.000	15:07:38
Subject Test	0.124*	15:10:46
Air Blank	0.000	15:11:18
Air Blank	0.000	15:11:46
Subject Test	XXX**	15:12:01
Air Blank	0.000	15:12:34

*Deficient Sample
**Invalid Sample
Begin new deprivation period

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002201
06/09/2011 15:13:00


5 of 9 sc

Standard Lot#: 20910100A1
Last Changed By: S. CIONE #50616

TS4

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	15:13:27
Cal Check	0.103	15:13:46
Air Blank	0.000	15:14:15

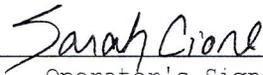

Operator's Signature

EXHIBIT G-2

6 of 9 sc

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

TSC
6-27-2011

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 6/9/11 TIME 1458
Intoxilyzer Serial # 80-002201 Location Crime Lab

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 104 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000

Location: ARIZONA

Serial Number: 80-002201

06/09/2011

14:58:04

7 of 9 SC

Standard Lot#: 20910100A1

Last Changed By: S. CIONE #50616

TS4

QAS: S. CIONE #50616

TUCSON PD

Test g/210L Time

Air Blank 0.000 14:58:31
Cal Check 0.104 14:58:50
Air Blank 0.000 14:59:19

Sarah Cione

Operator's Signature



TUCSON POLICE DEPARTMENT
CRIME LABORATORY

8 of 9 sc

INTOXILYZER SERVICE RECORD

TSY
6-27-2011

Date: 6/9/11

QAS: Sarah Cione

Intoxilyzer Serial Number: 80-002201

Reason For Service: Needs preventative maintenance. Calibration checks
(0.100) running consistently 3-4% high

Date Out of Service: 6/9/11

Calibration Check completed prior to removal from service

Date: 6/9/11

SQAP completed prior to removal from service

Date: 6/9/11

Corrective Action: Sent to CMI, Inc.



9 of 95c

Service Evaluation Form

This form **MUST** be completed and enclosed with instrument to be serviced.
Failure to complete and return this form may cause delays in service.

154

(Note: please ship items in their original shipping container(s) or a similar protective box.)

1. Contact information: Customer Number 105075 (contact Customer Service)
 Name Terry Gallegos Phone: (520) 837-7677
 Fax: (520) 791-4642 Email: Terry.Gallegos@tucsonaz.gov

2. Bill to Address: Ship to Address:
Finance / Accounts Payable Tucson Police Department Crime Lab
Po Box 27450 270 S. Stone Ave.
Tucson, AZ 85726-7450 Tucson, AZ 85701-1917

3. Serial Number: 80-002201 Instrument Model: I 8000

4. Detailed Description of Problem:
Calibration checks running 3-4% high. Preventative maintenance
needed.

Hazardous Material Warning! - DO NOT return gas cylinders with instrument!

I Authorize Repairs Up To: All \$250 \$500 \$750 Other \$ _____
 Purchase Order Number (attach a copy of P.O. if applicable) 34770

Authorized By:

Terry Gallegos Toxicology Coordinator
 Name (Please Print) Title
[Signature] 6-27-2011
 Signature Date

Ship to:
 CMI, Inc.
 Attn: Service Dept.
 316 East Ninth Street
 Owensboro, KY 42303

No, please send estimate before repairs are made.
 Note: An estimate will be faxed before performing any repairs and may cause delays in service.
 An evaluation fee (\$79.00 or actual costs) will apply to estimates that are not repaired.