



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

154
12-22-2011

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002201 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 12/9/11 Cal 1 0.101 AC Cal 2 0.101 AC

Between the dates of 11/10/11 and 12/9/11 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -

The instrument was returned to service on -

Sarah Cione

Sarah L. Cione


12/20/11

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 20, day of December, 20 11.

A.M. Hernandez

Notary Public

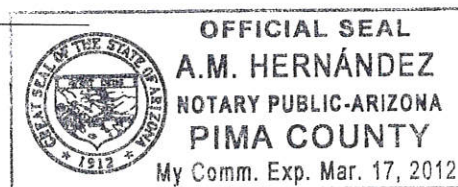

*Includes documents specified in A.A.C. R13-10-104.A.5

1 of 7 sc

TS4
12-22-2011

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 12/9/11 TIME 1610

INTOXILYZER SERIAL # 80-002201 LOCATION Crime Lab

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
(✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0.000 AC.
(✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
(✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
(✓) 4. Standard Calibration Check standard 0.100 AC
Result: 0.101 AC

Instrument is operating properly and accurately. YES ☒ NO ☐

COMMENTS

SIGNATURE Sarah Cione

2 of 7_{sc}

TS4

INTOXILYZER 8000

Location: TUCSON PD

Serial Number: 80-002201

Core Version: 8105.48

12/09/2011

16:12:18

Standard Lot#: 0491110A1

Last Changed By: S. CIONE

#50616

QAS: S. CIONE

#50616

TUCSON PD

Operator: S. CIONE

#50616

TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1980

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	16:13:08
Diagnostic Test	Pass	16:13:43
Air Blank	0.000	16:14:10
0.100 Cal Check	0.101	16:14:31
Air Blank	0.000	16:15:01
Subject Test	0.000	16:15:26
Air Blank	0.000	16:15:54
Wait	ABT*	
Air Blank	0.000	16:16:29

*Sequence Aborted

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

3 of 7_{sc}

BH

INTOXILYZER 8000

Location: TUCSON PD

Serial Number: 80-002201

Core Version: 8105.48

12/09/2011

16:16:41

Standard Lot#: 0491110A1

Last Changed By: S. CIONE

#50616

QAS: S. CIONE

#50616

TUCSON PD

Operator: S. CIONE

#50616

TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1980

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	16:17:09
Air Blank	IPS**	16:17:16

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

4 of 7_{sc}

B4

INTOXILYZER 8000
 Location: TUCSON PD
 Serial Number: 80-002201
 Core Version: 8105.48
 12/09/2011

16:17:22

Standard Lot#: 0491110A1
 Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
 TUCSON PD

Operator: S. CIONE #50616
 TUCSON PD

Subject: SQAP, SQAP
 DOB: 11/11/1980
 Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	16:18:06
Diagnostic Test	Pass	16:18:41
Air Blank	0.000	16:19:09
0.100 Cal Check	0.101	16:19:29
Air Blank	0.000	16:19:59
Subject Test	0.048*	16:23:07
Air Blank	0.000	16:23:38
Air Blank	0.000	16:24:06
Subject Test	XXX**	16:24:42
Air Blank	0.000	16:25:16

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000

Location: TUCSON PD

Serial Number: 80-002201

Core Version: 8105.48

12/09/2011

16:25:49

Standard Lot#: 0491110A1

Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616

TUCSON PD

Operator: S. CIONE #50616

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	16:26:18
0.100 Cal Check	0.101	16:26:39
Air Blank	0.000	16:27:09

Preventative Maintenance
Performed -- Timer Reset

5 of 7 sc

BH

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

BS4
12-22-2011

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 12/9/11 TIME 1610
Intoxilyzer Serial # 80-002201 Location Crime Lab

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0.100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0.101 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

7 of 7 sc

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002201
Core Version: 8105.48
12/09/2011

16:10:48

Standard Lot#: 0491110A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	16:11:18
0.100 Cal Check	0.101	16:11:39
Air Blank	0.000	16:12:08

Preventative Maintenance
Performed -- Timer Reset