



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

TS4
11-19-2012

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002201 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 11/9/12 Cal 1 0.102 AC Cal 2 0.101 AC

Between the dates of 10/12/12 and 11/9/12 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -
The instrument was returned to service on -

Sarah Cione

Sarah L. Cione


11/14/12

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 14th, day of November, 2012.

Elma Escarcega

Notary Public

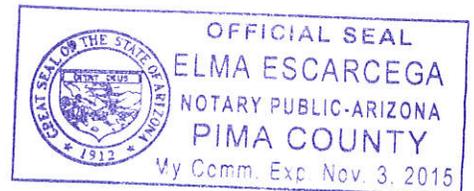

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

1 of 7sc

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

TSY
11-19-2012

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 11/9/12 TIME 1500

INTOXILYZER SERIAL # 80-002201 LOCATION Crime Lab

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

(✓) 1. Clock time check.

(✓) 2. Date check.

OPERATIONAL TESTS

(✓) 1. Alcohol-free subject test result 0. 000 AC.

(✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed

(✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.

(✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 101 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE

Sarah Cione

2 of 7 sc
B4

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002201
Core Version: 8105.48
11/09/2012

15:02:41

Standard Lot#: 23511100A2
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	15:03:25
Diagnostic Test	Pass	15:04:00
Air Blank	0.000	15:04:27
0.100 Cal Check	0.101	15:04:48
Air Blank	0.000	15:05:18
Subject Test	0.000	15:05:52
Air Blank	0.000	15:06:19
Wait	ABT*	
Air Blank	0.000	15:06:57

*Sequence Aborted

Not a Successfully Completed Test Sequence

3 of 7_{sc}
BH

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002201
Core Version: 8105.48
11/09/2012 . 15:07:07

Standard Lot#: 23511100A2
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	15:07:50
Diagnostic Test	Pass	15:08:25
Air Blank	0.000	15:08:53
0.100 Cal Check	0.102	15:09:14
Air Blank	0.000	15:09:44
Subject Test	RFI*	15:09:56
Air Blank	IPS**	15:10:03

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

4 of 7 sc

134

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002201
Core Version: 8105.48
11/09/2012

15:10:07

Standard Lot#: 23511100A2
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	15:10:49
Diagnostic Test	Pass	15:11:24
Air Blank	0.000	15:11:51
0.100 Cal Check	0.101	15:12:12
Air Blank	0.000	15:12:42
Subject Test	0.196*	15:15:49
Air Blank	0.000	15:16:21
Air Blank	0.000	15:16:48
Subject Test	XXX**	15:17:31
Air Blank	0.000	15:18:04

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

5 of 7_{sc}
BH

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002201
Core Version: 8105.48
11/09/2012 15:18:16

Standard Lot#: 23511100A2
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	15:18:44
0.100 Cal Check	0.101	15:19:05
Air Blank	0.000	15:19:35

Preventative Maintenance
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

6 of 7 sc

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

TS4
11-19-2012

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 11/9/12 TIME 1500
Intoxilyzer Serial # 80-002201 Location Crime Lab

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 102 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

7 of 7 sc

104

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002201
Core Version: 8105.48
11/09/2012 15:00:21

Standard Lot#: 23511100A2
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	15:00:50
0.100 Cal Check	0.102	15:01:11
Air Blank	0.000	15:01:41

Preventative Maintenance
Performed -- Timer Reset