



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
1306 W. Miracle Mile, Tucson, AZ 85705

TSy  
9.24.2013

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002201 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 9/12/13 Cal 1 0.098 AC Cal 2 0.097 AC

Between the dates of 8/15/13 and 9/12/13 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.  
The instrument was returned to service on -.

Sarah Cione  
Sarah L. Cione

9/16/13  
ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 16<sup>th</sup>, day of September, 2013.

Edma Escarcega  
Notary Public

\*Includes documents specified in A.A.C. R13-10-104.A.5

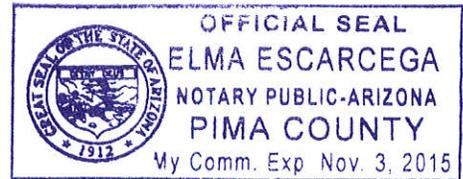


EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

1 of 7  
TS4  
9.24.2013

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 9/12/13 TIME 1328

INTOXILYZER SERIAL # 80-002201 LOCATION Crime Lab

(  ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (  ) 1. Clock time check.
- (  ) 2. Date check.

OPERATIONAL TESTS

- (  ) 1. Alcohol-free subject test result 0. 000 AC.
- (  ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- (  ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- (  ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 097 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

---



---



---



---

SIGNATURE Sarah Cione

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002201  
Core Version: 8105.48  
09/12/2013

13:30:36

2 of 7sc  
TS4

Standard Lot#: 30412100A1  
Last Changed By: R. SZELEWSKI #46927

QAS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1980  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	13:31:52
Diagnostic Test	Pass	13:32:27
Air Blank	0.000	13:32:54
0.100 Cal Check	0.098	13:33:15
Air Blank	0.000	13:33:45
Subject Test	0.000	13:34:09
Air Blank	0.000	13:34:37
Wait	ABT*	
Air Blank	0.000	13:35:12

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002201  
Core Version: 8105.48  
09/12/2013

13:35:19

3 of 7 sc  
T84

---

Standard Lot#: 30412100A1  
Last Changed By: R. SZELEWSKI #46927

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1980  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	13:36:03
Diagnostic Test	Pass	13:36:38
Air Blank	0.000	13:37:05
0.100 Cal Check	0.098	13:37:26
Air Blank	0.000	13:37:55
Subject Test	RFI*	13:38:19
Air Blank	IPS**	13:38:27

\*RFI Detect

\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002201  
Core Version: 8105.48  
09/12/2013

13:39:18

4 of 7<sub>sc</sub>  
TS4

---

Standard Lot#: 30412100A1  
Last Changed By: R. SZELEWSKI #46927

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1980  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	13:40:27
Diagnostic Test	Pass	13:41:02
Air Blank	0.000	13:41:30
0.100 Cal Check	0.098	13:41:50
Air Blank	0.000	13:42:20
Subject Test	0.179*	13:45:27
Air Blank	0.000	13:45:59
Air Blank	0.000	13:46:26
Subject Test	XXX**	13:46:46
Air Blank	0.000	13:47:20

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002201  
Core Version: 8105.48  
09/12/2013 13:47:45

5 of 7 sc  
TS4

-----  
Standard Lot#: 30412100A1  
Last Changed By: R. SZELEWSKI #46927  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Operator: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	13:48:14
0.100 Cal Check	0.097	13:48:35
Air Blank	0.000	13:49:04

Preventative Maintenance  
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6 of 7  
TSK  
9.24.2013

QA Specialist Sarah Cione AGENCY Tucson Police Department  
DATE 9/12/13 TIME 1328  
Intoxilyzer Serial # 80-002201 Location Crime Lab

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 098 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002201  
Core Version: 8105.48  
09/12/2013 13:28:35

7 of 7 SC  
TS

-----  
Standard Lot#: 30412100A1  
Last Changed By: R. SZELEWSKI #46927  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Operator: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	13:29:03
0.100 Cal Check	0.098	13:29:24
Air Blank	0.000	13:29:54

Preventative Maintenance  
Performed -- Timer Reset