



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

SC
3/20/14

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002201 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 3.11.2014 Cal 1 0.100 AC Cal 2 0.099 AC

Between the dates of N/A and N/A the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.

The instrument was returned to ^{TPD}service on 2.14.2014.

^{TS4}
On loan to UAPD from 1.22.2014 to 2.14.2014


Terry S. Gallegos

3.19.2014

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 19th, day of March, 2014.

Notary Public

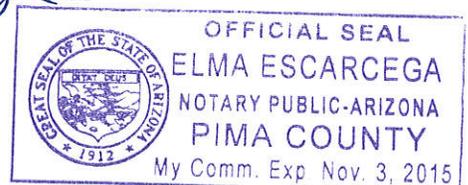

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

lof7
TS4
SC
3/20/14

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department
DATE 3-11-2014 TIME 1752
INTOXILYZER SERIAL # 80-002201 LOCATION Crime Lab

() 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- () 1. Clock time check.
- () 2. Date check.

OPERATIONAL TESTS

- () 1. Alcohol-free subject test result 0. 000 AC.
- () 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- () 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- () 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 099 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE 

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON PD U OF A
Serial Number: 80-002201
Core Version: 8105.48
03/11/2014 17:56:03

2 of 7
TS4
SC

Standard Lot#: 30412100A1
Last Changed By: R. SZELEWSKI #46927

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	17:56:55
Diagnostic Test	Pass	17:57:30
Air Blank	0.000	17:57:57
0.100 Cal Check	0.100	17:58:18
Air Blank	0.000	17:58:47
Subject Test	0.000	17:59:03
Air Blank	0.000	17:59:30
Wait	ABT*	
Air Blank	0.000	18:00:43

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD U OF A
Serial Number: 80-002201
Core Version: 8105.48
03/11/2014

18:00:54

30 of 7
TS4

SC

Standard Lot#: 30412100A1
Last Changed By: R. SZELEWSKI #46927

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	18:01:42
Diagnostic Test	Pass	18:02:17
Air Blank	0.000	18:02:44
0.100 Cal Check	0.100	18:03:05
Air Blank	RFI*	18:03:12
Air Blank	IPS**	18:03:20

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD U OF A
Serial Number: 80-002201
Core Version: 8105.48
03/11/2014 18:03:24

4 of 7
TS4

5C

Standard Lot#: 30412100A1
Last Changed By: R. SZELEWSKI #46927

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST,SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	18:04:18
Diagnostic Test	Pass	18:04:53
Air Blank	0.000	18:05:20
0.100 Cal Check	0.100	18:05:41
Air Blank	0.000	18:06:11
Subject Test	0.222*	18:09:18
Air Blank	0.000	18:09:51
Air Blank	0.000	18:10:19
Subject Test	XXX**	18:10:36
Air Blank	0.000	18:11:09

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD U OF A
Serial Number: 80-002201
Core Version: 8105.48
03/11/2014 18:11:25

507
TS4
SC

Standard Lot#: 30412100A1
Last Changed By: R. SZELEWSKI #46927

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	18:11:55
0.100 Cal Check	0.099	18:12:15
Air Blank	0.000	18:12:45

Preventative Maintenance
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

60f7
TBH

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

5C
3/20/14

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department
DATE 3.11.2014 TIME 1752
Intoxilyzer Serial # 80-002201 Location CrimeLab

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE 

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location: TUCSON PD U OF A
Serial Number: 80-002201
Core Version: 8105.48
03/11/2014 17:52:40

7 of 7
TS4
SC

Standard Lot#: 30412100A1
Last Changed By: R. SZELEWSKI #46927

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	17:53:31
0.100 Cal Check	0.100	17:53:52
Air Blank	0.000	17:54:22

Preventative Maintenance
Performed -- Timer Reset