



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

TS4
4.14.2014

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002201 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 3/28/14 Cal 1 0.100 AC Cal 2 0.100 AC

Between the dates of 3/11/14 and 3/28/14 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.
The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione


4/1/14

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 15th day of April, 2014.

Elma Escarcega

Notary Public

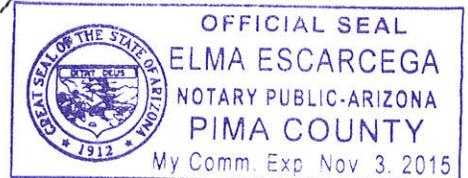

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

1 of 7 sc
TS4
4.14.2014

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 3/28/14 TIME 1514

INTOXILYZER SERIAL # 80-002201 LOCATION Crime Lab

() 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- () 1. Clock time check.
- () 2. Date check.

OPERATIONAL TESTS

- () 1. Alcohol-free subject test result 0. 000 AC.
- () 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- () 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- () 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 100 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Sarah Cione

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON PD U OF A
Serial Number: 80-002201
Core Version: 8105.48
03/28/2014

15:16:25

2 of 7_{sc}
TS4

Standard Lot#: 30412100A1
Last Changed By: R. SZELEWSKI #46927

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	15:17:23
Diagnostic Test	Pass	15:17:58
Air Blank	0.000	15:18:25
0.100 Cal Check	0.101	15:18:46
Air Blank	0.000	15:19:16
Subject Test	0.000	15:19:41
Air Blank	0.000	15:20:09
Wait	ABT*	
Air Blank	0.000	15:20:45

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD U OF A
Serial Number: 80-002201
Core Version: 8105.48
03/28/2014 15:21:08

3 of 7
sc
TSh

Standard Lot#: 30412100A1
Last Changed By: R. SZELEWSKI #46927

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	15:21:50
Diagnostic Test	Pass	15:22:25
Air Blank	0.000	15:22:53
0.100 Cal Check	0.100	15:23:13
Air Blank	RFI*	15:23:27
Air Blank	IPS**	15:23:35

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD U OF A
Serial Number: 80-002201
Core Version: 8105.48
03/28/2014

15:23:41

4 of 7^{sc}
Bh

Standard Lot#: 30412100A1
Last Changed By: R. SZELEWSKI #46927

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	15:24:29
Diagnostic Test	Pass	15:25:04
Air Blank	0.000	15:25:31
0.100 Cal Check	0.100	15:25:52
Air Blank	0.000	15:26:21
Subject Test	0.236*	15:29:29
Air Blank	0.000	15:30:01
Air Blank	0.000	15:30:28
Subject Test	XXX**	15:30:39
Air Blank	0.000	15:31:12

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD U OF A
Serial Number: 80-002201
Core Version: 8105.48
03/28/2014 15:31:33

5 of 7
sc
T84

Standard Lot#: 30412100A1
Last Changed By: R. SZELEWSKI #46927

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	15:32:02
0.100 Cal Check	0.100	15:32:22
Air Blank	0.000	15:32:52

Preventative Maintenance
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6 of 7 sc
TS4
4.14.2014

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 3/28/14 TIME 1514
Intoxilyzer Serial # 80-002201 Location Crime Lab

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location: TUCSON PD U OF A
Serial Number: 80-002201
Core Version: 8105.48
03/28/2014 15:14:23

7 of 7
sc
TS4

Standard Lot#: 30412100A1
Last Changed By: R. SZELEWSKI #46927

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	15:14:51
0.100 Cal Check	0.100	15:15:12
Air Blank	0.000	15:15:42

Preventative Maintenance
Performed -- Timer Reset