



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
1306 W. Miracle Mile, Tucson, AZ 85705

TJH  
10-29-2014

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002201 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 10/10/14 Cal 1 0.100 AC Cal 2 0.099 AC

Between the dates of 9/11/14 and 10/10/14 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.  
The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione


10/13/14

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 13<sup>th</sup>, day of October, 2014.

Elma Escarcega

Notary Public

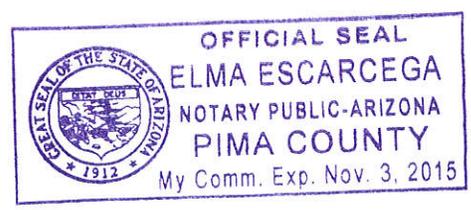

\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

1 of 7  
sc  
TJH  
10-29-2014

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 10/10/14 TIME 1427

INTOXILYZER SERIAL # 80-002201 LOCATION DUI HQ

( ✓ ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- ( ✓ ) 1. Clock time check.
- ( ✓ ) 2. Date check.

OPERATIONAL TESTS

- ( ✓ ) 1. Alcohol-free subject test result 0. 000 AC.
- ( ✓ ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- ( ✓ ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- ( ✓ ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 099 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE Sarah Cione

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002201  
Core Version: 8105.48  
10/10/2014

14:29:26

2 of 7<sub>sc</sub>

18h

---

Standard Lot#: 25913100A2  
Last Changed By: J. OLIVARES #36552

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1999  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	14:30:27
Diagnostic Test	Pass	14:31:02
Air Blank	0.000	14:31:29
0.100 Cal Check	0.100	14:31:50
Air Blank	0.000	14:32:20
Subject Test	0.000	14:32:50
Air Blank	0.000	14:33:17
Wait	ABT*	
Air Blank	0.000	14:33:53

\*Sequence Aborted

### Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002201  
Core Version: 8105.48  
10/10/2014

14:34:30

3 of 7  
50  
TS4

---

Standard Lot#: 25913100A2  
Last Changed By: J. OLIVARES #36552

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1999  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	14:35:43
Diagnostic Test	Pass	14:36:18
Air Blank	0.000	14:36:46
0.100 Cal Check	0.100	14:37:07
Air Blank	0.000	14:37:36
Subject Test	RFI*	14:38:30
Air Blank	IPS**	14:38:43

\*RFI Detect

\*\*Improper Sample

### Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002201  
Core Version: 8105.48  
10/10/2014

14:48:22

4 of 7  
sc  
184

---

Standard Lot#: 25913100A2  
Last Changed By: J. OLIVARES #36552

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1999  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	14:49:20
Diagnostic Test	Pass	14:49:55
Air Blank	0.000	14:50:22
0.100 Cal Check	0.100	14:50:43
Air Blank	0.000	14:51:12
Subject Test	0.196*	14:54:20
Air Blank	0.000	14:54:52
Air Blank	0.000	14:55:19
Subject Test	XXX**	14:55:45
Air Blank	0.000	14:56:18

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

### Not a Successfully Completed Test Sequence

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002201  
Core Version: 8105.48  
10/10/2014 14:56:33

5 of 7 sc  
Tdy

-----  
Standard Lot#: 25913100A2  
Last Changed By: J. OLIVARES #36552  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Operator: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	14:57:01
0.100 Cal Check	0.099	14:57:22
Air Blank	0.000	14:57:52

Preventative Maintenance  
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

607 7 SC  
TJ  
10-29-2014

QA Specialist Sarah Cione AGENCY Tucson Police Department  
DATE 10/10/14 TIME 1427  
Intoxilyzer Serial # 80-002201 Location DUI HQ

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 100 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002201  
Core Version: 8105.48  
10/10/2014 14:27:26

7 of 7  
50  
104

-----  
Standard Lot#: 25913100A2  
Last Changed By: J. OLIVARES #36552  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Operator: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	14:27:54
0.100 Cal Check	0.100	14:28:15
Air Blank	0.000	14:28:44

Preventative Maintenance  
Performed -- Timer Reset