



CERTIFICATION by CUSTODIAN of RECORDS  
 Tucson Police Department Crime Laboratory  
 1306 W. Miracle Mile, Tucson, AZ 85705

*AS*  
 10/23/15

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002201 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 10-8-2015 Cal 1 0.099 AC Cal 2 0.098 AC

Between the dates of 9-10-2015 and 10-8-2015 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on —.  
 The instrument was returned to service on —.

*Terry Gallegos*

Terry S. Gallegos


10-16-2015

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 16<sup>th</sup>, day of October, 2015.

*Elma Escarcega*

Notary Public

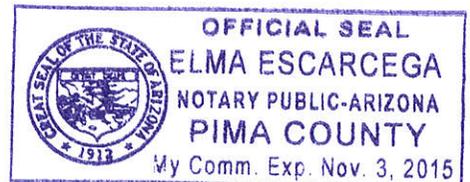

\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1 of 7  
TS4  
AB  
10/23/15

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department

DATE 10.8.2015 TIME 1206

INTOXILYZER SERIAL # 80-002201 LOCATION DUI-HQ

(  ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (  ) 1. Clock time check.
- (  ) 2. Date check.

OPERATIONAL TESTS

- (  ) 1. Alcohol-free subject test result 0. 000 AC.
- (  ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- (  ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- (  ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 098 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

---



---



---



---

SIGNATURE

Terry Gallegos

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002201  
Core Version: 8105.48  
10/08/2015

12:08:47

2 of 7  
TSy  
Jed

---

Standard Lot#: 25414100A4  
Last Changed By: T. GALLEGOS #28847

---

QAS: T. GALLEGOS #28847  
TUCSON PD

---

Operator: T. GALLEGOS #28847  
TUCSON PD

---

Subject: TEST, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	12:09:34
Diagnostic Test	Pass	12:10:09
Air Blank	0.000	12:10:36
0.100 Cal Check	0.099	12:10:57
Air Blank	0.000	12:11:27
Subject Test	0.000	12:11:56
Air Blank	0.000	12:12:23
Wait	ABT*	
Air Blank	0.000	12:13:07

\*Sequence Aborted

### Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002201  
Core Version: 8105.48  
10/08/2015

12:13:16

30 of 7  
TS4  
dd

Standard Lot#: 25414100A4  
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847  
TUCSON PD

Operator: T. GALLEGOS #28847  
TUCSON PD

Subject: TEST, SQAP

DOB: 11/11/1987

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	12:15:06
Diagnostic Test	Pass	12:15:41
Air Blank	0.000	12:16:09
0.100 Cal Check	RFI*	12:16:27
Air Blank,	IPS**	12:16:40

\*RFI Detect

\*\*Improper Sample

### Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002201  
Core Version: 8105.48  
10/08/2015

12:16:46

4 of 7  
By  
Me

Standard Lot#: 25414100A4  
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847  
TUCSON PD

Operator: T. GALLEGOS #28847  
TUCSON PD

Subject: TEST, SQAP

DOB: 11/11/1987

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	12:17:44
Diagnostic Test	Pass	12:18:19
Air Blank	0.000	12:18:47
0.100 Cal Check	0.098	12:19:07
Air Blank	0.000	12:19:37
Subject Test	0.114*	12:22:45
Air Blank	0.000	12:23:19
Air Blank	0.000	12:23:46
Subject Test	XXX**	12:24:22
Air Blank	0.000	12:24:56

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

### Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002201  
Core Version: 8105.48  
10/08/2015

12:25:21

*5 of 7*  
*TG*  
*add*

-----  
Standard Lot#: 25414100A4  
Last Changed By: T. GALLEGOS #28847  
-----

QAS: T. GALLEGOS #28847  
TUCSON PD  
-----

Operator: T. GALLEGOS #28847  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	12:25:50
0.100 Cal Check	0.098	12:26:11
Air Blank	0.000	12:26:41

Preventative Maintenance  
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6 of 7  
TS4  
dlt  
10/23/15

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department  
DATE 10-8-2015 TIME 1206  
Intoxilyzer Serial # 80 002201 Location DUI-HQ

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 099 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE 

DPS Form Exh G-2 (Rev. 05-01)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002201  
Core Version: 8105.48  
10/08/2015 12:06:13

7077  
TS4  
JLD

-----  
Standard Lot#: 25414100A4  
Last Changed By: T. GALLEGOS #28847  
-----

QAS: T. GALLEGOS #28847  
TUCSON PD  
-----

Operator: T. GALLEGOS #28847  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	12:06:42
0.100 Cal Check	0.099	12:07:03
Air Blank	0.000	12:07:33

Preventative Maintenance  
Performed -- Timer Reset