



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002261 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 12-28-2006 Cal 1 0.098 AC Cal 2 0.098 AC

Between the dates of 12-1-2006 and 12-28-2006 the instrument was operating properly and accurately.

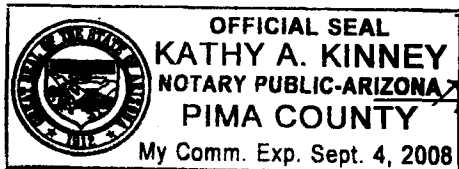
The instrument was taken out of service for repair on —.
The instrument was returned to service on —.

Terry S. Gallegos

1-4-2007

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 4, day of January, 2007.



Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104 (A)

1 of 9
TS4

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Terry Gallegos 20047 AGENCY Tucson Police DepartmentDATE 12-29-2006 TIME 1644INTOXILYZER SERIAL # 80-002261 LOCATION DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

(✓) 1. Clock time check.

(✓) 2. Date check.

OPERATIONAL TESTS

* (✓) 1. Alcohol-free subject test result 0. 000 AC.

(✓) 2. Error recognition logic system functioning.

Not a Successfully Completed Test Sequence printed ✓✓✓

(✓) 3. Proper sample recognition system.

Not a Successfully Completed Test Sequence printed

Deficient sample printed.

(✓) 4. Standard Calibration Check standard 0. 100 ACResult: 0. 098 ACInstrument is operating properly and accurately. YES ✓ NO

COMMENTS

Alcohol-free subject test performed twice due
to alcohol from a previous invalid sample check.

SIGNATURE

Terry Gallegos

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002261
12/28/2006

16:45:50

Standard Lot#: 502602I
Last Changed By: J. COOMBS #46160

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	16:46:57
Diagnostic Test	Pass	16:47:24
Air Blank	0.000	16:47:51
0.100 Cal Check	0.099	16:48:10
Air Blank	0.000	16:48:40
Subject Test	0.007	16:48:58
Air Blank	0.000	16:49:25
Five Minute Wait	ABT*	
Air Blank	0.000	16:50:03

*Sequence Aborted

Not a Successfully Completed Test Sequence

30f9
TS4

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002261
12/28/2006 16:50:43

Standard Lot#: 502602I
Last Changed By: J. COOMBS #46160

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	16:51:38
Diagnostic Test	Pass	16:52:06
Air Blank	0.000	16:52:33
0.100 Cal Check	0.099	16:52:52
Air Blank	0.000	16:53:21
Subject Test	0.000	16:53:43
Air Blank	RFI*	16:53:50
Air Blank	0.000	16:54:19

*RFI Detect

Not a Successfully Completed Test Sequence

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-002261

12/28/2006

16:54:31

40f9
TS4

Standard Lot#: 502602I

Last Changed By: J. COOMBS

#46160

QAS: T. GALLEGOS

#28847

TUCSON PD CRIME LAB

Operator: T. GALLEGOS

#28847

TUCSON PD CRIME LAB

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	IPS*	16:55:20
Air Blank	0.000	16:55:48

*Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-002261

12/28/2006

16:56:02

50f9
TS4

Standard Lot#: 502602I

Last Changed By: J. COOMBS

#46160

QAS: T. GALLEGOS

#28847

TUCSON PD CRIME LAB

Operator: T. GALLEGOS

#28847

TUCSON PD CRIME LAB

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	16:57:32
Diagnostic Test	Pass	16:57:59
Air Blank	0.000	16:58:26
0.100 Cal Check	0.098	16:58:45
Air Blank	0.000	16:59:15
Subject Test	0.043*	17:02:23
Air Blank	0.000	17:02:54
Air Blank	0.000	17:03:22
Subject Test	0.052	17:03:41
Air Blank	0.000	17:04:12
Five Minute Wait	ABT**	
Air Blank	0.000	17:04:56

*Deficient Sample

**Sequence Aborted

Invalid sample test accepted as 0.043^{52th} result.

Not a Successfully Completed Test Sequence

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-002261

12/28/2006

17:05:11

6089
TS4

Standard Lot#: 502602I

Last Changed By: J. COOMBS

#46160

QAS: T. GALLEGOS

#28847

TUCSON PD CRIME LAB

Operator: T. GALLEGOS

#28847

TUCSON PD CRIME LAB

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	17:06:01
Diagnostic Test	Pass	17:06:28
Air Blank	0.000	17:06:56
0.100 Cal Check	0.097	17:07:14
Air Blank	0.000	17:07:44
Subject Test	XXX*	17:10:36
Air Blank	0.000	17:11:12

*Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-002261

12/28/2006

17:11:30

Standard Lot#: 502602I

Last Changed By: J. COOMBS #46160

QAS: T. GALLEGOS #28847

TUCSON PD CRIME LAB

Test	g/210L	Time
Air Blank	0.000	17:11:57
Cal Check	0.098	17:12:16
Air Blank	0.000	17:12:45

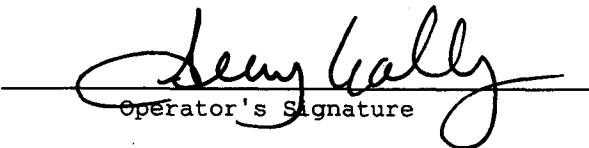

Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

8059
TS4

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Terry Callegos 28847 AGENCY Tucson Police Department
DATE 12-29-2006 TIME 1644
Intoxilyzer Serial # 80-002261 Location DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 098 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Terry Callegos

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-002261

12/28/2006

16:44:09

Standard Lot#: 502602I

Last Changed By: J. COOMBS #46160

QAS: T. GALLEGOS #28847

TUCSON PD CRIME LAB

Test	g/210L	Time
Air Blank	0.000	16:44:37
Cal Check	0.098	16:44:55
Air Blank	0.000	16:45:25


Operator's Signature