



CERTIFICATION by CUSTODIAN of RECORDS  
Tucson Police Department Crime Laboratory

77

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002261 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 3.26.2007 Cal 1 0.100 AC Cal 2 0.099 AC

Between the dates of 2.26.2007 and 3.26.2007 the instrument was operating properly and accurately.


The instrument was taken out of service for repair on \_\_\_\_\_.  
The instrument was returned to service on \_\_\_\_\_.

Terry S. Gallegos

April 6, 2007

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 6<sup>th</sup> day of April, 2007.

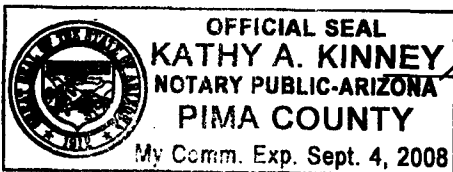

Notary Public

\*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

10F8  
TS4  
777

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department

DATE 3-26-2007 TIME 1453

INTOXILYZER SERIAL # 80-002261 LOCATION DUI HQ

( ✓ ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- ( ✓ ) 1. Clock time check.  
( ✓ ) 2. Date check.

OPERATIONAL TESTS

- ( ✓ ) 1. Alcohol-free subject test result 0. 000 AC.  
( ✓ ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed ✓✓✓  
( ✓ ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.  
( ✓ ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 099 AC

Instrument is operating properly and accurately. YES ✓ NO    

COMMENTS

\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE

Terry Gallegos

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-002261

03/26/2007

14:56:02

2 of 8  
Phm

---

Standard Lot#: 627101G

Last Changed By: T. GALLEGOS

#28847

---

QAS: T. GALLEGOS

#28847

TUCSON PD CRIME LAB

---

Operator: T. GALLEGOS

#28847

TUCSON PD CRIME LAB

---

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	14:56:59
Diagnostic Test	Pass	14:57:26
Air Blank	0.000	14:57:53
0.100 Cal Check	0.099	14:58:12
Air Blank	0.000	14:58:42
Subject Test	0.000	14:59:45
Air Blank	0.000	15:00:13
Five Minute Wait	ABT*	
Air Blank	0.000	15:00:54

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-002261  
03/26/2007

15:00:59

3058  
Bhmm

---

Standard Lot#: 627101G  
Last Changed By: T. GALLEGOS #28847

---

QAS: T. GALLEGOS #28847  
TUCSON PD CRIME LAB

---

Operator: T. GALLEGOS #28847  
TUCSON PD CRIME LAB

---

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	15:01:51
Diagnostic Test	Pass	15:02:19
Air Blank	RFI*	15:02:35
Air Blank	IPS**	15:02:45

\*RFI Detect

\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-002261  
03/26/2007

15:02:54

4088  
TShm

Standard Lot#: 627101G  
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847  
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847  
TUCSON PD CRIME LAB

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	15:03:45
Diagnostic Test	Pass	15:04:13
Air Blank	0.000	15:04:40
0.100 Cal Check	0.100	15:04:59
Air Blank	0.000	15:05:29
Subject Test	RGE*	15:07:18
Air Blank	0.000	15:07:55

\*Range Exceeded

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-002261  
03/26/2007

15:08:29

5088  
B4 m

Standard Lot#: 627101G  
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847  
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847  
TUCSON PD CRIME LAB

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	15:09:19
Diagnostic Test	Pass	15:09:46
Air Blank	0.000	15:10:14
0.100 Cal Check	0.099	15:10:32
Air Blank	0.000	15:11:02
Subject Test	0.048*	15:14:10
Air Blank	0.000	15:14:41
Air Blank	0.000	15:15:09
Subject Test	XXX**	15:15:39
Air Blank	0.000	15:16:13

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-002261

03/26/2007

15:16:42

60f8  
Dym

-----  
Standard Lot#: 627101G

Last Changed By: T. GALLEGOS #28847  
-----

QAS: T. GALLEGOS #28847

TUCSON PD CRIME LAB  
-----

Test	g/210L	Time
Air Blank	0.000	15:17:09
Cal Check	0.099	15:17:28
Air Blank	0.000	15:17:57

  
Operator's Signature

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY


STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

## STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Terry Gallegos 28047 AGENCY Tucson Police Department  
DATE 3-26-2007 TIME 1453  
Intoxilyzer Serial # 80-002261 Location DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach  $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

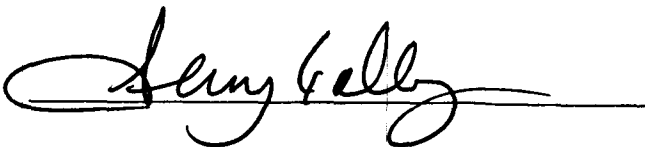

DPS Form Exh G-2 (Rev. 05-01)



INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-002261

03/26/2007

14:53:08

80f8  
Phm

-----  
Standard Lot#: 627101G

Last Changed By: T. GALLEGOS #28847  
-----

QAS: T. GALLEGOS #28847

TUCSON PD CRIME LAB  
-----

Test	g/210L	Time
Air Blank	0.000	14:53:36
Cal Check	0.100	14:53:54
Air Blank	0.000	14:54:24

  
Operator's Signature