



CERTIFICATION by CUSTODIAN of RECORDS  
Tucson Police Department Crime Laboratory <sup>WMP  
6-5-07</sup>

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002261 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 5-24-2007 Cal 1 0.101 AC Cal 2 0.100 AC

Between the dates of 4-26-2007 and 5-24-2007 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on                     .  
The instrument was returned to service on                     .

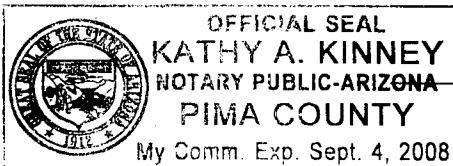
Terry S. Gallegos

Terry S. Gallegos

6-5-2007

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 5, day of June, 2007.



Kathy A. Kinney

Notary Public

\*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

10 of 10  
Tdy 6:50 PMQA SPECIALIST Terry Gallegos 2007 AGENCY Tucson Police DepartmentDATE 5.24.2007 TIME 1824INTOXILYZER SERIAL # 80-002261 LOCATION DUI HQ( ☒ ) 1. Display reads "PUSH BUTTON TO START"

## DIAGNOSTIC TESTS

- ( ☐ ) 1. Clock time check.  
( ☒ ) 2. Date check.

## OPERATIONAL TESTS

- ( ☒ ) 1. Alcohol-free subject test result 0. 000 AC.  
( ☒ ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed ✓✓✓  
( ☒ ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.  
( ☒ ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 100 AC

Instrument is operating properly and accurately. YES ☒ NO ☐

## COMMENTS

---

---

---

---

SIGNATURE

Terry Gallegos

DPS Form Exh G-4 (Rev. 05-01)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-002261

05/24/2007

18:26:20

20 of 10  
TS4<sup>mm</sup>  
6-5-11

Standard Lot#: 627101G

Last Changed By: T. GALLEGOS

#28847

QAS: T. GALLEGOS

#28847

TUCSON PD CRIME LAB

Operator: T. GALLEGOS

#28847

TUCSON PD CRIME LAB

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	18:27:10
Diagnostic Test	Pass	18:27:37
Air Blank	0.000	18:28:04
0.100 Cal Check	0.101	18:28:23
Air Blank	0.000	18:28:53
Subject Test	0.000	18:30:08
Air Blank	0.000	18:30:35
Five Minute Wait	ABT*	
Air Blank	0.000	18:31:10

\*Sequence Aborted

### Not a Successfully Completed Test Sequence

Pursuant to section 28-1527, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-002261

05/24/2007

18:31:54

3 of 10  
Tdy 7-5-07

Standard Lot#: 627101G

Last Changed By: T. GALLEGOS

#28847

QAS: T. GALLEGOS

#28847

TUCSON PD CRIME LAB

Operator: T. GALLEGOS

#28847

TUCSON PD CRIME LAB

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	18:32:33
Air Blank	IPS**	18:32:49

\*RFI Detect

\*\*Improper Sample

### Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-002261

05/24/2007

18:32:57

4 of 10  
TJ 6-5-07

Standard Lot#: 627101G

Last Changed By: T. GALLEGOS

#28847

QAS: T. GALLEGOS

#28847

TUCSON PD CRIME LAB

Operator: T. GALLEGOS

#28847

TUCSON PD CRIME LAB

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	AMB*	18:33:48
Air Blank	0.000	18:34:16

\*Ambient Fail

### Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-002261

05/24/2007

18:36:28

5 of 10  
TS4 6-5-07

Standard Lot#: 627101G

Last Changed By: T. GALLEGOS

#28847

QAS: T. GALLEGOS

#28847

TUCSON PD CRIME LAB

Operator: T. GALLEGOS

#28847

TUCSON PD CRIME LAB

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	18:37:21
Diagnostic Test	Pass	18:37:48
Air Blank	0.000	18:38:16
0.100 Cal Check	0.101	18:38:34
Air Blank	0.000	18:39:04
Subject Test	0.069*	18:42:12
Air Blank	0.000	18:42:45
Air Blank	0.000	18:43:12
Subject Test	XXX**	18:43:38
Air Blank	0.000	18:44:13

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

### Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-002261

05/24/2007

18:45:04

60510  
TS6<sup>mm</sup>  
6507

-----  
Standard Lot#: 627101G

Last Changed By: T. GALLEGOS #28847  
-----

QAS: T. GALLEGOS #28847

TUCSON PD CRIME LAB  
-----

Test	g/210L	Time
Air Blank	0.000	18:45:31
Cal Check	0.100	18:45:50
Air Blank	0.000	18:46:19

  
Operator's Signature

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 800040810  
TSG  
V-507

## STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department  
DATE 5-24-2007 TIME 1824  
Intoxilyzer Serial # 80-002261 Location DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach  $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 101 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Terry Gallegos

DPS Form Exh G-2 (Rev. 05-01)



INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-00226

05/24/2007

18:24:52

80 of 10  
TJG  
6-50-07

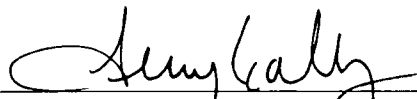
-----  
Standard Lot#: 627101G

Last Changed By: T. GALLEGOS #28847  
-----

QAS: T. GALLEGOS #28847

TUCSON PD CRIME LAB  
-----

Test	g/210L	Time
Air Blank	0.000	18:25:19
Cal Check	0.101	18:25:38
Air Blank	0.000	18:26:07

  
-----  
Operator's Signature



TUCSON POLICE DEPARTMENT  
CRIME LABORATORY

INTOXILYZER SERVICE RECORD  
FOR  
CALIBRATION GAS TANK CHANGE

Date: 5.24.2007 9 of 10 <sup>10:45</sup>  
PSY  
QAS: Terry Gallegos 28847  
Intoxilyzer Serial Number: 80-002261  
New Lot Number: 630601 I  
New Expiration Date: 11-6-2008

- (✓) 1. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 2. Remove the old gas tank and install the new one.
- (✓) 3. Hit the "ESC" key twice and swipe the QAS card.
- (✓) 4. Go to the "M-Maintenance" menu and select "S-Standard Change."  
Select "D" for dry gas and enter the following information:  
Enter Target Value: 0.100  
Enter Control Lot #: From the tank label  
Enter Exp Date: From the tank label (mm/dd/yyyy)  
Edit info Y/N? Y or N  
Update Std Info? Y ENTER  
"Please wait saving..."
- (✓) 5. Hit the "ESC" key twice.
- (✓) 6. Go to the "M-Maintenance" menu and select "T-Control Testing."  
Select "D" for dry control test. After selection is made press ENTER.
- (✓) 7. Air blank completed.
- (✓) 8. Calibration check completed. Test results 0. 101 AC.
- (✓) 9. Air blank completed.
- (✓) 10. Remove and sign printed record. Attach the record to the completed checklist.

SIGNATURE

Terry Gallegos

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-002261

05/24/2007

18:50:04

100f 10  
Ph<sup>mm</sup>  
6/5/07

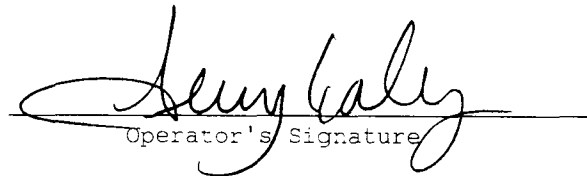
-----  
Standard Lot#: 630601I

Last Changed By: T. GALLEGOS #28847  
-----

QAS: T. GALLEGOS #28847

TUCSON PD CRIME LAB  
-----

Test	g/210L	Time
Air Blank	0.000	18:50:31
Cal Check	0.101	18:50:50
Air Blank	0.000	18:51:20

  
Operator's Signature