



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002261 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

TSY
3-17-08
BY

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 3/4/08 Cal 1 0.101 AC Cal 2 0.100 AC

Between the dates of 2/14/08 and 3/4/08 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on 3/4/08.
The instrument was returned to service on .


Sarah Cione

Sarah L. Cione

3/14/08

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 14, day of March, 2008.

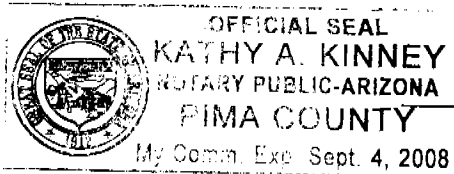

Kathy A Kinney

Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

1 of 75c

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

by
3.17.08

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 3/4/08 TIME 13:58
Intoxilyzer Serial # 80-002261 Location Crime Lab

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 101 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-002261

03/04/2008

13:58:39

Standard Lot#: 730901I

Last Changed By: C. NOLAN #38874

QAS: S. CIONE #50616

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	13:59:06
Cal Check	0.101	13:59:25
Air Blank	0.000	13:59:55



Operator's Signature

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

3 of 75c

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police DepartmentDATE 3/4/08 TIME 13:58INTOXILYZER SERIAL # 80-002261 LOCATION Crime Lab

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

(✓) 1. Clock time check.
(✓) 2. Date check.

OPERATIONAL TESTS

(✓) 1. Alcohol-free subject test result 0. 000 AC.
(✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
(✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
(✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 100 ACInstrument is operating properly and accurately. YES ✓ NO

COMMENTS

SIGNATURE

Sarah Cione

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002261
03/04/2008

14:00:54

4 of 75c
134

Standard Lot#: 730901I
Last Changed By: C. NOLAN #38874

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP,TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:01:52
Diagnostic Test	Pass	14:02:20
Air Blank	0.000	14:02:47
0.100 Cal Check	0.101	14:03:06
Air Blank	0.000	14:03:36
Subject Test	0.000	14:03:57
Air Blank	0.000	14:04:24
Five Minute Wait	ABT*	
Air Blank	0.000	14:05:03

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002261
03/04/2008

14:05:10

5 of 7sc
64

Standard Lot#: 730901I
Last Changed By: C. NOLAN #38874

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP,TEST

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:06:03
Diagnostic Test	Pass	14:06:30
Air Blank	RFI*	14:06:43
Air Blank	IPS**	14:06:54

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002261
03/04/2008

14:07:03

6 of 75c
BH

Standard Lot#: 730901I
Last Changed By: C. NOLAN #38874

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP,TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:07:49
Diagnostic Test	Pass	14:08:17
Air Blank	0.000	14:08:44
0.100 Cal Check	0.101	14:09:03
Air Blank	0.000	14:09:33
Subject Test	0.192*	14:12:41
Air Blank	0.000	14:13:14
Air Blank	0.000	14:13:41
Subject Test	XXX**	14:14:04
Air Blank	0.000	14:14:39

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-002261

03/04/2008

14:14:56


Standard Lot#: 730901I

Last Changed By: C. NOLAN #38874

QAS: S. CIONE #50616

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	14:15:23
Cal Check	0.100	14:15:42
Air Blank	0.000	14:16:12

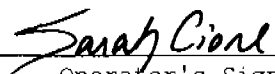

Operator's Signature



TUCSON POLICE DEPARTMENT
CRIME LABORATORY

INTOXILYZER SERVICE RECORD

1 of 25

TS4
3-13-08

Date: 3/4/08

QAS: Sarah Cione

Intoxilyzer Serial Number: 80-002261

Reason For Service: Could not complete data download on 2/14/08.
On 3/4/08 attempted to download w/ assistance from CMI, + still
could not complete data download.

Date Out of Service: 3/3/08

Calibration Check completed prior to removal from service

Date: 3/4/08

SQAP completed prior to removal from service

Date: 3/4/08

Corrective Action:

☒ Sent to CMI, Inc. (3/5/08)



Return Material Authorization Form

This form **MUST** be completed and enclosed with item(s) being returned.
Failure to complete and return this form may cause delays in crediting / repairing.

2 of 25c

BH

1. Item(s) Returned:



Instrument



Supplies



Other

(Note: please ship items in their original shipping container(s) or a similar protective box.)

2. Instrument Model:

Intox 8000

Serial Number 80-002261

3. Bill To Address:

Finance / Accounts Payable

P.O. Box 27450

TUCSON, AZ 85726-7450

Ship To Address: (NO PO BOXES)

Tucson Police Department

Crime Lab

270 S. Stone Ave

TUCSON, AZ 85701-1917

Customer # 105075

Name _____

Phone: (520) -791-4494

Fax: (520) 791-4642

Email: _____

4. Reason for Return: Note - If you are sending an item for repair, please give a detailed description of the problem. Please list any special instructions that you may have concerning this return. Note - Returns for credit may be subject to a 15% restocking charge.

Attn: Brian Faulkner - Please complete data download

5. Purchase Order Number (attach a copy of P.O. if applicable)

20669

6. I Authorize All Repairs:



Yes



No

7. I Authorize Repairs Up To:

\$ _____

8. I Require An Estimate Regardless of Cost



Yes



No

Note: An evaluation fee will apply to estimates that are not repaired.

9. Authorized By:

Terry Gallegos

Name (Please Print)

Terry Gallegos

Signature

3-4-2008

Date