



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

9-1-09

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002261 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 8-20-2009 Cal 1 0.099 AC Cal 2 0.098 AC

Between the dates of 7-23-2009 and 8-20-2009 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on _____.
The instrument was returned to service on _____.

Terry S. Gallegos

Terry S. Gallegos

9.1.2009


ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 15, day of September, 2009.

A. M. Hernández

Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

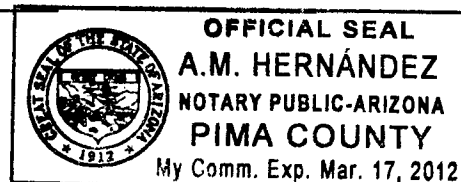

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

CU
9-1-09
107 7
TS4

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department

DATE 8-20-2009 TIME 1209

INTOXILYZER SERIAL # 80-002261 LOCATION DUI HQ

(☒) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (☒) 1. Clock time check.
(☒) 2. Date check.

OPERATIONAL TESTS

- (☒) 1. Alcohol-free subject test result 0. 000 AC.
(☒) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed ✓✓✓
(☒) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed ✓
Deficient sample printed.
(☒) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 098 AC

Instrument is operating properly and accurately. YES ☒ NO ☐

COMMENTS

SIGNATURE

Terry Gallegos

DPS Form Exh G-4 (Rev. 05-01)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department

INTOXILYZER 8000
Location: TUCSON PD
Serial Number : 80-002261
08/20/2009

12:13:11

CM
20f 7
TS4

Standard Lot#: 826001I
Last Changed By: C. NOLAN #38874

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	12:14:24
Diagnostic Test	Pass	12:14:52
Air Blank	0.000	12:15:19
0.100 Cal Check	0.099	12:15:38
Air Blank	0.000	12:16:08
Subject Test	0.000	12:16:25
Air Blank	0.000	12:16:52
Five Minute Wait	ABT*	
Air Blank	0.000	12:17:41

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number : 80-002261
08/20/2009

12:17:46

Standard Lot#: 826001I
Last Changed By: C. NOLAN #38874

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	12:18:39
Diagnostic Test	Pass	12:19:06
Air Blank	0.000	12:19:33
0.100 Cal Check	0.099	12:19:52
Air Blank	RFI*	12:20:07
Air Blank	IPS**	12:20:24

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number : 80-002261
08/20/2009

12:22:05

Standard Lot#: 826001I
Last Changed By: C. NOLAN #38874

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	12:22:55
Diagnostic Test	Pass	12:23:22
Air Blank	0.000	12:23:50
0.100 Cal Check	0.100	12:24:08
Air Blank	0.000	12:24:38
Subject Test	0.234*	12:27:46
Air Blank	0.000	12:28:18
Air Blank	0.000	12:28:45
Subject Test	XXX**	12:29:12
Air Blank	0.000	12:29:45

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000

Location: TUCSON PD

Serial Number: 80-002261

08/20/2009

12:30:37

Standard Lot#: 826001I

Last Changed By: C. NOLAN #38874

QAS: T. GALLEGOS #28847

TUCSON PD CRIME LAB

Test	g/210L	Time
Air Blank	0.000	12:31:04
Cal Check	0.098	12:31:23
Air Blank	0.000	12:31:53


Operator's Signature

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Terry Gallegos 2801A AGENCY Tucson Police Department
DATE 8.20.2009 TIME 1209
Intoxilyzer Serial # 80-002261 Location DUI HQ

- (☒) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (☒) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (☒) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (☒) 4. Air blank completed.
- (☒) 5. Calibration check completed. Test results 0. 099 AC.
- (☒) 6. Air blank completed.
- (☒) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Terry Gallegos

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000

Location: TUCSON PD

Serial Number: 80-002261

08/20/2009

12:09:48


Standard Lot#: 826001I

Last Changed By: C. NOLAN #38874

QAS: T. GALLEGOS #28847

TUCSON PD CRIME LAB

Test	g/210L	Time
Air Blank	0.000	12:10:15
Cal Check	0.099	12:10:34
Air Blank	0.000	12:11:03

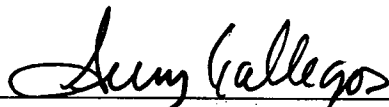

Operator's Signature



TUCSON POLICE DEPARTMENT
CRIME LABORATORY

INTOXILYZER SERVICE RECORD
FOR
CALIBRATION GAS TANK CHANGE

Date: 8/6/09
QAS: 110612
Intoxilyzer Serial Number: 80-002261
New Lot Number: 826001I
New Expiration Date: 24 Sep 2010

- (✓) 1. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 2. Remove the old gas tank and install the new one.
- (✓) 3. Hit the "ESC" key twice and swipe the QAS card.
- (✓) 4. Go to the "M-Maintenance" menu and select "S-Standard Change."
Select "D" for dry gas and enter the following information:
Enter Target Value: 0.100
Enter Control Lot #: From the tank label
Enter Exp Date: From the tank label (mm/dd/yyyy)
Edit info Y/N? Y or N
Update Std Info? Y ENTER
"Please wait saving..."
- (✓) 5. Hit the "ESC" key twice.
- (✓) 6. Go to the "M-Maintenance" menu and select "T-Control Testing."
Select "D" for dry control test. After selection is made press
ENTER.
- (✓) 7. Air blank completed.
- (✓) 8. Calibration check completed. Test results 0. 100 AC.
- (✓) 9. Air blank completed.
- (✓) 10. Remove and sign printed record. Attach the record to the completed checklist.

SIGNATURE

[Signature]

2 of 2
TSY
CM

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002261
08/06/2009 21:34:25

Standard Lot#: 8260011
Last Changed By: C. NOLAN #38874

OAS: C. NOLAN #38874
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	21:34:53
Cal Check	0.100	21:35:11
Air Blank	0.000	21:35:41


Operator's Signature