



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002261 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

3-10-2010

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 3-4-10 Cal 1 0.099 AC Cal 2 0.098 AC

Between the dates of 2-4-10 and 3-4-10 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.
The instrument was returned to service on N/A.

Clare Conley

Clare L. Conley


3-4-10

ADPS QAS Permit # 33314

SUBSCRIBED AND SWORN to before me this 4th day of March, 2010.

Susan A. Hitchye

Notary Public

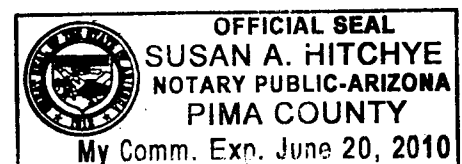

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1587
TS4
3.10.2010
Ch

QA SPECIALIST Clare Conley AGENCY Tucson Police Department

DATE 3-4-10 TIME 0902

INTOXILYZER SERIAL # 80-002261 LOCATION DDW DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

(✓) 1. Clock time check.

(✓) 2. Date check.

OPERATIONAL TESTS

(✓) 1. Alcohol-free subject test result 0. 000 AC.

(✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed

(✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.

(✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 078 AC

Instrument is operating properly and accurately. YES ☒ NO ☐

COMMENTS None.

SIGNATURE Clare Conley

INTOXILYZER 8000

Location: TUCSON PD

Serial Number : 80-002261

03/04/2010

09:03:37

Standard Lot#: 826001I

Last Changed By: C. NOLAN

#38874

QAS: C. CONLEY

#34517

TUCSON PD

Operator: C. CONLEY

#34517

TUCSON PD

Subject: TEST, SQAP

DOB: 11/11/1111

Sex: M

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:05:02
Diagnostic Test	Pass	09:05:29
Air Blank	0.000	09:05:56
0.100 Cal Check	0.100	09:06:15
Air Blank	0.000	09:06:44
Subject Test	0.000	09:07:14
Air Blank	ABT*	09:07:21
Air Blank	0.000	09:07:49

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number : 80-002261
03/04/2010

09:07:54

Standard Lot#: 826001I
Last Changed By: C. NOLAN #38874

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	α/210L	Time
Air Blank	RFI*	09:08:25
Air Blank	IPS**	09:08:32

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number : 80-002261
03/04/2010

09:09:42

457
CW

Standard Lot#: 826001I
Last Changed By: C. NOLAN #38874

BH

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP

DOB: 11/11/1111

Sex: M

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:10:31
Diagnostic Test	Pass	09:10:58
Air Blank	0.000	09:11:25
0.100 Cal Check	0.100	09:11:44
Air Blank	0.000	09:12:13
Subject Test	0.451*	09:15:22
Air Blank	0.000	09:15:54
Air Blank	0.000	09:16:21
Subject Test	XXX**	09:17:14
Air Blank	0.000	09:17:48

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000

Location: TUCSON PD

Serial Number: 80-002261

03/04/2010

09:18:10

Standard Lot#: 826001I

Last Changed By: C. NOLAN #38874

QAS: C. CONLEY #34517

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:18:37
Cal Check	0.098	09:18:56
Air Blank	0.000	09:19:25

Clare Conley

Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

687
154
3.10.2010

QA Specialist Clare Conley AGENCY Tucson Police Department
DATE 3-4-10 TIME 0902
Intoxilyzer Serial # 80-002261 Location ODW Du/40

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0.100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0.099 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Clare Conley

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000

Location: TUCSON PD

Serial Number: 80-002261

03/04/2010

09:02:14

Standard Lot#: 826001I

Last Changed By: C. NOLAN #38874

QAS: C. CONLEY #34517

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:02:41
Cal Check	0.099	09:03:00
Air Blank	0.000	09:03:29

Clare Conley

Operator's Signature

767
lu

154