



# CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
270 S. Stone Avenue Tucson, AZ 85701-1917

TS4  
7-22-2010

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002261 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 7-22-10 Cal 1 0.098 AC Cal 2 0.098 AC

Between the dates of 6-24-10 and 7-22-10  
the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.  
The instrument was returned to service on N/A.

Clare Conley

Clare L. Conley

7-22-10

ADPS QAS Permit # 33314


SUBSCRIBED AND SWORN to before me this 22, day of

July, 20 10.

AM Hernandez

Notary Public

\*Includes documents specified in A.A.C. R13-10-104.A.5

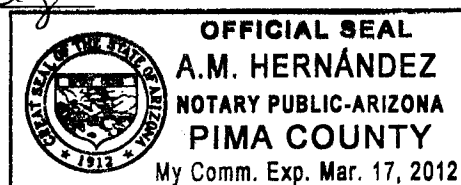

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

10/17  
eq  
BS  
7-22-2010

QA SPECIALIST Clare Conley AGENCY Tucson Police Department

DATE 7-22-10 TIME 1321

INTOXILYZER SERIAL # 80-002261 LOCATION Crime Lab (TPD)

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.  
(✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.  
(✓) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed  
(✓) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.  
(✓) 4. Standard Calibration Check standard 0. 107 AC  
Result: 0. 098 AC

Instrument is operating properly and accurately. YES ✓ NO     

COMMENTS None.  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE Clare Conley

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number : 80-002261  
07/22/2010

13:23:15

267  
cu

Ph

Standard Lot#: 826001I  
Last Changed By: C. NOLAN #38874

QAS: C. CONLEY #34517  
TUCSON PD

Operator: C. CONLEY #34517  
TUCSON PD

Subject: TEST, SQAP  
DOB: 11/11/1111  
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	13:23:59
Diagnostic Test	Pass	13:24:26
Air Blank	0.000	13:24:54
0.100 Cal Check	0.098	13:25:12
Air Blank	0.000	13:25:42
Subject Test	0.000	13:26:18
Air Blank	ABT*	13:26:24
Air Blank	0.000	13:26:52

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number : 80-002261  
07/22/2010

13:27:13

Standard Lot#: 826001I  
Last Changed By: C. NOLAN #38874

QAS: C. CONLEY #34517  
TUCSON PD

Operator: C. CONLEY #34517  
TUCSON PD

Subject: TEST, SQAP  
DOB: 11/11/1111  
Sex: M

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	13:28:11
Diagnostic Test	Pass	13:28:38
Air Blank	0.000	13:29:06
0.100 Cal Check	RFI*	13:29:20
Air Blank	IPS**	13:29:32

\*RFI Detect

\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number : 80-002261  
07/22/2010

13:29:50

4567  
cm

TS4

Standard Lot#: 826001I  
Last Changed By: C. NOLAN #38874

QAS: C. CONLEY #34517  
TUCSON PD

Operator: C. CONLEY #34517  
TUCSON PD

Subject: TEST, SQAP

DOB: 11/11/1111

Sex: M

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	13:30:39
Diagnostic Test	Pass	13:31:07
Air Blank	0.000	13:31:34
0.100 Cal Check	0.099	13:31:53
Air Blank	0.000	13:32:22
Subject Test	0.159*	13:35:31
Air Blank	0.000	13:36:02
Air Blank	0.000	13:36:29
Subject Test	XXX**	13:36:53
Air Blank	0.000	13:37:26

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000

Location: TUCSON PD

Serial Number: 80-002261

07/22/2010

13:37:46

-----  
Standard Lot#: 826001I

Last Changed By: C. NOLAN #38874  
-----

QAS: C. CONLEY #34517

TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	13:38:13
Cal Check	0.098	13:38:32
Air Blank	0.000	13:39:02

Clare Conley  
Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

b67  
on

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

PD  
7-22-2010

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Clare Conley AGENCY Tucson Police Department  
DATE 7-22-10 TIME 1321  
Intoxilyzer Serial # 80-002261 Location TPD Crime Lab

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach  $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 098 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Clare Conley

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000

Location: TUCSON PD

Serial Number: 80-002261

07/22/2010

13:21:47

-----  
Standard Lot#: 826001I

Last Changed By: C. NOLAN #38874  
-----

QAS: C. CONLEY #34517

TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	13:22:14
Cal Check	0.098	13:22:33
Air Blank	0.000	13:23:02

*Clare Conley*

\_\_\_\_\_  
Operator's Signature