



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

TS4
10-18-2010

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002261 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 10/14/10 Cal 1 0.099 AC Cal 2 0.099 AC

Between the dates of 9/16/10 and 10/14/10 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -
The instrument was returned to service on -

Sarah Cione

Sarah L. Cione


10/14/10

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 14, day of October, 2010.

Ann Hernandez

Notary Public

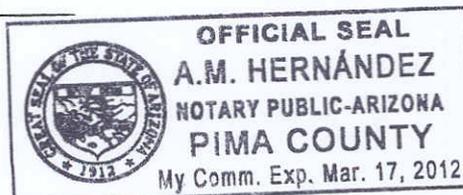

*Includes documents specified in A.A.C. R13-10-104.A.5

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

1 of 8 sc

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

TSH
10-18-2010

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 10/14/10 TIME 942

INTOXILYZER SERIAL # 80-002261 LOCATION DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 099 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Sarah Cione

2 of 8_{sc}

104

INTOXILYZER 8000
Location: TUCSON PD
Serial Number : 80-002261
10/14/2010 09:43:38

Standard Lot#: 719919
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:44:22
Diagnostic Test	Pass	09:44:49
Air Blank	0.000	09:45:16
0.100 Cal Check	0.099	09:45:35
Air Blank	0.000	09:46:04
Subject Test	0.000	09:46:35
Air Blank	0.000	09:47:02
Five Minute Wait	ABT*	
Air Blank	0.000	09:47:41

*Sequence Aborted

Not a Successfully Completed Test Sequence

3 of 8 SC

134

INTOXILYZER 8000
Location: TUCSON PD
Serial Number : 80-002261
10/14/2010

09:47:54

Standard Lot#: 719919
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:48:36
Diagnostic Test	Pass	09:49:03
Air Blank	0.000	09:49:31
0.100 Cal Check	0.099	09:49:49
Air Blank	0.000	09:50:19
Subject Test	RFI*	09:50:31
Air Blank	IPS**	09:50:40

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

4 of 8_{sc}

TS4

INTOXILYZER 8000
Location: TUCSON PD
Serial Number : 80-002261
10/14/2010 09:52:51

Standard Lot#: 719919
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:53:46
Diagnostic Test	Fail*	09:53:50
Air Blank	0.000	09:54:18

* Voltage/Current Test

Not a Successfully Completed Test Sequence

5 of 8 sc

INTOXILYZER 8000
Location: TUCSON PD
Serial Number : 80-002261
10/14/2010

09:55:01

134

Standard Lot#: 719919
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:55:53
Diagnostic Test	Pass	09:56:21
Air Blank	0.000	09:56:48
0.100 Cal Check	0.099	09:57:07
Air Blank	0.000	09:57:37
Subject Test	0.025*	10:00:45
Air Blank	0.000	10:01:15
Air Blank	0.000	10:01:42
Subject Test	XXX**	10:02:01
Air Blank	0.000	10:02:33

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

6 of 8 sc


INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002261
10/14/2010 10:03:05

BH

Standard Lot#: 719919
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	10:03:32
Cal Check	0.099	10:03:51
Air Blank	0.000	10:04:20

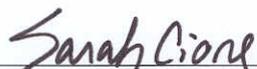

Operator's Signature

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

TS4
10.18.2010

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 10/14/10 TIME 942
Intoxilyzer Serial # 80-002261 Location DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 099 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

8 of 8 sc

154

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002261
10/14/2010 09:42:06

Standard Lot#: 719919
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:42:33
Cal Check	0.099	09:42:52
Air Blank	0.000	09:43:21

Sarah Cione
Operator's Signature