



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
270 S. Stone Avenue Tucson, AZ 85701-1917

TSy  
3.21.2011

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002261 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 3-3-11 Cal 1 0.099 AC Cal 2 0.099 AC

Between the dates of 2-3-11 and 3-3-11 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.  
The instrument was returned to service on N/A.

Clare Conley

Clare L. Conley


3-16-11

ADPS QAS Permit # 33314

SUBSCRIBED AND SWORN to before me this 16th, day of March, 20 11.

Isabel C. Trujillo

Notary Public

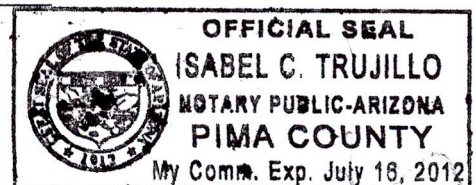

\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1067  
CM

TSY  
3-21-2011

QA SPECIALIST Clare Conley AGENCY Tucson Police Department

DATE 3-3-11 TIME 0749

INTOXILYZER SERIAL # 80-002261 LOCATION ODW DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.  
(✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.  
(✓) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed  
(✓) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.  
(✓) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 099 AC

Instrument is operating properly and accurately. YES ✓ NO       

COMMENTS None.  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE Clare Conley

2867  
cu

TS4

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002261  
03/03/2011 07:50:44

-----  
Standard Lot#: 719919  
Last Changed By: S. CIONE #50616  
-----

OAS: C. CONLEY #34517  
TUCSON PD  
-----

Operator: C. CONLEY #34517  
TUCSON PD  
-----

Subject: TEST,SOAP  
DOB: 11/11/1111  
Sex: M Weight: 111  
-----

15 Minute Deprivation Period? Yes

-----  
Test g/210L Time  
-----  
Air Blank 0.000 07:51:29  
Diagnostic Test Pass 07:51:58  
Air Blank 0.000 07:52:25  
0.100 Cal Check 0.100 07:52:44  
Air Blank 0.000 07:53:13  
Subject Test 0.000 07:53:32  
Air Blank ABT\* 07:53:41  
Air Blank 0.000 07:54:09

\*Sequence Aborted

Not a Successfully  
Completed  
Test Sequence

3567  
CW

TS4

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002261  
03/03/2011 07:54:24

Standard Lot#: 719919  
Last Changed By: S. CIONE #50616

QAS: C. CONLEY #34517  
TUCSON PD

Operator: C. CONLEY #34517  
TUCSON PD

Subject: TEST, SQAP  
DOB: 11/11/1111  
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

| Test      | g/210L | Time     |
|-----------|--------|----------|
| Air Blank | RFI*   | 07:55:03 |
| Air Blank | IPS**  | 07:55:10 |

\*RFI Detect  
\*\*Improper Sample

Not a Successfully  
Completed  
Test Sequence

4567  
Cu

734

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002261  
03/03/2011 07:56:07

Standard Lot#: 719919  
Last Changed By: S. CIONE #50616

QAS: C. CONLEY #34517  
TUCSON PD

Operator: C. CONLEY #34517  
TUCSON PD

Subject: TEST, SQAP  
DOB: 11/11/1111  
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

| Test            | g/210L | Time     |
|-----------------|--------|----------|
| Air Blank       | 0.000  | 07:56:58 |
| Diagnostic Test | Pass   | 07:57:27 |
| Air Blank       | 0.000  | 07:57:54 |
| 0.100 Cal Check | 0.100  | 07:58:13 |
| Air Blank       | 0.000  | 07:58:42 |
| Subject Test    | 0.088* | 08:01:51 |
| Air Blank       | 0.000  | 08:02:22 |
| Air Blank       | 0.000  | 08:02:50 |
| Subject Test    | XXX**  | 08:03:31 |
| Air Blank       | 0.000  | 08:04:04 |

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

Not a Successfully  
Completed  
Test Sequence

5067  
Cu  
154

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002261  
03/03/2011 08:05:03

-----  
Standard Lot#: 719919  
Last Changed By: S. CLONE #50616  
-----

QAS: C. CONLEY #34517  
TUCSON PD  
-----

| Test      | g/210L | Time     |
|-----------|--------|----------|
| Air Blank | 0.000  | 08:05:30 |
| Cal Check | 0.099  | 08:05:49 |
| Air Blank | 0.000  | 08:06:18 |

-----  
*Clare Conley*  
Operator's Signature



EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Clare Conley AGENCY Tucson Police Department  
DATE 3-3-11 TIME 0749  
Intoxilyzer Serial # 80-002261 Location ODW Day HB

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0.100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach  $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0.099 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Clare Conley

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000

Location: TUCSON PD

Serial Number: 80-002261

03/03/2011

07:49:19

-----  
Standard Lot#: 719919

Last Changed By: S. CIONE #50616  
-----

QAS: C. CONLEY #34517

TUCSON PD  
-----

| Test      | g/210L | Time     |
|-----------|--------|----------|
| Air Blank | 0.000  | 07:49:47 |
| Cal Check | 0.099  | 07:50:05 |
| Air Blank | 0.000  | 07:50:35 |

*Clare Conley*

\_\_\_\_\_  
Operator's Signature