



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
270 S. Stone Avenue Tucson, AZ 85701-1917

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002261 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

TS4  
3-31-2011

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 3-31-11 Cal 1 0.099 AC Cal 2 0.098 AC

Between the dates of 3-3-11 and 3-31-11 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.  
The instrument was returned to service on N/A.

Clare Conley

Clare L. Conley


3-31-11

ADPS QAS Permit # 33314

SUBSCRIBED AND SWORN to before me this 31, day of March, 20 11.

Laura B. Hodgson

Notary Public

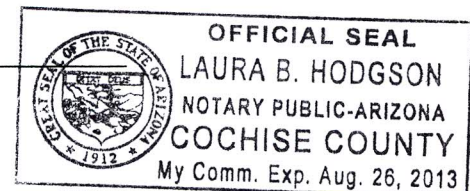

\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1867  
CM  
By  
3-31-2011

QA SPECIALIST Clare Conley AGENCY Tucson Police Department

DATE 3-31-11 TIME 0835

INTOXILYZER SERIAL # 50-002261 LOCATION ODW DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

(✓) 1. Clock time check.

(✓) 2. Date check.

OPERATIONAL TESTS

(✓) 1. Alcohol-free subject test result 0. 000 AC.

(✓) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed

(✓) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.

(✓) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 098 AC

Instrument is operating properly and accurately. YES ✓ NO     

COMMENTS None.  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE Clare Conley

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number : 80-002261  
03/31/2011

08:36:48

Standard Lot#: 20910100A1  
Last Changed By: S. CIONE #50616

QAS: C. CONLEY #34517  
TUCSON PD

Operator: C. CONLEY #34517  
TUCSON PD

Subject: TEST, SQAP  
DOB: 11/11/1111  
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:37:41
Diagnostic Test	Pass	08:38:09
Air Blank	0.000	08:38:36
0.100 Cal Check	0.099	08:38:55
Air Blank	0.000	08:39:24
Subject Test	0.000	08:39:56
Air Blank	ABT*	08:40:02
Air Blank	0.000	08:40:31

\*Sequence Aborted

## Not a Successfully Completed Test Sequence

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number : 80-002261  
03/31/2011

08:40:39

Standard Lot#: 20910100A1  
Last Changed By: S. CIONE #50616

QAS: C. CONLEY #34517  
TUCSON PD

Operator: C. CONLEY #34517  
TUCSON PD

Subject: TEST, SQAP  
DOB: 11/11/1111  
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	08:41:13
Air Blank	IPS**	08:41:20

\*RFI Detect

\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number : 80-002261  
03/31/2011

08:43:18

457  
CM  
TS4

Standard Lot#: 20910100A1  
Last Changed By: S. CIONE #50616

QAS: C. CONLEY #34517  
TUCSON PD

Operator: C. CONLEY #34517  
TUCSON PD

Subject: TEST, SQAP  
DOB: 11/11/1111  
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:44:05
Diagnostic Test	Pass	08:44:33
Air Blank	0.000	08:45:00
0.100 Cal Check	0.099	08:45:19
Air Blank	0.000	08:45:49
Subject Test	0.106*	08:48:57
Air Blank	0.000	08:49:28
Air Blank	0.000	08:49:56
Subject Test	XXX**	08:50:32
Air Blank	0.000	08:51:05

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**



INTOXILYZER 8000

Location: TUCSON PD

Serial Number: 80-002261

03/31/2011

08:51:24

Standard Lot#: 20910100A1

Last Changed By: S. CIONE #50616

QAS: C. CONLEY #34517

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	08:51:51
Cal Check	0.098	08:52:10
Air Blank	0.000	08:52:40

*Claire Conley*

Operator's Signature

50617  
CW

TS4

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Clare Conley AGENCY Tucson Police Department  
DATE 3-31-11 TIME 0835  
Intoxilyzer Serial # 80-002261 Location ODW DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach  $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 099 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Clare Conley

DPS Form Exh G-2 (Rev. 05-01)

↓

INTOXILYZER 8000

Location: TUCSON PD

Serial Number: 80-002261

03/31/2011

08:35:09

757  
cm

-----  
Standard Lot#: 20910100A1

Last Changed By: S. CIONE #50616

Bh

-----  
QAS: C. CONLEY #34517

TUCSON PD

-----  
Test g/210L Time  
-----  
Air Blank 0.000 08:35:36  
Cal Check 0.099 08:35:55  
Air Blank 0.000 08:36:25

*Clare Conley*

\_\_\_\_\_  
Operator's Signature