



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
270 S. Stone Avenue Tucson, AZ 85701-1917

BH  
5.31.2011

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002261 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 5/26/11 Cal 1 0.099 AC Cal 2 0.098 AC

Between the dates of 4/29/11 and 5/26/11 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on —  
The instrument was returned to service on —


Sarah Cione

Sarah L. Cione

5/30/11

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 30, day of MAY, 20 11.

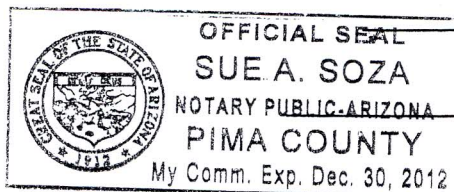

Sue A Soza

Notary Public

\*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

154  
5.31.2011

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department  
DATE 5/26/11 TIME 1035  
INTOXILYZER SERIAL # 80-002261 LOCATION DUI HQ

( / ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- ( / ) 1. Clock time check.  
( / ) 2. Date check.

OPERATIONAL TESTS

- ( / ) 1. Alcohol-free subject test result 0. 000 AC.  
( / ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed  
( / ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.  
( / ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 098 AC

Instrument is operating properly and accurately. YES / NO     

COMMENTS

\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE

Sarah Cione

2 of 75c

TS4

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number : 80-002261  
05/26/2011 10:39:08

Standard Lot#: 20910100A1  
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:39:53
Diagnostic Test	Pass	10:40:20
Air Blank	0.000	10:40:48
0.100 Cal Check	0.099	10:41:06
Air Blank	0.000	10:41:36
Subject Test	0.000	10:42:01
Air Blank	0.000	10:42:28
Five Minute Wait	ABT*	
Air Blank	0.000	10:43:31

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

3 of 7sc

TS4

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number : 80-002261  
05/26/2011 10:43:53

---

Standard Lot#: 20910100A1  
Last Changed By: S. CIONE #50616

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	RFI*	10:44:22
Air Blank	IPS**	10:44:30

\*RFI Detect  
\*\*Improper Sample

**Not a Successfully Completed Test Sequence**



4 of 7sc

INTOXILYZER 8000

Location: TUCSON PD

Serial Number : 80-002261

05/26/2011

10:44:37

1/24

Standard Lot#: 20910100A1

Last Changed By: S. CIONE

#50616

QAS: S. CIONE

#50616

TUCSON PD

Operator: S. CIONE

#50616

TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:45:19
Diagnostic Test	Pass	10:45:46
Air Blank	0.000	10:46:14
0.100 Cal Check	0.098	10:46:33
Air Blank	0.000	10:47:02
Subject Test	0.082*	10:50:10
Air Blank	0.000	10:50:42
Air Blank	0.000	10:51:09
Subject Test	XXX**	10:51:39
Air Blank	0.000	10:52:13

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

5 of 7<sub>sc</sub>

INTOXILYZER 8000

Location: TUCSON PD

Serial Number: 80-002261

05/26/2011

10:52:30

-----  
Standard Lot#: 20910100A1

Last Changed By: S. CIONE #50616  
-----

QAS: S. CIONE #50616

TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	10:52:57
Cal Check	0.098	10:53:16
Air Blank	0.000	10:53:45

*Sarah Cione*

\_\_\_\_\_  
Operator's Signature

6 of 7 sc

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY  
STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

156  
5.31.2011

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department  
DATE 5/26/11 TIME 1035  
Intoxilyzer Serial # 80-002261 Location DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach  $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 099 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000

Location: TUCSON PD

Serial Number: 80-002261

05/26/2011

10:35:27

-----  
Standard Lot#: 20910100A1

Last Changed By: S. CIONE #50616  
-----

QAS: S. CIONE #50616

TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	10:35:54
Cal Check	0.099	10:36:13
Air Blank	0.000	10:36:42

Sarah Cione

Operator's Signature

7 of 7<sub>sc</sub>

Sh