



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002261 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

TS4
7-26-2011

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 7/21/11 Cal 1 0.097 AC Cal 2 0.097 AC

Between the dates of 6/23/11 and 7/21/11
the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -
The instrument was returned to service on -

Sarah Cione

Sarah L. Cione


7/25/11

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 25 day of July, 20 11.

A.M. Hernandez

Notary Public

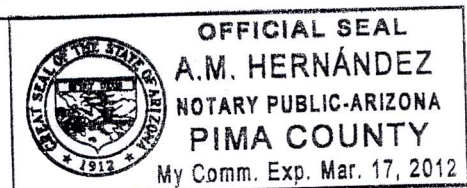

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

1 of 8 SC

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

TS4
7-26-2011

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department
DATE 7/21/11 TIME 950
INTOXILYZER SERIAL # 80-002261 LOCATION DUI HQ

(☒) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (☒) 1. Clock time check.
(☒) 2. Date check.

OPERATIONAL TESTS

- (☒) 1. Alcohol-free subject test result 0. 000 AC.
(☒) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
(☒) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
(☒) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 097 AC

Instrument is operating properly and accurately. YES ☒ NO ☐

COMMENTS

SIGNATURE Sarah Cione

2 of 8 sc

B6

INTOXILYZER 8000

Location: TUCSON PD

Serial Number : 80-002261

07/21/2011

09:52:08

Standard Lot#: 20910100A1

Last Changed By: S. CIONE

#50616

QAS: S. CIONE

#50616

TUCSON PD

Operator: S. CIONE

#50616

TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:53:16
Diagnostic Test	Pass	09:53:43
Air Blank	0.000	09:54:11
0.100 Cal Check	0.098	09:54:29
Air Blank	0.000	09:54:59
Subject Test	0.000	09:55:18
Air Blank	0.000	09:55:45
Five Minute Wait	ABT*	
Air Blank	0.000	09:56:49

*Sequence Aborted

Not a Successfully Completed Test Sequence

3 of 8 sc

PD4

INTOXILYZER 8000

Location: TUCSON PD

Serial Number : 80-002261

07/21/2011

09:57:50

Standard Lot#: 20910100A1

Last Changed By: S. CIONE

#50616

QAS: S. CIONE

#50616

TUCSON PD

Operator: S. CIONE

#50616

TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	09:58:23
Air Blank	IPS**	09:58:31

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000

Location: TUCSON PD

Serial Number : 80-002261

07/21/2011

09:58:59

Standard Lot#: 20910100A1

Last Changed By: S. CIONE

#50616

QAS: S. CIONE

#50616

TUCSON PD

Operator: S. CIONE

#50616

TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:59:49
Diagnostic Test	Pass	10:00:17
Air Blank	0.000	10:00:44
0.100 Cal Check	0.097	10:01:03
Air Blank	0.000	10:01:32
Subject Test	0.046*	10:04:41
Air Blank	0.000	10:05:12
Air Blank	0.000	10:05:39
Subject Test	0.055**	10:06:23
Air Blank	0.000	10:07:01
Five Minute Wait		
Air Blank	0.000	10:11:51
Subject Test	0.000**	10:12:48
Air Blank	0.000	10:13:16

*Deficient Sample

**No .020 agreement -- Registered actual reading during mouth alcohol check. Waited + performed another subject test to demonstrate no alcohol present.

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number : 80-002261
07/21/2011

10:13:29

5 of 8_{sc}

BH

Standard Lot#: 20910100A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:14:12
Diagnostic Test	Pass	10:14:40
Air Blank	0.000	10:15:07
0.100 Cal Check	0.097	10:15:26
Air Blank	0.000	10:15:55
Subject Test	XXX*	10:16:22
Air Blank	0.000	10:16:55

*Invalid Sample
Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002261
07/21/2011

10:17:34

6 of 8 sc

Standard Lot#: 20910100A1
Last Changed By: S. CIONE #50616

TS4

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	10:18:01
Cal Check	0.097	10:18:20
Air Blank	0.000	10:18:49



Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 7/21/11 TIME 950
Intoxilyzer Serial # 80-002261 Location DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 097 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002261
07/21/2011

09:50:21

Standard Lot#: 20910100A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:50:48
Cal Check	0.097	09:51:07
Air Blank	0.000	09:51:37

Sarah Cione

Operator's Signature

8 of 8 sc

154