



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002261 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

TS4
8.22.2011

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 8-18-11 Cal 1 0.098 AC Cal 2 0.098 AC

Between the dates of 7.21-11 and 8-18-11 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.
The instrument was returned to service on N/A.

Clare Conley

Clare L. Conley


8-18-11

ADPS QAS Permit # 33314

SUBSCRIBED AND SWORN to before me this 18th, day of August, 20 11.

Susan A. Hitchye

Notary Public

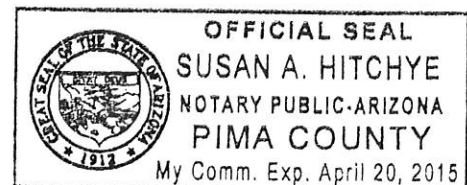

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1610
Ch
B4
8.22.2011

QA SPECIALIST Clare Conley AGENCY Tucson Police Department

DATE 8-18-11 TIME 0816

INTOXILYZER SERIAL # 80-002261 LOCATION ODW Det HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

(✓) 1. Clock time check.

(✓) 2. Date check.

OPERATIONAL TESTS

(✓) 1. Alcohol-free subject test result 0. 000 AC.

(✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed

(✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.

(✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 098 AC

Instrument is operating properly and accurately. YES ✓ NO

COMMENTS None.

SIGNATURE Clare Conley

INTOXILYZER 8000
Location: TUCSON PD
Serial Number : 80-002261
08/18/2011

08:17:58

2 of 10
CW
B4

Standard Lot#: 20910100A1
Last Changed By: S. CIONE #50616

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP

DOB: 11/11/1111

Sex: M

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:18:43
Diagnostic Test	Pass	08:19:10
Air Blank	0.000	08:19:38
0.100 Cal Check	0.098	08:19:56
Air Blank	0.000	08:20:26
Subject Test	0.000	08:21:04
Air Blank	ABT*	08:21:11
Air Blank	0.000	08:21:39

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number : 80-002261
08/18/2011

08:23:20

Standard Lot#: 20910100A1
Last Changed By: S. CIONE #50616

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP

DOB: 11/11/1111

Sex: M

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	08:23:46
Air Blank	IPS**	08:23:56

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number : 80-002261
08/18/2011

08:24:38

Standard Lot#: 20910100A1
Last Changed By: S. CIONE #50616

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP

DOB: 11/11/1111

Sex: M

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:25:23
Diagnostic Test	Pass	08:25:51
Air Blank	0.000	08:26:18
0.100 Cal Check	0.098	08:26:37
Air Blank	0.000	08:27:06
Subject Test	0.114*	08:30:15
Air Blank	0.000	08:30:46
Air Blank	0.000	08:31:13
Subject Test	RGE**	08:31:48
Air Blank	0.000	08:32:23

*Deficient Sample

**Range Exceeded

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number : 80-002261
08/18/2011

08:32:32

Standard Lot#: 20910100A1
Last Changed By: S. CIONE #50616

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP

DOB: 11/11/1111

Sex: M

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:33:18
Diagnostic Test	Pass	08:33:45
Air Blank	0.000	08:34:13
0.100 Cal Check	0.098	08:34:31
Air Blank	0.000	08:35:01
Subject Test	XXX*	08:35:36
Air Blank	0.000	08:36:10

*Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002261
08/18/2011 08:37:01

b610
cu

Standard Lot#: 20910100A1
Last Changed By: S. CIONE #50616

QAS: C. CONLEY #34517
TUCSON PD

134

Test	g/210L	Time
Air Blank	0.000	08:37:28
Cal Check	0.098	08:37:47
Air Blank	0.000	08:38:16

Claire Conley
Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

78610
Ch

TS4
8.22.2011

QA Specialist Clare Conley AGENCY Tucson Police Department
DATE 8-18-11 TIME 0816
Intoxilyzer Serial # 80-802261 Location ODW DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 098 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Clare Conley

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000

Location: TUCSON PD

Serial Number: 80-002261

08/18/2011

08:16:38


Standard Lot#: 20910100A1

Last Changed By: S. CIONE #50616

QAS: C. CONLEY #34517

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	08:17:05
Cal Check	0.098	08:17:24
Air Blank	0.000	08:17:53

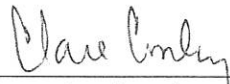

Operator's Signature

85610
CIC

134



TUCSON POLICE DEPARTMENT
CRIME LABORATORY

9810
Cw

INTOXILYZER SERVICE RECORD
FOR
CALIBRATION GAS TANK CHANGE

BS4
8.22.2011

Date: 8-18-11
QAS: Clare Conley
Intoxilyzer Serial Number: 80-802261
New Lot Number: 04911100A1
New Expiration Date: 03/01/2013

- (✓) 1. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 2. Remove the old gas tank and install the new one.
- (✓) 3. Hit the "ESC" key twice and swipe the QAS card.
- (✓) 4. Go to the "M-Maintenance" menu and select "S-Standard Change."
Select "D" for dry gas and enter the following information:
Enter Target Value: 0.100
Enter Control Lot #: From the tank label
Enter Exp Date: From the tank label (mm/dd/yyyy)
Edit info Y/N? Y or N
Update Std Info? Y ENTER
"Please wait saving..."
- (✓) 5. Hit the "ESC" key twice.
- (✓) 6. Go to the "M-Maintenance" menu and select "Control Testing."
Select "D" for dry control test. After selection is made press ENTER.
- (✓) 7. Air blank completed.
- (✓) 8. Calibration check completed. Test results 0. 097 AC.
- (✓) 9. Air blank completed.
- (✓) 10. Remove printed record. Attach the record to the completed checklist.

SIGNATURE


Clare Conley

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002261
08/18/2011 08:45:06

Standard Lot#: 04911100A1
Last Changed By: C. CONLEY #34517

QAS: C. CONLEY #34517
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	08:45:33
Cal Check	0.097	08:45:52
Air Blank	0.000	08:46:22

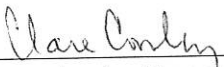

Operator's Signature