



CERTIFICATION by CUSTODIAN of RECORDS
 Tucson Police Department Crime Laboratory
 1306 W. Miracle Mile, Tucson, AZ 85705

*TSK
 11-7-2011*

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002261 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 10/14/11 Cal 1 0.097 AC Cal 2 0.096 AC

Between the dates of 9/15/11 and 10/14/11 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.
 The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione


10/27/11

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 27, day of October, 2011.

A.M. Hernandez

Notary Public

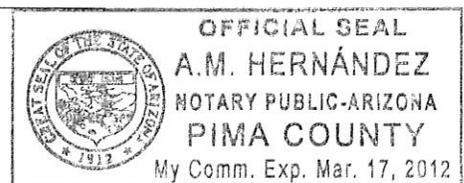

*Includes documents specified in A.A.C. R13-10-104.A.5

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

B4
11.7.2011

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 10/14/11 TIME 1326

INTOXILYZER SERIAL # 80-002261 LOCATION DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 096 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Sarah Cione

2 of 7 sc

INTOXILYZER 8000
Location: TUCSON PD
Serial Number : 80-002261
10/14/2011 13:27:48

TSY

Standard Lot#: 04911100A1
Last Changed By: C. CONLEY #34517

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|------------------|--------|----------|
| Air Blank | 0.000 | 13:28:38 |
| Diagnostic Test | Pass | 13:29:05 |
| Air Blank | 0.000 | 13:29:32 |
| 0.100 Cal Check | 0.097 | 13:29:51 |
| Air Blank | 0.000 | 13:30:21 |
| Subject Test | 0.000 | 13:30:39 |
| Air Blank | 0.000 | 13:31:06 |
| Five Minute Wait | ABT* | |
| Air Blank | 0.000 | 13:31:44 |

*Sequence Aborted

Not a Successfully Completed Test Sequence

3 of 7_{sc}

134

INTOXILYZER 8000
Location: TUCSON PD
Serial Number : 80-002261
10/14/2011 13:32:02

Standard Lot#: 04911100A1
Last Changed By: C. CONLEY #34517

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------|--------|----------|
| Air Blank | RFI* | 13:32:29 |
| Air Blank | IPS** | 13:32:36 |

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

4 of 7 sc

B4

INTOXILYZER 8000
 Location: TUCSON PD
 Serial Number : 80-002261
 10/14/2011 13:32:55

Standard Lot#: 04911100A1
 Last Changed By: C. CONLEY #34517

QAS: S. CIONE #50616
 TUCSON PD

Operator: S. CIONE #50616
 TUCSON PD

Subject: SQAP, SQAP
 DOB: 11/11/1111
 Sex: F Weight: 111

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 13:33:45 |
| Diagnostic Test | Pass | 13:34:12 |
| Air Blank | 0.000 | 13:34:40 |
| 0.100 Cal Check | 0.098 | 13:34:59 |
| Air Blank | 0.000 | 13:35:28 |
| Subject Test | 0.119* | 13:38:36 |
| Air Blank | 0.000 | 13:39:09 |
| Air Blank | 0.000 | 13:39:36 |
| Subject Test | XXX** | 13:40:10 |
| Air Blank | 0.000 | 13:40:43 |

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002261
10/14/2011

13:41:00

5 of 7sc

TS4

Standard Lot#: 04911100A1
Last Changed By: C. CONLEY #34517

QAS: S. CIONE #50616
TUCSON PD

| Test | g/210L | Time |
|-----------|--------|----------|
| Air Blank | 0.000 | 13:41:27 |
| Cal Check | 0.096 | 13:41:46 |
| Air Blank | 0.000 | 13:42:16 |

Sarah Cione

Operator's Signature

6 of 7 sc

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

T54
11.7.2011

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 10/14/11 TIME 1326
Intoxilyzer Serial # ~~80-400~~ 5C Location DUI HQ
80-002261

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 097 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002261
10/14/2011 13:26:20

7 of 7 sc

TS4

Standard Lot#: 04911100A1
Last Changed By: C. CONLEY #34517

QAS: S. CIONE #50616
TUCSON PD

| Test | g/210L | Time |
|-----------|--------|----------|
| Air Blank | 0.000 | 13:26:47 |
| Cal Check | 0.097 | 13:27:06 |
| Air Blank | 0.000 | 13:27:35 |

Sarah Cione

Operator's Signature