



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
1306 W. Miracle Mile, Tucson, AZ 85705

CH  
11-30-11

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002261 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 11/10/11 Cal 1 0.098 AC Cal 2 0.097 AC

Between the dates of 10/14/11 and 11/10/11 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on —.  
The instrument was returned to service on —.

Sarah Cione

Sarah L. Cione


11/16/11

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 16, day of November, 20—.

A.M. Hernandez

Notary Public

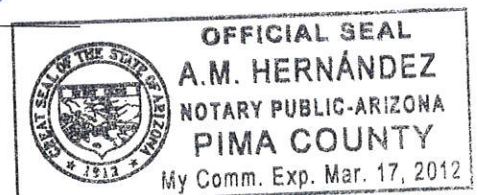

\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department  
DATE 11/10/11 TIME 924  
INTOXILYZER SERIAL # 80-002264 LOCATION DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.  
(✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.  
(✓) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed  
(✓) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.  
(✓) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 097 AC

Instrument is operating properly and accurately. YES ✓ NO       

COMMENTS

SIGNATURE Sarah Cione

2 of 7<sub>sc</sub>

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number : 80-002261  
11/10/2011

09:26:50

CM

---

Standard Lot#: 04911100A1  
Last Changed By: C. CONLEY #34517

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	09:27:37
Diagnostic Test	Pass	09:28:05
Air Blank	0.000	09:28:32
0.100 Cal Check	0.098	09:28:51
Air Blank	0.000	09:29:20
Subject Test	0.000	09:29:45
Air Blank	0.000	09:30:12
Five Minute Wait	ABT*	
Air Blank	0.000	09:30:54

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

3 of 7<sub>sc</sub>

W

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number : 80-002261  
11/10/2011

09:34:04

---

Standard Lot#: 04911100A1  
Last Changed By: C. CONLEY #34517

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

---

Test	g/210L	Time
Air Blank	0.000	09:35:01
Diagnostic Test	Pass	09:35:29
Air Blank	RFI*	09:35:37
Air Blank	IPS**	09:35:45

---

\*RFI Detect

\*\*Improper Sample

**Not a Successfully Completed Test Sequence**



4 of 7<sub>sc</sub>

CM

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number : 80-002261  
11/10/2011

09:36:42

---

Standard Lot#: 04911100A1  
Last Changed By: C. CONLEY #34517

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	09:37:29
Diagnostic Test	Pass	09:37:56
Air Blank	0.000	09:38:24
0.100 Cal Check	0.098	09:38:42
Air Blank	0.000	09:39:12
Subject Test	0.126*	09:42:20
Air Blank	0.000	09:42:52
Air Blank	0.000	09:43:20
Subject Test	XXX**	09:45:15
Air Blank	0.000	09:45:49

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000

Location: TUCSON PD

Serial Number: 80-002261

11/10/2011

09:47:54

-----  
Standard Lot#: 04911100A1

Last Changed By: C. CONLEY #34517  
-----

QAS: S. CIONE #50616

TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	09:48:21
Cal Check	0.097	09:48:40
Air Blank	0.000	09:49:09

  
Operator's Signature

EXHIBIT G-2

6 of 7<sub>sc</sub>

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

DM  
11-30-11

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department  
DATE 11/10/11 TIME 925  
Intoxilyzer Serial # 80-002261 Location DUT HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 098 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000

Location: TUCSON PD

Serial Number: 80-002261

11/10/2011

09:25:04


-----  
Standard Lot#: 04911100A1

Last Changed By: C. CONLEY #34517  
-----

QAS: S. CIONE #50616

TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	09:25:31
Cal Check	0.098	09:25:50
Air Blank	0.000	09:26:20

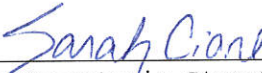  
\_\_\_\_\_  
Operator's Signature