



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

TS4
2-11-2014

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002261 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 1/29/14 Cal 1 0.101 AC Cal 2 0.100 AC

Between the dates of 1/2/14 and 1/29/14 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.
The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione


2/11/14

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 11th, day of February, 20 14.

Elma Escarcega

Notary Public

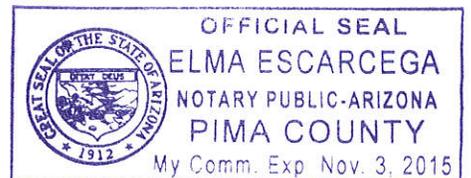

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

1 of 8 sc
TS4
2.11.2014

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 1/29/14 TIME 1549

INTOXILYZER SERIAL # 80-002261 LOCATION Crime Lab

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 100 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Sarah Cione

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002261
Core Version: 8105.48
01/29/2014

15:51:10

2 of 8
sc
Tsh

Standard Lot#: 30412100A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	15:51:50
Diagnostic Test	Pass	15:52:25
Air Blank	0.000	15:52:53
0.100 Cal Check	0.101	15:53:13
Air Blank	0.000	15:53:43
Subject Test	0.000	15:54:01
Air Blank	0.000	15:54:28
Wait	ABT*	
Air Blank	0.000	15:55:02

*Sequence Aborted

Not a Successfully Completed Test Sequence

3 of 8 sc
134

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002261
Core Version: 8105.48
01/29/2014 15:55:06

Standard Lot#: 30412100A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	15:55:49
Diagnostic Test	Pass	15:56:24
Air Blank	0.000	15:56:51
0.100 Cal Check	RFI*	15:57:06
Air Blank	IPS**	15:57:13

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

4 of 8 sec
T54

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002261
Core Version: 8105.48
01/29/2014 15:57:18

Standard Lot#: 30412100A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	15:58:05
Diagnostic Test	Pass	15:58:40
Air Blank	0.000	15:59:07
0.100 Cal Check	0.101	15:59:28
Air Blank	0.000	15:59:58
Subject Test	0.044*	16:03:05
Air Blank	0.000	16:03:36
Air Blank	0.000	16:04:03
Subject Test	XXX**	16:04:19
Air Blank	0.000	16:04:52

*Deficient Sample
**Invalid Sample
Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002261
Core Version: 8105.48
01/29/2014 16:05:03

5 of 8 sc
TS4

Standard Lot#: 30412100A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	16:05:32
0.100 Cal Check	0.100	16:05:53
Air Blank	0.000	16:06:22

Preventative Maintenance
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6 of 8 sc
TS4
2.11.2014

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 1/29/14 TIME 1549
Intoxilyzer Serial # 80-002261 Location Crime Lab

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 101 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002261
Core Version: 8105.48
01/29/2014

15:49:39

7 of 8
TS4 SC

Standard Lot#: 30412100A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	15:50:08
0.100 Cal Check	0.101	15:50:28
Air Blank	0.000	15:50:58

Preventative Maintenance
Performed -- Timer Reset



TUCSON POLICE DEPARTMENT
CRIME LABORATORY

INTOXILYZER SERVICE RECORD

8 of 8 sc
TBH
2-11-2014

Date: 1/29/14

QAS: Sarah Cione

Intoxilyzer Serial Number: 80-002261

Reason For Service: Officer Nolan notified me that this intox would not go to "Ready Mode"

Date Out of Service: N/A

Calibration Check completed prior to removal from service

Date: N/A

SQAP completed prior to removal from service

Date: N/A

Corrective Action: Sent to CMI, Inc.

Breath hose would not heat up. Replaced breath hose.