



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

TS4
6.3.2014

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002261 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 5/22/14 Cal 1 0.100 AC Cal 2 0.100 AC

Between the dates of 4/24/14 and 5/22/14 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.
The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione

5/29/14


ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 29th, day of May, 2014.

Elma Escarcega

Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

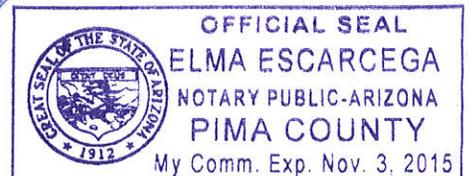

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

1 of 7
TSG
6.3.2014

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 5/22/14 TIME 1044

INTOXILYZER SERIAL # 80-002261 LOCATION DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 100 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Sarah Cione

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002261
Core Version: 8105.48
05/22/2014 10:48:40

2 of 7
TS4

Standard Lot#: 25913100A2
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:49:38
Diagnostic Test	Pass	10:50:13
Air Blank	0.000	10:50:40
0.100 Cal Check	0.101	10:51:01
Air Blank	0.000	10:51:31
Subject Test	0.000	10:51:54
Air Blank	0.000	10:52:21
Wait	ABT*	
Air Blank	0.000	10:53:30

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002261
Core Version: 8105.48
05/22/2014

10:55:39

3 of 7_{sc}
TS₄

Standard Lot#: 25913100A2
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:56:37
Diagnostic Test	Pass	10:57:12
Air Blank	0.000	10:57:39
0.100 Cal Check	0.100	10:58:00
Air Blank	0.000	10:58:30
Subject Test	RFI*	10:58:38
Air Blank	IPS**	10:58:47

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002261
Core Version: 8105.48
05/22/2014

11:02:31

4 of 7_{sc}
TS4

Standard Lot#: 25913100A2
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:03:35
Diagnostic Test	Pass	11:04:10
Air Blank	0.000	11:04:38
0.100 Cal Check	0.100	11:04:58
Air Blank	0.000	11:05:28
Subject Test	0.311*	11:08:35
Air Blank	0.000	11:09:07
Air Blank	0.000	11:09:35
Subject Test	XXX**	11:09:54
Air Blank	0.000	11:10:28

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002261
Core Version: 8105.48
05/22/2014 11:10:53

5 of 7 sc

TJ4

Standard Lot#: 25913100A2
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	11:11:22
0.100 Cal Check	0.100	11:11:43
Air Blank	0.000	11:12:12

Preventative Maintenance
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

6 of 7
sc

ARIZONA DEPARTMENT OF PUBLIC SAFETY

TS4

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

6.3.2014

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 5/22/14 TIME 1044
Intoxilyzer Serial # 80-002261 Location DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002261
Core Version: 8105.48
05/22/2014 10:44:51

7 of 75c
TS4

Standard Lot#: 25913100A2
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	10:45:20
0.100 Cal Check	0.100	10:45:40
Air Blank	0.000	10:46:10

Preventative Maintenance
Performed -- Timer Reset