



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
1306 W. Miracle Mile, Tucson, AZ 85705

TSG  
9.29.2014

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002261 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 9/18/14 Cal 1 0.098 AC Cal 2 0.097 AC

Between the dates of 9/9/14 and 9/18/14 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.  
The instrument was returned to service on -.


Sarah Cione  
Sarah L. Cione

9/23/14  
ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 23<sup>rd</sup>, day of September, 2014.

Elma Escarcega  
Notary Public

\*Includes documents specified in A.A.C. R13-10-104.A.5

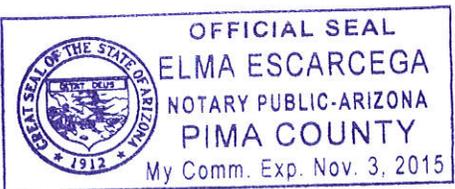

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

1 of 7  
TB4  
9-29-2014

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 9/18/14 TIME 907

INTOXILYZER SERIAL # 80-002261 LOCATION Crime Lab

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 097 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS Intox checked 9 days after 9/9/14 to get it on different QAS  
schedule. Intox re-located to Team location.

SIGNATURE Sarah Cione

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002261  
Core Version: 8105.48  
09/18/2014

09:08:53

2 of 7<sub>sc</sub>  
T34

---

Standard Lot#: 25913100A2  
Last Changed By: S. CIONE #50616

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1990  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

| Test            | g/210L | Time     |
|-----------------|--------|----------|
| Air Blank       | 0.000  | 09:09:37 |
| Diagnostic Test | Pass   | 09:10:12 |
| Air Blank       | 0.000  | 09:10:39 |
| 0.100 Cal Check | 0.098  | 09:11:00 |
| Air Blank       | 0.000  | 09:11:29 |
| Subject Test    | 0.000  | 09:11:45 |
| Air Blank       | 0.000  | 09:12:12 |
| Wait            | ABT*   |          |
| Air Blank       | 0.000  | 09:12:47 |

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002261  
Core Version: 8105.48  
09/18/2014

09:12:51

3 of 7 sc  
154

---

Standard Lot#: 25913100A2  
Last Changed By: S. CIONE #50616

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1999  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

| Test            | g/210L | Time     |
|-----------------|--------|----------|
| Air Blank       | 0.000  | 09:13:51 |
| Diagnostic Test | Pass   | 09:14:26 |
| Air Blank       | 0.000  | 09:14:53 |
| 0.100 Cal Check | 0.098  | 09:15:14 |
| Air Blank       | 0.000  | 09:15:44 |
| Reference       | RFI*   | 09:15:48 |
| Air Blank       | IPS**  | 09:16:00 |

\*RFI Detect

\*\*Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002261  
Core Version: 8105.48  
09/18/2014

09:16:06

4 of 7  
SC  
64

Standard Lot#: 25913100A2  
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1999  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

| Test            | g/210L | Time     |
|-----------------|--------|----------|
| Air Blank       | 0.000  | 09:16:48 |
| Diagnostic Test | Pass   | 09:17:23 |
| Air Blank       | 0.000  | 09:17:50 |
| 0.100 Cal Check | 0.097  | 09:18:11 |
| Air Blank       | 0.000  | 09:18:40 |
| Subject Test    | 0.241* | 09:21:48 |
| Air Blank       | 0.000  | 09:22:19 |
| Air Blank       | 0.000  | 09:22:47 |
| Subject Test    | XXX**  | 09:23:02 |
| Air Blank       | 0.000  | 09:23:35 |

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

### Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002261  
Core Version: 8105.48  
09/18/2014 09:24:00

5 of 7  
SC  
T84

-----  
Standard Lot#: 25913100A2  
Last Changed By: S. CIONE #50616  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Operator: S. CIONE #50616  
TUCSON PD  
-----

| Test            | g/210L | Time     |
|-----------------|--------|----------|
| Air Blank       | 0.000  | 09:24:28 |
| 0.100 Cal Check | 0.097  | 09:24:49 |
| Air Blank       | 0.000  | 09:25:19 |

Preventative Maintenance  
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6 of 7<sub>sc</sub>

TJG  
9.29.2014

QA Specialist Sarah Cione AGENCY Tucson Police Department  
DATE 9/18/14 TIME 907  
Intoxilyzer Serial # 80-002261 Location Crime Lab

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 098 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002261  
Core Version: 8105.48  
09/18/2014 09:07:24

7 of 7 sc  
Dy

-----  
Standard Lot#: 25913100A2  
Last Changed By: S. CIONE #50616  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Operator: S. CIONE #50616  
TUCSON PD  
-----

| Test            | g/210L | Time     |
|-----------------|--------|----------|
| Air Blank       | 0.000  | 09:07:52 |
| 0.100 Cal Check | 0.098  | 09:08:13 |
| Air Blank       | 0.000  | 09:08:43 |

Preventative Maintenance  
Performed -- Timer Reset