



CERTIFICATION by CUSTODIAN of RECORDS
 Tucson Police Department Crime Laboratory
 1306 W. Miracle Mile, Tucson, AZ 85705

TS4
 2-10-2015

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002261 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 2-6-2015 Cal 1 0.100 AC Cal 2 0.099 AC

Between the dates of 1-9-2015 and 2-6-2015 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on —.
 The instrument was returned to service on —.

[Handwritten signature of Amy L. Autret]

Amy L. Autret


2-10-2015

ADPS QAS Permit # 33718

SUBSCRIBED AND SWORN to before me this 10th, day of February, 2015.

[Handwritten signature of Elma Escarcega]

Notary Public

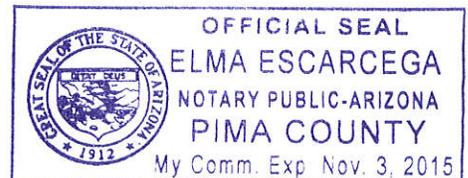

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

1 of 7
MS
TJG
2-10-2015

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Amy Autret 53546 AGENCY Tucson Police Department

DATE 2-6-2015 TIME 1105

INTOXILYZER SERIAL # 80-002261 LOCATION ODE

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS


- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 099 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

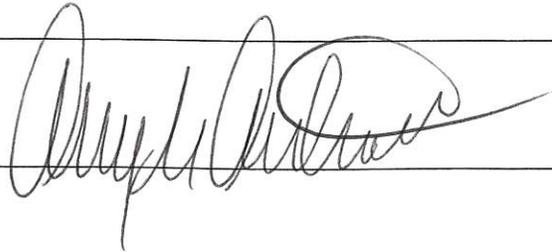SIGNATURE 

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

2 of 7
JLS
154

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002261
Core Version: 8105.48
02/06/2015

11:08:55

Standard Lot#: 25913100A2
Last Changed By: S. CIONE #50616

QAS: A. AUTRET #53546
TUCSON PD

Operator: A. AUTRET #53546
TUCSON PD

Subject: TEST,TEST
DOB: 11/11/1987
Sex: F Weight: 135

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 11:10:15 |
| Diagnostic Test | Pass | 11:10:50 |
| Air Blank | 0.000 | 11:11:17 |
| 0.100 Cal Check | 0.099 | 11:11:38 |
| Air Blank | 0.000 | 11:12:07 |
| Subject Test | 0.000 | 11:12:30 |
| Air Blank | 0.000 | 11:12:57 |
| Wait | ABT* | |
| Air Blank | 0.000 | 11:13:35 |

*Sequence Aborted

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

3 of 7
JLJ
B4

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002261
Core Version: 8105.48
02/06/2015 11:14:55

Standard Lot#: 25913100A2
Last Changed By: S. CIONE #50616

QAS: A. AUTRET #53546
TUCSON PD

Operator: A. AUTRET #53546
TUCSON PD

Subject: TEST,SQAP
DOB: 11/12/1987
Sex: F Weight: 135

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 11:16:13 |
| Diagnostic Test | Pass | 11:16:48 |
| Air Blank | 0.000 | 11:17:16 |
| 0.100 Cal Check | RFI* | 11:17:28 |
| Air Blank | IPS** | 11:17:38 |

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

4 of 7
old
TS4

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002261
Core Version: 8105.48
02/06/2015

11:18:24

Standard Lot#: 25913100A2
Last Changed By: S. CIONE #50616

QAS: A. AUTRET #53546
TUCSON PD

Operator: A. AUTRET #53546
TUCSON PD

Subject: TEST,SQAP
DOB: 11/12/1987
Sex: F Weight: 135

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 11:19:35 |
| Diagnostic Test | Pass | 11:20:10 |
| Air Blank | 0.000 | 11:20:37 |
| 0.100 Cal Check | 0.100 | 11:20:58 |
| Air Blank | 0.000 | 11:21:27 |
| Subject Test | 0.188* | 11:24:35 |
| Air Blank | 0.000 | 11:25:07 |
| Air Blank | 0.000 | 11:25:34 |
| Subject Test | XXX** | 11:26:06 |
| Air Blank | 0.000 | 11:26:39 |

*Deficient Sample
**Invalid Sample
Begin new deprivation period

Not a Successfully Completed Test Sequence

Pursuant to section 23-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002261
Core Version: 8105.48
02/06/2015 11:28:34

5 of 7
AS
B4

Standard Lot#: 25913100A2
Last Changed By: S. CIONE #50616

QAS: A. AUTRET #53546
TUCSON PD

Operator: A. AUTRET #53546
TUCSON PD

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 11:29:07 |
| 0.100 Cal Check | 0.099 | 11:29:28 |
| Air Blank | 0.000 | 11:29:57 |

Preventative Maintenance
Performed -- Timer Reset

EXHIBIT G-2

6 of 7
JLH
TBH
2.10.2015

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY


STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Amy Autret 53546 AGENCY Tucson Police Department
DATE 2-6-2015 TIME 1105
Intoxilyzer Serial # 80-002261 Location ODE

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

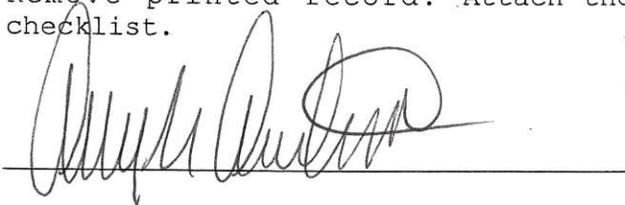

DPS Form Exh G-2 (Rev. 05-01)

7 of 7
all
Bh

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002261
Core Version: 8105.48
02/06/2015 11:05:32

Standard Lot#: 25913100A2
Last Changed By: S. CIONE #50616

QAS: A. AUTRET #53546
TUCSON PD

Operator: A. AUTRET #53546
TUCSON PD

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 11:06:05 |
| 0.100 Cal Check | 0.100 | 11:06:26 |
| Air Blank | 0.000 | 11:06:56 |

Preventative Maintenance
Performed -- Timer Reset