



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
1306 W. Miracle Mile, Tucson, AZ 85705

TSy  
3.11.2015

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002261 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 3-6-15 Cal 1 0.100 AC Cal 2 0.098 AC

Between the dates of 2-6-15 and 3-6-15 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.  
The instrument was returned to service on N/A.

[Signature]  
Amy L. Autret


3-9-15

ADPS QAS Permit # 33718

SUBSCRIBED AND SWORN to before me this 9<sup>th</sup>, day of March, 2015.

[Signature]

Notary Public

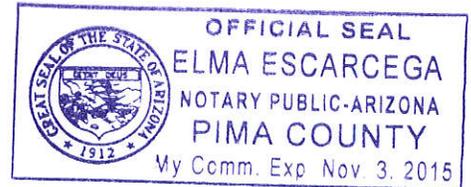

\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

1 of 7  
MA  
TSh  
3.11.2015

ARIZONA DEPARTMENT OF PUBLIC SAFETY  
STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000  
STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Amy Autret 53546 AGENCY Tucson Police Department

DATE 3-6-15 TIME 0902

INTOXILYZER SERIAL # 80-002261 LOCATION ODE

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 098 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE



Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

2 of 7  
ad

TS4

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002261  
Core Version: 8105.48  
03/06/2015

09:07:23

Standard Lot#: 25913100A2  
Last Changed By: S. CIONE #50616

QAS: A. AUTRET #53546  
TUCSON PD

Operator: A. AUTRET #53546  
TUCSON PD

Subject: SOAP, SOAP  
DOB: 11/11/1987  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:08:35
Diagnostic Test	Pass	09:09:10
Air Blank	0.000	09:09:37
0.100 Cal Check	0.100	09:09:58
Air Blank	0.000	09:10:27
Subject Test	0.000	09:10:48
Air Blank	0.000	09:11:15
Wait	ABT*	
Air Blank	0.000	09:11:50

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

3 of 7  
JLH  
TSh

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002261  
Core Version: 8105.48  
03/06/2015 09:12:14

Standard Lot#: 25913100A2  
Last Changed By: S. CIONE #50616

QAS: A. AUTRET #53546  
TUCSON PD

Operator: A. AUTRET #53546  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:13:17
Diagnostic Test	Pass	09:13:52
Air Blank	0.000	09:14:19
0.100 Cal Check	RFI*	09:14:32
Air Blank	IPS**	09:14:40

\*RFI Detect  
\*\*Improper Sample

### Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

4 of 7  
[Signature]  
B4

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002261  
Core Version: 8105.48  
03/06/2015 09:15:21

Standard Lot#: 25913100A2  
Last Changed By: S. CIONE #50616

QAS: A. AUTRET #53546  
TUCSON PD

Operator: A. AUTRET #53546  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:16:27
Diagnostic Test	Pass	09:17:02
Air Blank	0.000	09:17:29
0.100 Cal Check	0.099	09:17:50
Air Blank	0.000	09:18:20
Subject Test	0.268*	09:21:27
Air Blank	0.000	09:21:59
Air Blank	0.000	09:22:26
Subject Test	XXX**	09:22:45
Air Blank	0.000	09:23:19

\*Deficient Sample  
\*\*Invalid Sample  
Begin new deprivation period

**Not a Successfully Completed Test Sequence**

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

5 of 7  
MS  
734

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002261  
Core Version: 8105.48  
03/06/2015 09:24:21

-----  
Standard Lot#: 25913100A2  
Last Changed By: S. CIONE #50616  
-----

QAS: A. AUTRET #53546  
TUCSON PD  
-----

Operator: A. AUTRET #53546  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	09:24:51
0.100 Cal Check	0.098	09:25:12
Air Blank	0.000	09:25:42

Preventative Maintenance  
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

6 of 7  
JLB  
TS4  
3.11.2015

ARIZONA DEPARTMENT OF PUBLIC SAFETY


STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Amy Autret 53546 AGENCY Tucson Police Department  
DATE 3-6-15 TIME 0902  
Intoxilyzer Serial # 80-002261 Location ODE

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 100 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

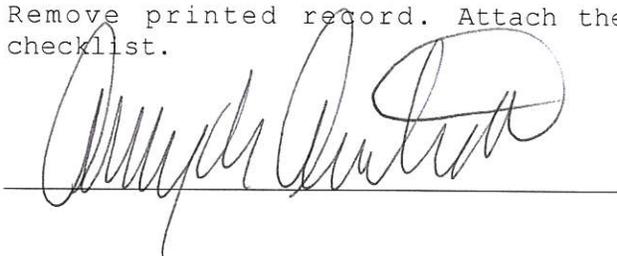

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002261  
Core Version: 8105.48  
03/06/2015 09:02:15

7 of 7  
old  
TSK

-----  
Standard Lot#: 25913100A2  
Last Changed By: S. CIONE #50616  
-----

QAS: A. AUTRET #53546  
TUCSON PD  
-----

Operator: A. AUTRET #53546  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	09:02:53
0.100 Cal Check	0.100	09:03:14
Air Blank	0.000	09:03:44

Preventative Maintenance  
Performed -- Timer Reset