



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002261 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

TS4
4.10.2015

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 4/2/2015 Cal 1 0.101 AC Cal 2 0.099 AC

Between the dates of 3/6/2015 and 4/2/2015 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on
The instrument was returned to service on

William L Loomis

William L. Loomis


4/2/2015

ADPS QAS Permit # 33720

SUBSCRIBED AND SWORN to before me this 2nd day of April, 2015.

Elma Escarcega

Notary Public

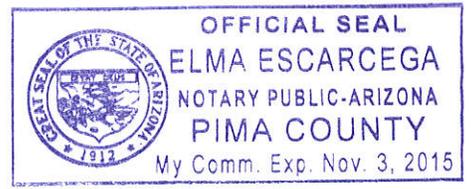

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1 of 7
WLL
TS4 4.10.2015

QA SPECIALIST W. LOOMIS 53945 AGENCY Tucson Police Department

DATE 4/2/2015 TIME 09:50

INTOXILYZER SERIAL # 80-002261 LOCATION ODE

() 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- () 1. Clock time check.
- () 2. Date check.

OPERATIONAL TESTS

- () 1. Alcohol-free subject test result 0. 000 AC.
- () 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- () 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- () 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 099 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE

Will Loomis

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002261
Core Version: 8105.48
04/02/2015 09:52:32

2 of 7
WLV
Bh

Standard Lot#: 25913100A2
Last Changed By: S. CIONE #50616

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Subject: SQAP,TEST
DOB: 11/11/1987
Sex: M Weight: 222

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:53:21
Diagnostic Test	Pass	09:53:56
Air Blank	0.000	09:54:23
0.100 Cal Check	0.100	09:54:44
Air Blank	0.000	09:55:14
Subject Test	0.000	09:55:33
Air Blank	0.000	09:56:01
Wait	ABT*	
Air Blank	0.000	09:56:37

*Sequence Aborted

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002261
Core Version: 8105.48
04/02/2015

09:56:42

3 of 7
WLL
TS4

Standard Lot#: 25913100A2
Last Changed By: S. CIONE #50616

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Subject: SQAP,TEST
DOB: 11/11/1987
Sex: M Weight: 222

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:57:31
Diagnostic Test	Pass	09:58:06
Air Blank	0.000	09:58:34
0.100 Cal Check	RFI*	09:58:47
Air Blank	IPS**	09:58:56

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002261
Core Version: 8105.48
04/02/2015

09:59:00

4 of 7
wll
Bh

Standard Lot#: 25913100A2
Last Changed By: S. CIONE #50616

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Subject: SQAP,TEST
DOB: 11/11/1987
Sex: M Weight: 222

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:59:49
Diagnostic Test	Pass	10:00:24
Air Blank	0.000	10:00:51
0.100 Cal Check	0.100	10:01:12
Air Blank	0.000	10:01:41
Subject Test	0.215*	10:04:49
Air Blank	0.000	10:05:22
Air Blank	0.000	10:05:49
Subject Test	XXX**	10:06:06
Air Blank	0.000	10:06:40

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002261
Core Version: 8105.48
04/02/2015 10:07:05

Standard Lot#: 25913100A2
Last Changed By: S. CIONE #50616

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	10:07:41
0.100 Cal Check	0.099	10:08:02
Air Blank	0.000	10:08:31

Preventative Maintenance
Performed -- Timer Reset

5 of 7
WLL
Bh

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6 of 7
WU
By 4.10.2015

QA Specialist W. LOOMIS 53545 AGENCY Tucson Police Department
DATE 4/2/2015 TIME 09:50
Intoxilyzer Serial # 80-002261 Location ODE

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 101 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Will Loomis

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002261
Core Version: 8105.48
04/02/2015 09:50:00

Standard Lot#: 25913100A2
Last Changed By: S. CIONE #50616

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:50:52
0.100 Cal Check	0.101	09:51:13
Air Blank	0.000	09:51:42

Preventative Maintenance
Performed -- Timer Reset

7 of 7
wll
BJ