



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory

64

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002262 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 12/1/06 Cal 1 0.099 AC Cal 2 0.100 AC


The instrument was placed into service on 12/1/06.

Jennifer L. Coombs
Jennifer L. Coombs

12/27/06

ADPS QAS Permit # 33120

SUBSCRIBED AND SWORN to before me this 27, day of December, 2006.

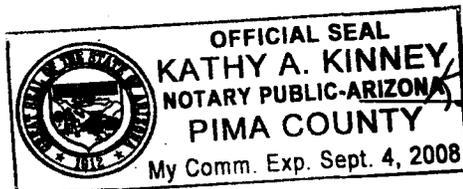

Kathy A. Kinney
Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

1 of 8
file

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE *SM*

QA SPECIALIST Jennifer Coombs 46160 AGENCY Tucson Police Department

DATE December 1, 2006 TIME 1455

INTOXILYZER SERIAL # 80-002262 LOCATION LAB

() 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- () 1. Clock time check.
- () 2. Date check.

OPERATIONAL TESTS

- () 1. Alcohol-free subject test result 0. 000 AC.
- () 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- () 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- () 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 099 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE

Jennifer L. Coombs

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002262
12/01/2006

14:56:35

BH

2 of 8
file

Standard Lot#: 502602I
Last Changed By: J. COOMBS #46160

QAS: J. COOMBS #46160
TUCSON PD CRIME LAB

Operator: J. COOMBS #46160
TUCSON PD CRIME LAB

Subject: SQAP,TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:57:36
Diagnostic Test	Pass	14:58:03
Air Blank	0.000	14:58:30
0.100 Cal Check	0.099	14:58:49
Air Blank	0.000	14:59:19
Subject Test	0.000	14:59:52
Air Blank	0.000	15:00:20
Five Minute Wait	ABT*	
Air Blank	0.000	15:00:56

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002262
12/01/2006

15:01:30

Ph

3 of 8
file

Standard Lot#: 502602I
Last Changed By: J. COOMBS #46160

QAS: J. COOMBS #46160
TUCSON PD CRIME LAB

Operator: J. COOMBS #46160
TUCSON PD CRIME LAB

Subject: SQAP, TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	IPS*	15:02:04
Air Blank	0.000	15:02:32

*Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002262
12/01/2006

15:03:12

154

4 of 8
jlc

Standard Lot#: 502602I
Last Changed By: J. COOMBS #46160

QAS: J. COOMBS #46160
TUCSON PD CRIME LAB

Operator: J. COOMBS #46160
TUCSON PD CRIME LAB

Subject: SQAP, TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	15:04:06
Diagnostic Test	Pass	15:04:34
Air Blank	0.000	15:05:01
0.100 Cal Check	RFI*	15:05:18
Air Blank	0.000	15:05:48

*RFI Detect

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002262
12/01/2006

15:06:08

BA

5 of 8
JL

Standard Lot#: 502602I
Last Changed By: J. COOMBS #46160

QAS: J. COOMBS #46160
TUCSON PD CRIME LAB

Operator: J. COOMBS #46160
TUCSON PD CRIME LAB

Subject: SQAP,TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	15:06:57
Diagnostic Test	Pass	15:07:24
Air Blank	0.000	15:07:52
0.100 Cal Check	0.099	15:08:10
Air Blank	0.000	15:08:40
Subject Test	0.137*	15:11:48
Air Blank	0.000	15:12:22
Air Blank	0.000	15:12:50
Subject Test	XXX**	15:13:23
Air Blank	0.000	15:13:59

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-002262

12/01/2006

15:15:17

6 of 8
file

54

Standard Lot#: 502602I

Last Changed By: J. COOMBS #46160

QAS: J. COOMBS #46160

TUCSON PD CRIME LAB

Test	g/210L	Time
Air Blank	0.000	15:15:44
Cal Check	0.099	15:16:03
Air Blank	0.000	15:16:33



Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

7 of 8
JL

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

SM

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Jennifer Coombs 46160 AGENCY Tucson Police Department
DATE December 1, 2006 TIME 1455
Intoxilyzer Serial # 80-002262 Location LAB

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Jennifer L. Coombs


DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-002262
12/01/2006 14:55:01

Standard Lot#: 502602I
Last Changed By: J. COOMBS #46160

QAS: J. COOMBS #46160
TUCSON PD CRIME LAB

Test	g/210L	Time
Air Blank	0.000	14:55:28
Cal Check	0.100	14:55:47
Air Blank	0.000	14:56:17

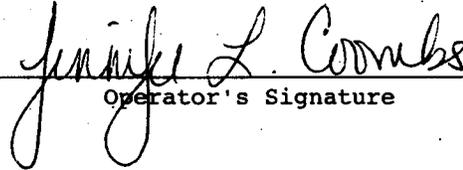
Operator's Signature



TUCSON POLICE DEPARTMENT
CRIME LABORATORY

INTOXILYZER SERVICE RECORD

1 of 1
file

PH

Intoxilyzer Serial Number: 80-002262

Initial Service Date: 12-1-06

CMI, Inc. Certificate of Calibration Date: 11-10-06