



CERTIFICATION by CUSTODIAN of RECORDS  
Tucson Police Department Crime Laboratory

ck  
5-27-08

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002262 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 5.7.2008 Cal 1 0.101 AC Cal 2 0.101 AC

Between the dates of N/A and N/A the instrument was operating properly and accurately.

The instrument was taken out of service for repair on —.  
The instrument was returned to service on 5.7.2008.

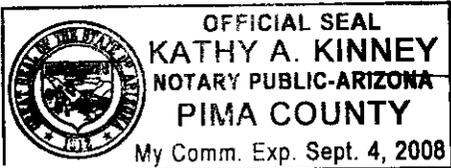
Terry S. Gallegos

Terry S. Gallegos

5.12.2008

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 12, day of May, 2008.



Kathy A. Kinney  
Notary Public

\*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

Ch  
5.27.08  
10 of 11  
154

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department

DATE 5.7.2008 TIME 1347

INTOXILYZER SERIAL # 80-002262 LOCATION Lab

(  ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (  ) 1. Clock time check.
- (  ) 2. Date check.

OPERATIONAL TESTS

- (  ) 1. Alcohol-free subject test result 0. 000 AC.
- (  ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed ✓✓✓
- (  ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- (  ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 101 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE

Terry Gallegos

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-002262  
05/07/2008

13:48:35

CU  
2 of 11  
154

Standard Lot#: 725401I  
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847  
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847  
TUCSON PD CRIME LAB

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	13:49:23
Diagnostic Test	Pass	13:49:50
Air Blank	0.000	13:50:18
0.100 Cal Check	0.101	13:50:36
Air Blank	0.000	13:51:06
Subject Test	0.000	13:51:22
Air Blank	0.000	13:51:50
Five Minute Wait	ABT*	
Air Blank	0.000	13:52:25

\*Sequence Aborted

### Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-002262  
05/07/2008

13:52:33

Ch  
3 of 11  
T84

Standard Lot#: 725401I  
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847  
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847  
TUCSON PD CRIME LAB

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	13:53:20
Diagnostic Test	Pass	13:53:47
Air Blank	IPS*	13:53:56
Air Blank	RFI**	13:54:07

\*Improper Sample

\*\*RFI Detect

### Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-002262  
05/07/2008

13:54:41

CM  
4 of 11  
134

Standard Lot#: 725401I  
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847  
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847  
TUCSON PD CRIME LAB

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	13:55:32
Diagnostic Test	Pass	13:55:59
Air Blank	0.000	13:56:27
0.100 Cal Check	0.101	13:56:45
Air Blank	0.000	13:57:15
Subject Test	0.294*	14:00:23
Air Blank	0.000	14:00:55
Air Blank	0.000	14:01:23
Subject Test	0.127	14:01:42
Air Blank	IPS**	14:02:08
Air Blank	0.000	14:02:37

\*Deficient Sample

\*\*Improper Sample

### Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-002262  
05/07/2008

14:03:18

CM  
50f 11  
TS4

Standard Lot#: 725401I  
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847  
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847  
TUCSON PD CRIME LAB

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:04:07
Diagnostic Test	Pass	14:04:34
Air Blank	0.000	14:05:02
0.100 Cal Check	0.100	14:05:20
Air Blank	0.000	14:05:50
Subject Test	XXX*	14:07:56
Air Blank	0.000	14:08:30

\*Invalid Sample  
Begin new deprivation period

### Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

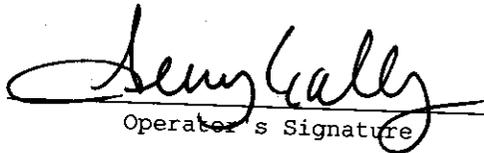
INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-002262  
05/07/2008

14:08:53

-----  
Standard Lot#: 725401I  
Last Changed By: T. GALLEGOS #28847  
-----

QAS: T. GALLEGOS #28847  
TUCSON PD CRIME LAB  
-----

Test	g/210L	Time
Air Blank	0.000	14:09:21
Cal Check	0.101	14:09:39
Air Blank	0.000	14:10:09

  
Operator's Signature

CA  
607.11  
TS4

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

CA  
5-27-08  
Fof II  
TS4

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department  
DATE 5.7.2008 TIME 1347  
Intoxilyzer Serial # 80-002262 Location Lab

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 101 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE 

DPS Form Exh G-2 (Rev. 05-01)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-002262  
05/07/2008

13:47:09

-----  
Standard Lot#: 725401I  
Last Changed By: T. GALLEGOS #28847  
-----

QAS: T. GALLEGOS #28847  
TUCSON PD CRIME LAB  
-----

Test	g/210L	Time
Air Blank	0.000	13:47:36
Cal Check	0.101	13:47:55
Air Blank	0.000	13:48:24

Ch  
80f 11  
TS4

  
Operator's Signature



TUCSON POLICE DEPARTMENT  
CRIME LABORATORY

INTOXILYZER SERVICE RECORD  
FOR  
CALIBRATION GAS TANK CHANGE

CU  
5-27-08  
9 of 11  
154

Date: 5.7.2008  
 QAS: Terry Gallegos 28947  
 Intoxilyzer Serial Number: 80-00 2262  
 New Lot Number: 725401I  
 New Expiration Date: 9.12.2009

- (✓) 1. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 2. Remove the old gas tank and install the new one.
- (✓) 3. Hit the "ESC" key twice and swipe the QAS card.
- (✓) 4. Go to the "M-Maintenance" menu and select "S-Standard Change."  
 Select "D" for dry gas and enter the following information:  
     Enter Target Value: 0.100  
     Enter Control Lot #: From the tank label  
     Enter Exp Date: From the tank label (mm/dd/yyyy)  
     Edit info Y/N? Y or N  
     Update Std Info? Y ENTER  
     "Please wait saving..."
- (✓) 5. Hit the "ESC" key twice.
- (✓) 6. Go to the "M-Maintenance" menu and select "T-Control Testing."  
 Select "D" for dry control test. After selection is made press ENTER.
- (✓) 7. Air blank completed.
- (✓) 8. Calibration check completed. Test results 0. 100 AC.
- (✓) 9. Air blank completed.
- (✓) 10. Remove and sign printed record. Attach the record to the completed checklist.

SIGNATURE Terry Gallegos

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

f

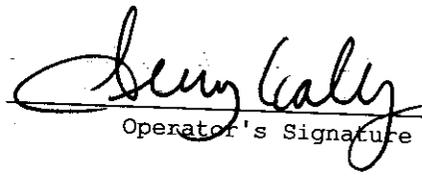
INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-002262  
05/07/2008 13:41:23

OK  
100f11  
TS4

-----  
Standard Lot#: 725401I  
Last Changed By: T. GALLEGOS #28847  
-----

QAS: T. GALLEGOS #28847  
TUCSON PD CRIME LAB  
-----

Test	g/210L	Time
Air Blank	0.000	13:41:50
Cal Check	0.100	13:42:09
Air Blank	0.000	13:42:39

  
Operator's Signature



TUCSON POLICE DEPARTMENT  
CRIME LABORATORY

INTOXILYZER SERVICE RECORD

CH  
5.27.08

11 of 11  
734

QAS: Terry Gallegos 28847

Intoxilyzer Serial Number: 80-002262

Date Returned to Service: 5.7.2008



# CERTIFICATE OF CALIBRATION

## **INTOXILYZER® 8000**

*This is to certify that the calibration of INTOXILYZER serial number 80-002262, manufactured by CMI, subsidiary of MPD, Inc. of Owensboro, Kentucky, was tested and found to be in compliance with the National Safety Administration Standard for Devices to Measure Breath Alcohol ( F.R., No. 179 48705-48710 Sept. 17, 1993) and with International Recommendation OIML R126 Edition 1998 (E). Calibration solutions are traceable to NIST (NBS) standard material 1828.*

Date 4-24-08 Signed D. Howard  
Technician

**CMI** INC.

316 EAST 9TH STREET  
OWENSBORO, KY. 42303

Pursuant to Section 28.1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.