



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory

704
8.29.08

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002262 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 8-28-08 Cal 1 0.099 AC Cal 2 0.098 AC

Between the dates of 7.31-08 and 8.28-08 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A
The instrument was returned to service on N/A


Clare Conley

Clare L. Conley

8.29.08

ADPS QAS Permit # 33314

SUBSCRIBED AND SWORN to before me this 29th day of August, 2008.

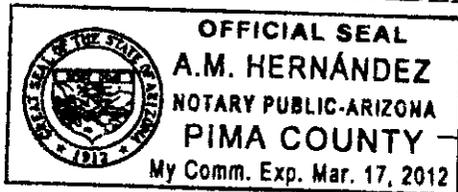

A.M. Hernandez

Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000
STANDARD QUALITY ASSURANCE PROCEDURE

1 of 7
cu
154 8.29.08

QA SPECIALIST Clare Conley AGENCY Tucson Police Department

DATE 8-28-08 TIME 0750

INTOXILYZER SERIAL # 80-002262 LOCATION ODW DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 098 AC

Instrument is operating properly and accurately. YES NO

COMMENTS None

SIGNATURE Clare Conley

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002262
08/28/2008

08:05:48

2007
cu

Standard Lot#: 735401G
Last Changed By: S. CIONE #50616

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:06:47
Diagnostic Test	Pass	08:07:14
Air Blank	0.000	08:07:42
0.100 Cal Check	0.099	08:08:00
Air Blank	0.000	08:08:30
Subject Test	0.000	08:08:48
Air Blank	ABT*	08:08:54
Air Blank	0.000	08:09:22

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002262
08/28/2008

08:09:29

Standard Lot#: 735401G
Last Changed By: S. CIONE #50616

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP

DOB: 11/11/1111

Sex: M

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	08:10:04
Air Blank	IPS**	08:10:13

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

3867
cu

SM

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002262
08/28/2008

08:10:17

Standard Lot#: 735401G
Last Changed By: S. CIONE #50616

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:11:08
Diagnostic Test	Pass	08:11:35
Air Blank	0.000	08:12:03
0.100 Cal Check	0.099	08:12:21
Air Blank	0.000	08:12:51
Subject Test	0.151*	08:15:59
Air Blank	0.000	08:16:32
Air Blank	0.000	08:17:00
Subject Test	XXX**	08:17:26
Air Blank	0.000	08:18:00

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

4567
cu
154

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002262
08/28/2008 08:18:21

Standard Lot#: 735401G
Last Changed By: S. CIONE #50616

QAS: C. CONLEY #34517
TUCSON PD

SAB
cu
104

Test	g/210L	Time
Air Blank	0.000	08:18:48
Cal Check	0.098	08:19:07
Air Blank	0.000	08:19:36

Clare Conley

Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000
STANDARD CALIBRATION CHECK PROCEDURE

6887
cu
BY 8.29.08

QA Specialist Clare Conley AGENCY Tucson Police Department
DATE 8-28-08 TIME 0750
Intoxilyzer Serial # 80-002262 Location ODW DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 099 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Clare Conley

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002262
08/28/2008


07:50:52

Standard Lot#: 735401G
Last Changed By: S. CIONE #50616

QAS: C. CONLEY #34517
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	07:51:20
Cal Check	0.099	07:51:39
Air Blank	0.000	07:52:08

767
lu
134

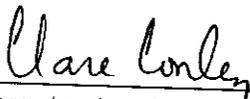

Operator's Signature

Reprinted From Database

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.