



# CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
270 S. Stone Avenue Tucson, AZ 85701-1917

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002262 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

*Handwritten:* 6-8-2009

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 6/4/09 Cal 1 0.101 AC Cal 2 0.101 AC

Between the dates of 5/7/09 and 6/4/09 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -  
The instrument was returned to service on -

Sarah Cione

Sarah L. Cione


6/4/09

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 4th day of June, 2009.

A.M. Hernandez

Notary Public

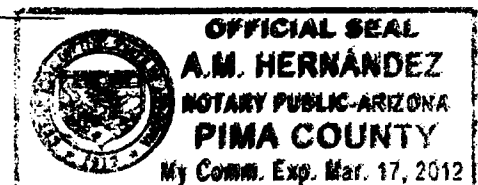

\*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

1 of 7sc

## ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000TS4  
6.8.2009

## STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police DepartmentDATE 6/4/09 TIME 10:31INTOXILYZER SERIAL # 80-002262 LOCATION DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

## DIAGNOSTIC TESTS

(✓) 1. Clock time check.

(✓) 2. Date check.

## OPERATIONAL TESTS

(✓) 1. Alcohol-free subject test result 0. 000 AC.(✓) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed(✓) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.(✓) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 101 ACInstrument is operating properly and accurately. YES ✓ NO     

## COMMENTS

---

---

---

---

SIGNATURE Sarah Cione

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-002262  
06/04/2009

10:32:42

2 of 7sc

Standard Lot#: 826001I  
Last Changed By: C. MORAND #24219

QAS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SQAP,TEST  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:33:32
Diagnostic Test	Pass	10:33:59
Air Blank	0.000	10:34:27
0.100 Cal Check	0.102	10:34:46
Air Blank	0.000	10:35:15
Subject Test	0.000	10:36:06
Air Blank	0.000	10:36:33
Five Minute Wait	ABT*	
Air Blank	0.000	10:37:19

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-002262  
06/04/2009

10:37:30

3 of 7<sub>sc</sub>

Standard Lot#: 826001I  
Last Changed By: C. MORAND #24219

QAS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SQAP, TEST  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	10:38:10
Air Blank	IPS**	10:38:18

\*RFI Detect

\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-002262  
06/04/2009

10:38:43

Standard Lot#: 826001I  
Last Changed By: C. MORAND #24219

QAS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SQAP,TEST  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:39:28
Diagnostic Test	Pass	10:39:55
Air Blank	0.000	10:40:23
0.100 Cal Check	0.101	10:40:41
Air Blank	0.000	10:41:11
Subject Test	0.101*	10:44:19
Air Blank	0.000	10:44:51
Air Blank	0.000	10:45:18
Subject Test	XXX**	10:45:36
Air Blank	0.000	10:46:10

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-002262  
06/04/2009

10:46:53

-----  
Standard Lot#: 826001I  
Last Changed By: C. MORAND #24219  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	10:47:20
Cal Check	0.101	10:47:39
Air Blank	0.000	10:48:08

*Sarah Cione*

Operator's Signature

6 of 7<sub>sc</sub>THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000134  
6.9.09

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department  
DATE 6/4/09 TIME 10:31  
Intoxilyzer Serial # 80-002262 Location DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 101 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-002262  
06/04/2009

10:31:00

-----  
Standard Lot#: 826001I  
Last Changed By: C. MORAND #24219  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	10:31:28
Cal Check	0.101	10:31:46
Air Blank	0.000	10:32:16

  
Operator's Signature



TUCSON POLICE DEPARTMENT  
CRIME LABORATORY

INTOXILYZER SERVICE RECORD  
FOR  
CALIBRATION GAS TANK CHANGE

1 of 2<sub>sc</sub>

134  
6.8.2009

Date: MAY 15 2009  
QAS: MORAND #24219  
Intoxilyzer Serial Number: 80.002262  
New Lot Number: 826001 I  
New Expiration Date: 24 SEPT 2010

- (X) 1. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (X) 2. Remove the old gas tank and install the new one.
- (X) 3. Hit the "ESC" key twice and swipe the QAS card.
- (X) 4. Go to the "M-Maintenance" menu and select "S-Standard Change." Select "D" for dry gas and enter the following information:
  - Enter Target Value: 0.100
  - Enter Control Lot #: From the tank label
  - Enter Exp Date: From the tank label (mm/dd/yyyy)
  - Edit info Y/N? Y or N
  - Update Std Info? Y ENTER
  - "Please wait saving..."
- (X) 5. Hit the "ESC" key twice.
- (X) 6. Go to the "M-Maintenance" menu and select "T-Control Testing." Select "D" for dry control test. After selection is made press ENTER.
- (X) 7. Air blank completed.
- (X) 8. Calibration check completed. Test results 0. 101 AC.
- (X) 9. Air blank completed.
- (X) 10. Remove and sign printed record. Attach the record to the completed checklist.

SIGNATURE

Cho Rymann

2 of 2sc

154

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-002262  
05/15/2009 21:54:21

Standard Lot#: 8260011  
Last Changed By: C. MORAND #24219

OAS: C. MORAND #24219  
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	21:54:48
Cal Check	0.101	21:55:07
Air Blank	0.000	21:55:37

  
Operator's Signature