



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
270 S. Stone Avenue Tucson, AZ 85701-1917

I hereby certify that these documents are a true and correct copy of the  
periodic maintenance\* for the Intoxilyzer 8000 Serial Number  
80-002262 maintained by the Tucson Police Department  
Crime Laboratory pursuant to the requirements of the Arizona Department  
of Public Safety for a quality assurance program.

TSK  
5.5.2010

A quality assurance specialist (QAS) tested the instrument on the following  
dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2,  
for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 4/28/10 Cal 1 0.100 AC Cal 2 0.099 AC

Between the dates of N/A and N/A  
the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.  
The instrument was returned to service on 4/28/10.

Sarah Cione

Sarah L. Cione

4/30/10


ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 30<sup>th</sup>, day of  
April, 20 10.

A.M. Hernandez

Notary Public

\*Includes documents specified in A.A.C. R13-10-104.A.5

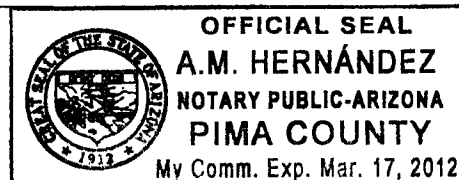



TUCSON POLICE DEPARTMENT  
CRIME LABORATORY

1 of 12 SC

INTOXILYZER SERVICE RECORD

By  
5-5-2010

QAS: Sarah Cione

Intoxilyzer Serial Number: 80-002262

Date Returned to Service: 4/28/2010

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

2 of 12 sc

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

by  
5.5.2010

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 4/28/10 TIME 1:35b

INTOXILYZER SERIAL # 80-002262 LOCATION Crime Lab

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.  
(✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.  
(✓) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed  
(✓) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.  
(✓) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 099 AC

Instrument is operating properly and accurately. YES ☒ NO ☐

COMMENTS

\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE Sarah Cione

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-002262  
04/28/2010

13:57:53

3 of 12<sub>sc</sub>

154

---

Standard Lot#: 642845  
Last Changed By: S. CIONE #50616

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	13:58:39
Diagnostic Test	Pass	13:59:06
Air Blank	0.000	13:59:33
0.100 Cal Check	0.100	13:59:52
Air Blank	0.000	14:00:22
Subject Test	0.000	14:01:00
Air Blank	0.000	14:01:27
Five Minute Wait	ABT*	
Air Blank	0.000	14:02:01

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

4 of 12<sub>SC</sub>

64

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-002262  
04/28/2010

14:02:17

---

Standard Lot#: 642845  
Last Changed By: S. CIONE #50616

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	14:03:05
Diagnostic Test	Pass	14:03:32
Air Blank	0.000	14:03:59
0.100 Cal Check	RFI*	14:04:15
Air Blank	IPS**	14:04:36

\*RFI Detect  
\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-002262  
04/28/2010

14:04:46

5 of 12<sub>sc</sub>

by

Standard Lot#: 642845  
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:05:29
Diagnostic Test	Pass	14:05:57
Air Blank	0.000	14:06:24
0.100 Cal Check	0.094	14:06:43
Air Blank	0.000	14:07:12
Subject Test	0.045*	14:10:20
Air Blank	0.000	14:10:51
Air Blank	0.000	14:11:19
Subject Test	XXX**	14:11:45
Air Blank	0.000	14:12:18

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-002262  
04/28/2010

14:12:34

-----  
Standard Lot#: 642845  
Last Changed By: S. CIONE #50616  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	14:13:01
Cal Check	0.099	14:13:20
Air Blank	0.000	14:13:49

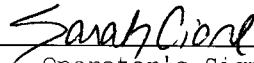
  
Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department  
DATE 4/28/10 TIME 1356  
Intoxilyzer Serial # 80-002262 Location Crime Lab

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach  $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)



INTOXILYZER 8000

Location: ARIZONA

Serial Number: 80-002262

04/28/2010

13:56:32

-----  
Standard Lot#: 642845

Last Changed By: S. CIONE #50616  
-----

QAS: S. CIONE #50616

TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	13:57:00
Cal Check	0.100	13:57:18
Air Blank	0.000	13:57:48

*Sarah Cione*

\_\_\_\_\_  
Operator's Signature



TUCSON POLICE DEPARTMENT  
CRIME LABORATORY

9 of 12 SC

INTOXILYZER SERVICE RECORD <sup>54</sup>  
FOR  
CALIBRATION GAS TANK CHANGE

Date: 4/28/2010

QAS: Sarah Cione

Intoxilyzer Serial Number: 80-002262

New Lot Number: 642845

New Expiration Date: 7/23/2011

- (✓) 1. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 2. Remove the old gas tank and install the new one.
- (✓) 3. Hit the "ESC" key twice and swipe the QAS card.
- (✓) 4. Go to the "M-Maintenance" menu and select "S-Standard Change."  
Select "D" for dry gas and enter the following information:  
Enter Target Value: 0.100  
Enter Control Lot #: From the tank label  
Enter Exp Date: From the tank label (mm/dd/yyyy)  
Edit info Y/N? Y or N  
Update Std Info? Y ENTER  
"Please wait saving..."
- (✓) 5. Hit the "ESC" key twice.
- (✓) 6. Go to the "M-Maintenance" menu and select "Control Testing."  
Select "D" for dry control test. After selection is made press ENTER.
- (✓) 7. Air blank completed.
- (✓) 8. Calibration check completed. Test results 0. 100 AC.
- (✓) 9. Air blank completed.
- (✓) 10. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

INTOXILYZER 8000

Location: ARIZONA

Serial Number: 80-002262

04/28/2010

13:55:07

10 of 12 SC  
PH

-----  
Standard Lot#: 642845

Last Changed By: S. CIONE #50616  
-----

QAS: S. CIONE #50616

TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	13:55:35
Cal Check	0.100	13:55:53
Air Blank	0.000	13:56:23



Operator's Signature

Seal # 63703		SHIPPING DATE/CHARGES	
HOURS:		SPECIAL SHIPPING INSTRUCTIONS	
TECH.	DATE	CAL. TECH.	DATE
FORM # CMI 8WO 1 REV (-)		<input type="checkbox"/> COD <input type="checkbox"/> AMOUNT \$	

12 of 12 sc  
154

## Certificate of Calibration

### **Intoxilyzer® 8000**

*This is to certify that the calibration of **Intoxilyzer®** serial number 80-002262, manufactured by CMI Inc., a subsidiary of MPD, Inc. of Owensboro, Kentucky, was tested and found to be in compliance with the National Highway Traffic Safety Administration Standard for Devices to Measure Breath Alcohol (F.R., No. 179 48705-48710 Sept. 17, 1993) and with International Recommendation OIML R126 Edition 1998 (E). Calibration solutions are traceable to NIST standard material 1828.*

Date 4/14/10 Signed Kristi Robinson  
Technician

**CMI** INC.

316 East 9<sup>th</sup> Street  
Owensboro, KY 42303  
USA

Part No. 650519 CMI 7/9/09

Pursuant to Section 28.1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.