



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

TS4
8-23-2010

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002262 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 8/20/10 Cal 1 0.099 AC Cal 2 0.099 AC

Between the dates of 7/22/10 and 8/20/10 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -
The instrument was returned to service on -

Sarah Cione

Sarah L. Cione


8/23/10

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 23 day of August, 2010.

A.M. Hernandez

Notary Public

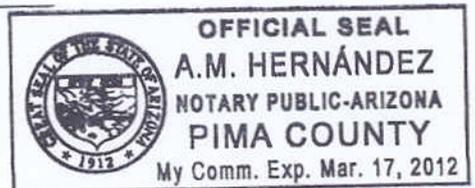

*Includes documents specified in A.A.C. R13-10-104.A.5

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

1 of 7 sc

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000
STANDARD QUALITY ASSURANCE PROCEDURE

TS4
8.23.2010

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 8/20/10 TIME 1311

INTOXILYZER SERIAL # 80-002262 LOCATION Crime Lab

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0.000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0.100 AC
Result: 0.099 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Sarah Cione

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

2 of 7sc

TS4

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002262
08/20/2010 13:12:33

Standard Lot#: 642845
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	13:13:14
Diagnostic Test	Pass	13:13:42
Air Blank	0.000	13:14:09
0.100 Cal Check	0.099	13:14:28
Air Blank	0.000	13:14:57
Subject Test	0.000	13:15:12
Air Blank	0.000	13:15:39
Five Minute Wait	ABT*	
Air Blank	0.000	13:16:14

*Sequence Aborted

Not a Successfully Completed Test Sequence

3 of 7_{sc}

TS4

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002262
08/20/2010 13:16:22

Standard Lot#: 642845
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	13:17:04
Diagnostic Test	Pass	13:17:32
Air Blank	RFI*	13:17:57
Air Blank	IPS**	13:18:06

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

4 of 7sc

TS4

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002262
08/20/2010 13:18:11

Standard Lot#: 642845
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	13:18:53
Diagnostic Test	Pass	13:19:20
Air Blank	0.000	13:19:48
0.100 Cal Check	0.099	13:20:06
Air Blank	0.000	13:20:36
Subject Test	0.021*	13:23:44
Air Blank	0.000	13:24:15
Air Blank	0.000	13:24:42
Subject Test	XXX**	13:25:07
Air Blank	0.000	13:25:41

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

5 of 7 sc
B4

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002262
08/20/2010 13:26:01

Standard Lot#: 642845
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	13:26:28
Cal Check	0.099	13:26:47
Air Blank	0.000	13:27:16

Sarah Cione
Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

6 of 7 sc

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000
STANDARD CALIBRATION CHECK PROCEDURE

TS4
8.23.2010

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 8/20/10 TIME 1311
Intoxilyzer Serial # 80-002262 Location Crime Lab

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 099 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

7 of 7 sc

AS4

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002262
08/20/2010 13:11:07

Standard Lot#: 642845
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	13:11:34
Cal Check	0.099	13:11:53
Air Blank	0.000	13:12:22

Sarah Cione

Operator's Signature