



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

TS4
12-13-2010

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002262 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 12-9-10 Cal 1 0.100 AC Cal 2 0.100 AC

Between the dates of 11-11-10 and 12-9-10 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.
The instrument was returned to service on N/A.


Clare Conley
Clare L. Conley

12-10-10

ADPS QAS Permit # 33314

SUBSCRIBED AND SWORN to before me this 10th day of December, 2010.

Rosa Reyna
Notary Public

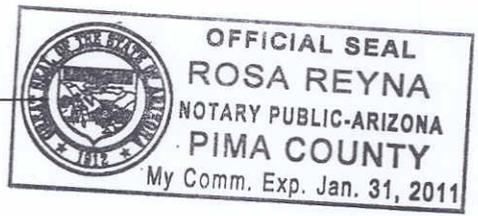

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1567
CM

PS
12-13-2010

QA SPECIALIST Clare Conley AGENCY Tucson Police Department

DATE 12-9-10 TIME 0642

INTOXILYZER SERIAL # 80-002262 LOCATION Crime Lab

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 100 AC

Instrument is operating properly and accurately. YES NO

COMMENTS None.

SIGNATURE Clare Conley

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002262
12/09/2010

06:44:28

Standard Lot#: 642845
Last Changed By: S. CIONE #50616

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	06:45:15
Diagnostic Test	Pass	06:45:42
Air Blank	0.000	06:46:09
0.100 Cal Check	0.100	06:46:28
Air Blank	0.000	06:46:58
Subject Test	0.000	06:47:44
Air Blank	ABT*	06:47:50
Air Blank	0.000	06:48:18

*Sequence Aborted

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

2567
m
754

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002262
12/09/2010

06:48:33

357
cu
PH

Standard Lot#: 642845
Last Changed By: S. CIONE #50616

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	06:49:04
Air Blank	IPS**	06:49:16

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002262
12/09/2010

06:49:30

487
cu

Standard Lot#: 642845
Last Changed By: S. CIONE #50616

124

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	06:50:19
Diagnostic Test	Pass	06:50:46
Air Blank	0.000	06:51:14
0.100 Cal Check	0.101	06:51:33
Air Blank	0.000	06:52:02
Subject Test	0.145*	06:55:10
Air Blank	0.000	06:55:42
Air Blank	0.000	06:56:09
Subject Test	XXX**	06:56:26
Air Blank	0.000	06:57:01

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002262
12/09/2010 06:57:28

Standard Lot#: 642845
Last Changed By: S. CIONE #50616

QAS: C. CONLEY #34517
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	06:57:56
Cal Check	0.100	06:58:15
Air Blank	0.000	06:58:44

50617
cu
TS4

Clare Conley
Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6587
lu
154
12-13-2010

QA Specialist Clare Conley AGENCY Tucson Police Department
DATE 12-9-10 TIME 0642
Intoxilyzer Serial # 80-002262 Location Crime Lab

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Clare Conley

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002262
12/09/2010 06:42:40

Standard Lot#: 642845
Last Changed By: S. CIONE #50616

QAS: C. CONLEY #34517
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	06:43:08
Cal Check	0.100	06:43:26
Air Blank	0.000	06:43:56

Clare Conley

Operator's Signature

7567
OK