



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
270 S. Stone Avenue Tucson, AZ 85701-1917

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002262 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

B4  
3.31.2011

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 3-31-11 Cal 1 0.101 AC Cal 2 0.100 AC

Between the dates of 3-3-11 and 3-31-11 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.  
The instrument was returned to service on N/A.

Clare Conley

Clare L. Conley


3-31-11

ADPS QAS Permit # 33314

SUBSCRIBED AND SWORN to before me this 31, day of March, 2011.

Laura B. Hodgson

Notary Public

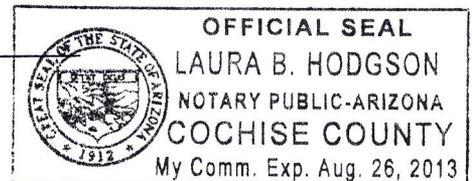

\*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

1567  
lu

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

By  
3.31.2011

QA SPECIALIST Clare Conky AGENCY Tucson Police Department

DATE 3-31-11 TIME 0808

INTOXILYZER SERIAL # 80-002262 LOCATION ODW DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 100 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS None.  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE Clare Conky

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-002262  
03/31/2011

08:10:17

2867  
en

Standard Lot#: 642845  
Last Changed By: S. CIONE #50616

B4

QAS: C. CONLEY #34517  
TUCSON PD

Operator: C. CONLEY #34517  
TUCSON PD

Subject: TEST, SQAP  
DOB: 11/11/1111  
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:11:05
Diagnostic Test	Pass	08:11:32
Air Blank	0.000	08:12:00
0.100 Cal Check	0.100	08:12:18
Air Blank	0.000	08:12:48
Subject Test	0.000	08:13:23
Air Blank	ABT*	08:13:30
Air Blank	0.000	08:13:58

\*Sequence Aborted

### Not a Successfully Completed Test Sequence

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-002262  
03/31/2011 08:14:46

3567  
lu

Standard Lot#: 642845  
Last Changed By: S. CIONE #50616

B4

QAS: C. CONLEY #34517  
TUCSON PD

Operator: C. CONLEY #34517  
TUCSON PD

Subject: TEST, SQAP  
DOB: 11/11/1111  
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	08:15:15
Air Blank	IPS**	08:15:23

\*RFI Detect  
\*\*Improper Sample

### Not a Successfully Completed Test Sequence

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-002262  
03/31/2011

08:15:43

4567  
Cu

Standard Lot#: 642845  
Last Changed By: S. CIONE #50616

TS4

QAS: C. CONLEY #34517  
TUCSON PD

Operator: C. CONLEY #34517  
TUCSON PD

Subject: TEST, SQAP  
DOB: 11/11/1111  
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:16:29
Diagnostic Test	Pass	08:16:57
Air Blank	0.000	08:17:24
0.100 Cal Check	0.099	08:17:43
Air Blank	0.000	08:18:12
Subject Test	0.087*	08:21:21
Air Blank	0.000	08:21:52
Air Blank	0.000	08:22:19
Subject Test	XXX**	08:22:50
Air Blank	0.000	08:23:23

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

### Not a Successfully Completed Test Sequence

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-002262  
03/31/2011 08:23:43

5067  
cu

-----  
Standard Lot#: 642845  
Last Changed By: S. CIONE #50616  
-----

BH

-----  
QAS: C. CONLEY #34517  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	08:24:10
Cal Check	0.100	08:24:29
Air Blank	0.000	08:24:58

*Claire Conley*

\_\_\_\_\_  
Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6567  
Cu

TBY  
3.31.2011

QA Specialist Clare Conley AGENCY Tucson Police Department  
DATE 3-31-11 TIME 0808  
Intoxilyzer Serial # SD-002262 Location ODW DW #

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 101 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Clare Conley

DPS Form Exh G-2 (Rev. 05-01)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-002262  
03/31/2011 08:08:30

787  
en

-----  
Standard Lot#: 642845  
Last Changed By: S. CIONE #50616  
-----

104

-----  
QAS: C. CONLEY #34517  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	08:08:57
Cal Check	0.101	08:09:16
Air Blank	0.000	08:09:45

*Clare Conley*

\_\_\_\_\_  
Operator's Signature