



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

TS4
9.27.2011

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002262 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 9-15-11 Cal 1 0.098 AC Cal 2 0.096 AC

Between the dates of 8-18-11 and 9-15-11 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.
The instrument was returned to service on N/A.


Clare Conley

Clare L. Conley

9-23-11

ADPS QAS Permit # 33314

SUBSCRIBED AND SWORN to before me this 23, day of

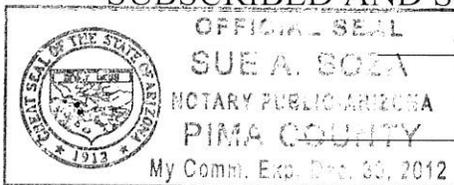

September, 20 11.

Sue A Soza

Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

1 of 8
em

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

TS4
9-27-2011

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Clare Conley AGENCY Tucson Police Department

DATE 9-15-11 TIME 0814

INTOXILYZER SERIAL # 80-002262 LOCATION ODW DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

(✓) 1. Clock time check.

(✓) 2. Date check.

OPERATIONAL TESTS

(✓) 1. Alcohol-free subject test result 0. 000 AC.

(✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed

(✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.

(✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 096 AC

Instrument is operating properly and accurately. YES NO

COMMENTS None.

SIGNATURE Clare Conley

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002262
09/15/2011

08:16:26

2/8
lu
154

Standard Lot#: 20910100A1
Last Changed By: S. CIONE #50616

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST,SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:17:17
Diagnostic Test	Pass	08:17:44
Air Blank	0.000	08:18:12
0.100 Cal Check	0.098	08:18:31
Air Blank	0.000	08:19:00
Subject Test	0.000	08:19:34
Air Blank	ABT*	08:19:40
Air Blank	0.000	08:20:08

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002262
09/15/2011

08:20:34

3.18
em
TS4

Standard Lot#: 20910100A1
Last Changed By: S. CIONE #50616

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:21:22
Diagnostic Test	Pass	08:21:49
Air Blank	0.000	08:22:17
0.100 Cal Check	0.098	08:22:35
Air Blank	0.000	08:23:05
Subject Test	RFI*	08:23:36
Air Blank	IPS**	08:23:44

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002262
09/15/2011

08:24:11

468
en

Standard Lot#: 20910100A1
Last Changed By: S. CIONE #50616

B4

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:24:57
Diagnostic Test	Pass	08:25:24
Air Blank	0.000	08:25:51
0.100 Cal Check	0.098	08:26:10
Air Blank	0.000	08:26:39
Subject Test	0.302*	08:29:48
Air Blank	0.000	08:30:21
Air Blank	0.000	08:30:48
Subject Test	RGE**	08:31:42
Air Blank	0.000	08:32:18

*Deficient Sample

**Range Exceeded

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002262
09/15/2011 08:32:26

5068
CM

TS4

Standard Lot#: 20910100A1
Last Changed By: S. CIONE #50616

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:33:14
Diagnostic Test	Pass	08:33:41
Air Blank	0.000	08:34:08
0.100 Cal Check	0.097	08:34:27
Air Blank	0.000	08:34:56
Subject Test	XXX*	08:35:25
Air Blank	0.000	08:36:00

*Invalid Sample
Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002262
09/15/2011 08:36:41

6.78
cu
754

Standard Lot#: 20910100A1
Last Changed By: S. CIONE #50616

QAS: C. CONLEY #34517
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	08:37:08
Cal Check	0.096	08:37:27
Air Blank	0.000	08:37:56



Operator's Signature

EXHIBIT G-2

7568
Cu

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

TS4
9.27.2011

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Clare Conley AGENCY Tucson Police Department
DATE 9-15-11 TIME 0814
Intoxilyzer Serial # 80-002262 Location ODW DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 098 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Clare Conley

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002262
09/15/2011 08:14:52

8/28
en

Standard Lot#: 20910100A1
Last Changed By: S. CIONE #50616

154

QAS: C. CONLEY #34517
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	08:15:20
Cal Check	0.098	08:15:38
Air Blank	0.000	08:16:08

Clare Conley

Operator's Signature