



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

TS4
11.7-2011

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002262 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 10/14/11 Cal 1 0.098 AC Cal 2 0.097 AC

Between the dates of 9/15/11 and 10/14/11 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on —.
The instrument was returned to service on —.

Sarah Cione

Sarah L. Cione


10/27/11

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 27, day of October, 20 11.

A.M. Hernández

Notary Public

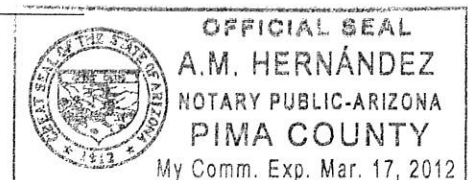

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

TS4
11.7.2011

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 10/14/11 TIME 1008

INTOXILYZER SERIAL # 80-002262 LOCATION DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
(✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
(✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
(✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
(✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 097 AC

Instrument is operating properly and accurately. YES ✓ NO

COMMENTS

SIGNATURE

Sarah Cione

2 of 7_{sc}

TS4

INTOXILYZER 8000
 Location: ARIZONA
 Serial Number : 80-002262
 10/14/2011

10:10:54

Standard Lot#: 20910100A1
 Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
 TUCSON PD

Operator: S. CIONE #50616
 TUCSON PD

Subject: SQAP, SQAP
 DOB: 11/11/1111
 Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:11:41
Diagnostic Test	Pass	10:12:08
Air Blank	0.000	10:12:36
0.100 Cal Check	0.099	10:12:54
Air Blank	0.000	10:13:24
Subject Test	0.000	10:13:46
Air Blank	0.000	10:14:13
Five Minute Wait	ABT*	
Air Blank	0.000	10:14:49

*Sequence Aborted

Not a Successfully Completed Test Sequence

3 of 7_{SC}

B4

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002262
10/14/2011 10:16:15

Standard Lot#: 20910100A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111
15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	10:16:49
Air Blank	IPS**	10:16:57

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

4 of 7sc

154

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002262
10/14/2011

10:18:32

Standard Lot#: 20910100A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:19:20
Diagnostic Test	Pass	10:19:47
Air Blank	0.000	10:20:14
0.100 Cal Check	0.098	10:20:33
Air Blank	0.000	10:21:03
Subject Test	0.034*	10:24:11
Air Blank	0.000	10:24:43
Air Blank	0.000	10:25:10
Subject Test	XXX**	10:25:57
Air Blank	0.000	10:26:32

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000

Location: ARIZONA

Serial Number: 80-002262

10/14/2011

10:26:57

Standard Lot#: 20910100A1

Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	10:27:24
Cal Check	0.097	10:27:43
Air Blank	0.000	10:28:12

Sarah Cione

Operator's Signature

5 of 7sc

TS4

EXHIBIT G-2

6 of 7sc

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

B4
11.7.2011

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 10/14/11 TIME 1008
Intoxilyzer Serial # 80-002262 Location DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 098 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002262
10/14/2011

10:08:37

Standard Lot#: 20910100A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	10:09:04
Cal Check	0.098	10:09:23
Air Blank	0.000	10:09:52

Sarah Cione

Operator's Signature