



CERTIFICATION by CUSTODIAN of RECORDS
 Tucson Police Department Crime Laboratory
 1306 W. Miracle Mile, Tucson, AZ 85705

CM
 11-30-11

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number PC-002262 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 11/10/11 Cal 1 0.099 AC Cal 2 0.098 AC

Between the dates of 10/14/11 and 11/10/11 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.
 The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione


11/16/11

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 16, day of November, 2011.

A.M. Hernandez

Notary Public

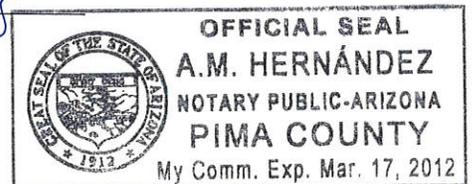

*Includes documents specified in A.A.C. R13-10-104.A.5

1077 SC

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

CM
11-30-11

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 11/10/11 TIME 1550

INTOXILYZER SERIAL # 80-002262 LOCATION ODS-Vehicle #944

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 098 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Sarah Cione

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

2 of 7 sc

OM

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002262
11/10/2011 15:52:07

Standard Lot#: 20910100A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: TEST,TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	15:53:05
Diagnostic Test	Pass	15:53:34
Air Blank	0.000	15:54:01
0.100 Cal Check	0.100	15:54:20
Air Blank	0.000	15:54:49
Subject Test	0.000	15:55:06
Air Blank	0.000	15:55:33
Five Minute Wait	ABT*	
Air Blank	0.000	15:56:07

*Sequence Aborted

Not a Successfully
Completed
Test Sequence

3 of 7 sc

CM

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002262
11/10/2011 15:56:14

Standard Lot#: 20910100A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RF1*	15:56:48
Air Blank	IPS**	15:56:56

*RF1 Detect
**Improper Sample

Not a Successfully
Completed
Test Sequence

4 of 7sc

CM

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002262
11/10/2011 15:57:04

Standard Lot#: 20910100A1
Last Changed By: S. CIONE #50616

DAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SOAP, SOAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	15:57:58
Diagnostic Test	Pass	15:58:26
Air Blank	0.000	15:58:54
0.100 Cal Check	0.100	15:59:13
Air Blank	0.000	15:59:42
Subject Test	0.139*	16:02:50
Air Blank	0.000	16:03:23
Air Blank	0.000	16:03:50
Subject Test	XXX**	16:04:07
Air Blank	0.000	16:04:41

*Deficient Sample
**Invalid Sample
Begin new deprivation period

Not a Successfully
Completed
Test Sequence

5 of 7 sc

CH

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002262
11/10/2011 16:04:59

Standard Lot#: 20910100A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	16:05:27
Cal Check	0.098	16:05:45
Air Blank	0.000	16:06:15

Sarah Cione
Operator's Signature

EXHIBIT G-2

6 of 7sc

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

CM
11-30-11

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 11/10/11 TIME 1550
Intoxilyzer Serial # 80-002262 Location ODS - Vehicle #944

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 099 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

7 of 7sc

cm

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002262
11/10/2011 15:50:43

Standard Lot#: 20910100A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	15:51:10
Cal Check	0.099	15:51:29
Air Blank	0.000	15:51:58

Sarah Cione
Operator's Signature