



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002262 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 3.1.2012 Cal 1 0.098 AC Cal 2 0.100 AC

Between the dates of 2-3.2012 and 3.1.2012 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on [blank]
The instrument was returned to service on [blank]

[Handwritten signature of Terry S. Gallegos]

Terry S. Gallegos


3.26.2012

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 26th, day of March, 2012.

[Handwritten signature of Edna Escarcega]

Notary Public

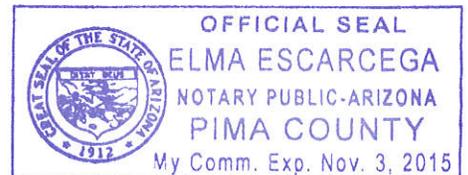

\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

10 of 7  
TS4

SC  
3/27/12

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department  
DATE 3.1.2012 TIME 1053  
INTOXILYZER SERIAL # 80-002262 LOCATION DUI HQ

( ✓ ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- ( ✓ ) 1. Clock time check.
- ( ✓ ) 2. Date check.

OPERATIONAL TESTS

- ( ✓ ) 1. Alcohol-free subject test result 0. 000 AC.
- ( ✓ ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- ( ✓ ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- ( ✓ ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 098 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE

Terry Gallegos

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-002262  
03/01/2012 10:58:26

2077  
B4  
SC

Standard Lot#: 20910100A1  
Last Changed By: S. CIONE #50616

QAS: T. GALLEGOS #28847  
TUCSON PD

Operator: T. GALLEGOS #28847  
TUCSON PD

Subject: SQAP,TEST  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:59:26
Diagnostic Test	Pass	10:59:53
Air Blank	0.000	11:00:21
0.100 Cal Check	0.099	11:00:39
Air Blank	0.000	11:01:09
Subject Test	0.000	11:01:25
Air Blank	ABT*	11:01:43
Air Blank	0.000	11:02:11

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-002262  
03/01/2012 11:03:53

3087  
T54  
SC

---

Standard Lot#: 20910100A1  
Last Changed By: S. CIONE #50616

---

QAS: T. GALLEGOS #28847  
TUCSON PD

---

Operator: T. GALLEGOS #28847  
TUCSON PD

---

Subject: SQAP,TEST  
DOB: 11/11/1111  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	11:04:44
Diagnostic Test	Pass	11:05:11
Air Blank	RFI*	11:05:18
Air Blank	IPS**	11:05:29

\*RFI Detect  
\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-002262  
03/01/2012 11:08:26

4077

B4

5C

Standard Lot#: 20910100A1  
Last Changed By: S. CIONE #50616

QAS: T. GALLEGOS #28847  
TUCSON PD

Operator: T. GALLEGOS #28847  
TUCSON PD

Subject: SQAP, TEST  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:09:17
Diagnostic Test	Pass	11:09:44
Air Blank	0.000	11:10:11
0.100 Cal Check	0.100	11:10:30
Air Blank	0.000	11:10:59
Subject Test	0.123*	11:14:08
Air Blank	0.000	11:14:39
Air Blank	0.000	11:15:07
Subject Test	XXX**	11:15:28
Air Blank	0.000	11:16:03

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-002262  
03/01/2012 11:16:44

Sof 7  
TS4

-----  
Standard Lot#: 20910100A1  
Last Changed By: S. CIONE #50616  
-----

SC

-----  
QAS: T. GALLEGOS #28847  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	11:17:11
Cal Check	0.098	11:17:30
Air Blank	0.000	11:17:59

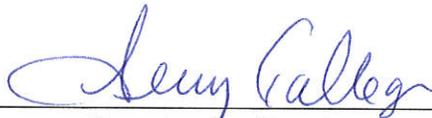
  
-----  
Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY


STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6 of 7  
TS4  
SC  
3/27/12

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department  
DATE 3.1.2012 TIME 1053  
Intoxilyzer Serial # 80-002262 Location DUI HQ

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 100 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

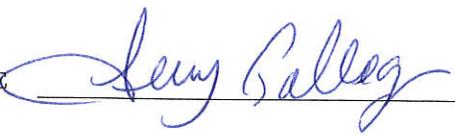SIGNATURE 

DPS Form Exh G-2 (Rev. 05-01)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-002262  
03/01/2012 10:53:48

7077  
TS4

SC

-----  
Standard Lot#: 20910100A1  
Last Changed By: S. CIONE #50616  
-----

QAS: T. GALLEGOS #28847  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	10:54:15
Cal Check	0.100	10:54:34
Air Blank	0.000	10:55:04



Operator's Signature